

Integration and Ruggedization of a Commercially Available Gas Chromatograph and Mass Spectrometer (GCMS) for the Resource Prospector Mission (RPM)

Presentation to HEMS
September , 2013



NASA
KSC
JPL

Inficon
OI Analytical
Creare



K. Loftin, J.. Captain, T. Griffin, NASA/KSC
R.Kidd, E.Niedholdt NASA/JPL
G. Kibelka OI Analytical
K. Wright. Inficon



Presentation Overview



- **Mission Summary**
- **Preliminary Payload Requirements**
- **Analog field tests Engineering Test Unit Design**
- **Challenges (Budget, Intellectual Property, Technical)**
- **Flight Forward Plans**



Resource Prospector Mission (RPM)



RPM is a potential low cost mission to the moon comprising a soft lander, and a rover carrying a payload designed to detect and map volatiles and for demonstration of in situ resource utilization (production of oxygen from regolith). The concept of operations includes landing in an area of limited solar illumination in the vicinity of the lunar south pole where volatiles may be trapped in the sub-surface.



Why do we need a Resource Prospector Mission (RPM)

To maintain a long term human presence in space , we must learn to use the resources that are available and not rely completely on transporting all of our supplies in other words”

We have to learn to “Live off the Land”

Imagine how far Lewis and Clark would have gotten if they had to rely on all the water, food and shelter they packed in.



What is RESOLVE?

Regolith & Environment Science and Oxygen & Lunar Volatile Extraction



RESOLVE is an internationally developed payload (NASA and CSA) that that can perform two important missions for Science and Human Exploration of the Moon

Resource Prospecting Mission: (Polar site)

- ✓ **Verify the existence of and characterize the constituents and distribution of water and other volatiles in lunar polar surface materials**
 - Map the surface distribution of hydrogen rich materials
 - Determine the mineral/chemical properties of polar regolith
 - Measure bulk properties & extract core sample from selected sites
 - To a depth of 1m with minimal loss of volatiles
 - Heat multiple samples from each core to drive off volatiles for analysis
 - From <100K to 423 K (150°C)
 - From 0 up to 100 psia (reliably seal in aggressively abrasive lunar environment)
 - Determine the constituents and quantities of the volatiles extracted
 - Quantify important volatiles: H_2 , He, CO, CO_2 , CH_4 , H_2O , N_2 , NH_3 , H_2S , SO_2 ,
 - Survive limited exposure to HF, HCl, and Hg

ISRU Processing Demonstration Mission: (Equatorial and/or Polar Site)

- ✓ **Demonstrate the Hydrogen Reduction process to extract oxygen from lunar regolith**
 - Heat sample to reaction temperature
 - From 423 K (150°C) to 1173 K (900°C)
 - Flow H_2 through regolith to extract oxygen in the form of water
 - Capture, quantify, and display the water generated



RESOLVE Top-Level Science/ISRU Requirements



- **Measure the water and hydrogen bearing volatiles content in a lunar subsurface**
 - **Determine the horizontal and vertical distribution of hydrogen and hydrogen bearing volatiles**
 - Measure the spatial distribution of hydrogen and hydrogen bearing volatiles with a min. horizontal resolution of 2.0 meters and a minimum vertical resolution of 0.25 m.
 - Measure neutron flux consistent with at least 0.5 wt% water equivalent hydrogen,
 - Measure 0.5 wt% water-equivalent layer at 1 meter depth under a dry overburden, with 25 cm depth resolution
 - Measure spatial OH and H₂O in the Near IR spectrum
 - **Extract subsurface material**
 - Extract a subsurface sample up to a depth of 0.75 m. (Goal of 1 m)
 - Maintain a minimum of 1 subsurface core segment per core below 175 deg K (-98 deg C)
 - Selectively accept sections of an acquired subsurface sample
 - Obtain augured cuttings from a depth of 0.5 meters
 - **Measure the abundance of water and hydrogen bearing volatiles in the lunar subsurface**
 - Quantify water in the lunar regolith when water concentrations are between 0.5% to 8.0% (95% TBR) by mass
 - Process a minimum of 40 subsurface core segments.
 - Heat lunar regolith samples to a minimum of 425 deg K (150 deg C) for volatile extraction.
 - Identify and measure the relative abundance of the volatile constituents of the lunar regolith below 70 amu
 - Measure the isotope ratio of Deuterium/Hydrogen and Oxygen 16/18
- **Measure geotechnical characteristics of the lunar highlands and cold traps**
 - Measure the distribution of grains in the lunar regolith with respect to size and shape. (GOAL)
 - Measure bulk characteristics of lunar regolith
 - Determining geotechnical parameters of the drilling media during the sample acquisition phase
 - Identifying mineralogical features in the lunar regolith
- **Demonstrate oxygen extraction from regolith using the Hydrogen Reduction process**
 - Heat samples to 1175 K (902 C) to hydrogen reduction
 - Measure water vapor produced
 - Image water condensate/droplets produced during volatile analysis and H₂ reduction



RESOLVE Development Toward Flight



Internal Call for Proposal Awarded 2006

Performed Lunar Polar Design Ref Mission

Lab. Tests

1st & 2nd ISRU
Analog Test

3rd ISRU
Analog Test

Lunar Env.
Chamber Test

Gen I

2006-2007

Demonstrate Feasibility & Subsystem Performance

- Hardware designed to demonstrate functions needed for RESOLVE
- Minimal integration between functions
- Minimal software and autonomous control development
- No mission operation considerations
- Not considered in design:
 - Flight environment
 - Mass, power, and volume
 - Mission operations



Gen II

2007-2008

Demonstrate Integration & Operations

- Hardware designed to demonstrate functions needed for RESOLVE in one 'flight like' package
- Flight mass and volume for RESOLVE functions considered in design
- Start of software and control development
- Start of mission operation considerations
- Not considered in design:
 - Flight environment
 - Flight-like avionics and power conditioning, and ground support hardware



Gen IIIA

2010-2012

Develop 'flight like' unit for mission simulation

- Hardware designed to address lunar polar ice/volatile mission requirements
- Software and control development
- Focus on mission operations
- Design to operate under lab/analog conditions with path to lunar env.



Gen IIIB

2012-2014

Develop flight prototype for vacuum operation

- Hardware designed to operate under lunar conditions
- Focus on
 - Flight design for all RESOLVE hardware
 - Software & control of hardware operation
 - Mission operation timeline and power profile
 - Environment: vacuum, temperatures, EMI, materials




Test Site on Hawaii Very Much Like the Moon!





RESOLVE Analog Field Tests

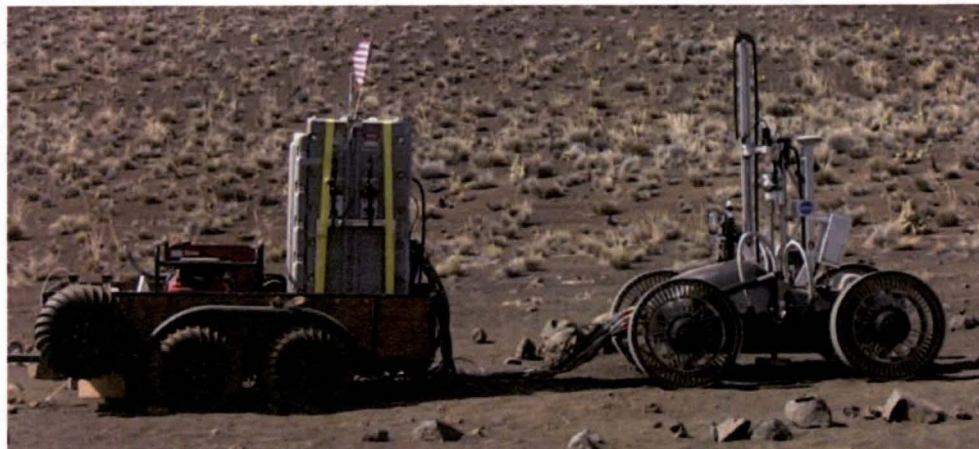


Nov. 2008

- **RESOLVE Gen II on Scarab Rover**
- Power, avionics, and ground support equipment on separate trailer




FEB. 2010

- **RESOLVE Gen II+ on CSA Juno Rover**
- Power, avionics, and ground support equipment on separate Juno

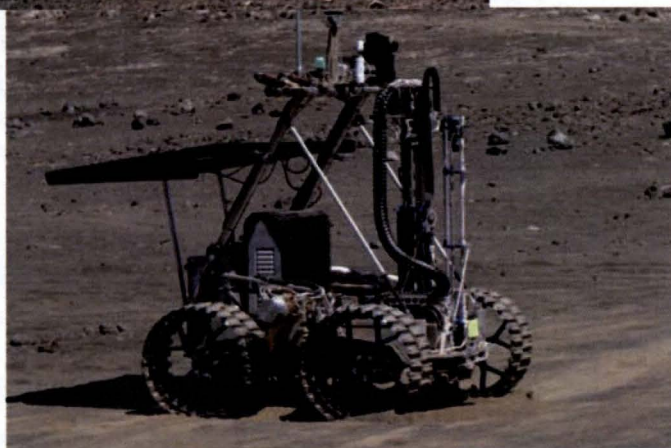

July 2012

- **RESOLVE Gen IIIA on CSA Artemis Jr. Rover**
- Everything on single rover platform

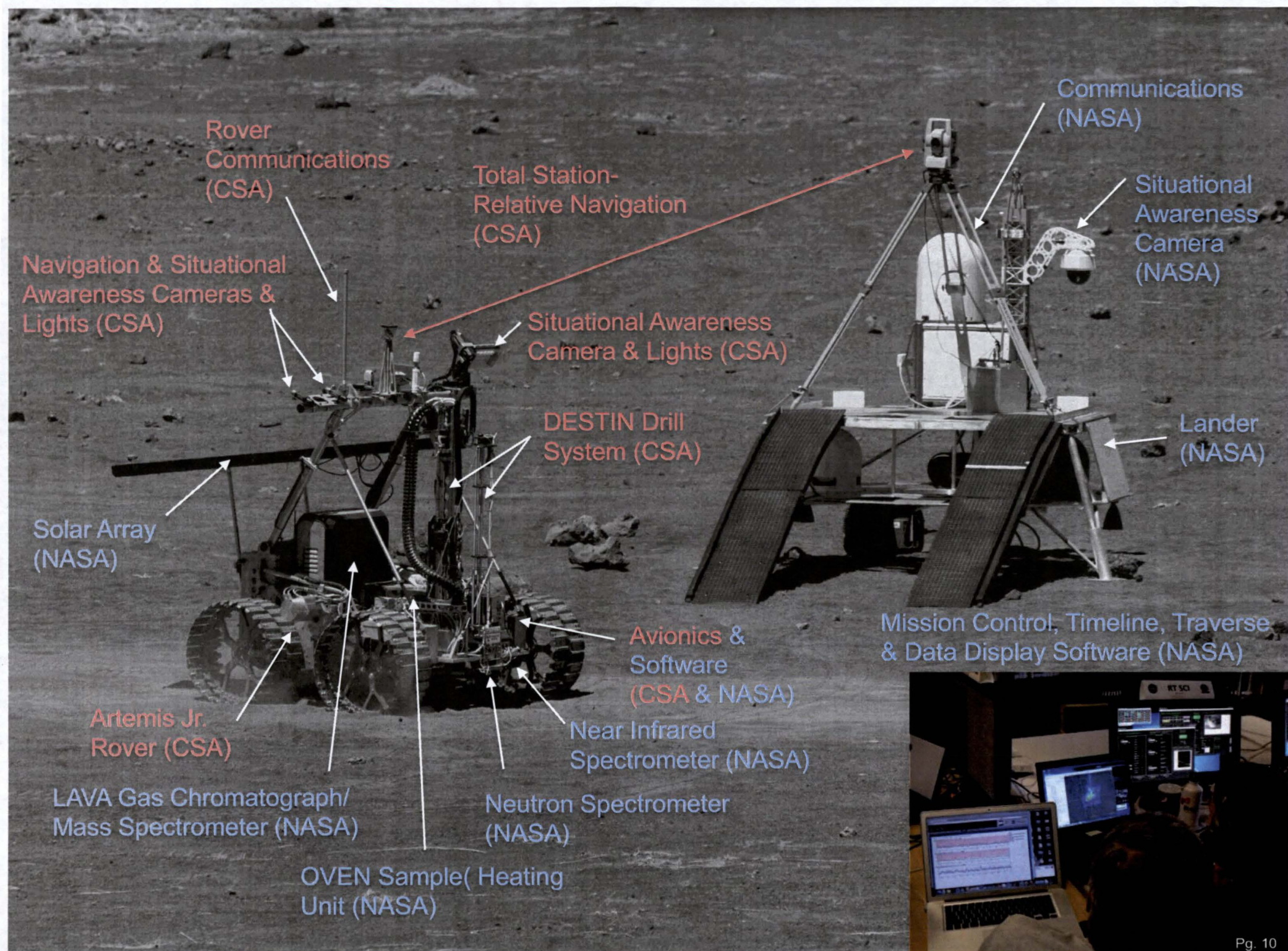



Neutron Spectrometer Subsystem (NSS) Functions & Design Constraints




NSS Functions

- Determine the horizontal and vertical distribution of hydrogen and hydrogen bearing compounds/minerals
 - Map to 1 meter depth and ~ 1 m wide path
 - Map to depth at rover speed $\leq 10 \pm 1$ cm/s
- Detect water at a minimum abundance of 0.5% by mass with $<10\%$ uncertainty
- Operate ≥ 6 hrs in permanently shadowed area

NSS Design Constraints


- Mount ~1 m above the surface aimed in front of the rover
- Operate -30 to $+40$ °C
- Max temperature change rate: 20 °C /hour
- Instrument Mass: 1.85 kg
- Power: 2 W ave.; <4 W max for heaters

NSS Sensor Module (SM)
On boom or fixed structure

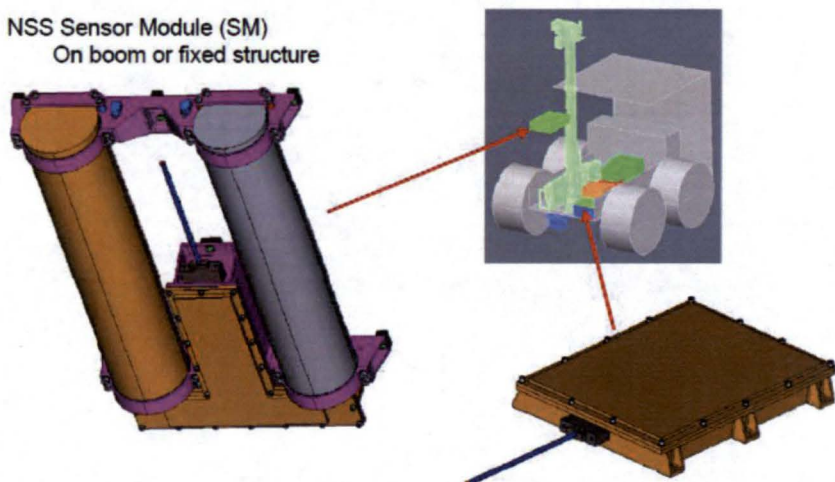

NSS Data Processing Module (DPM)
Integrated within main payload

NSS Brassboard

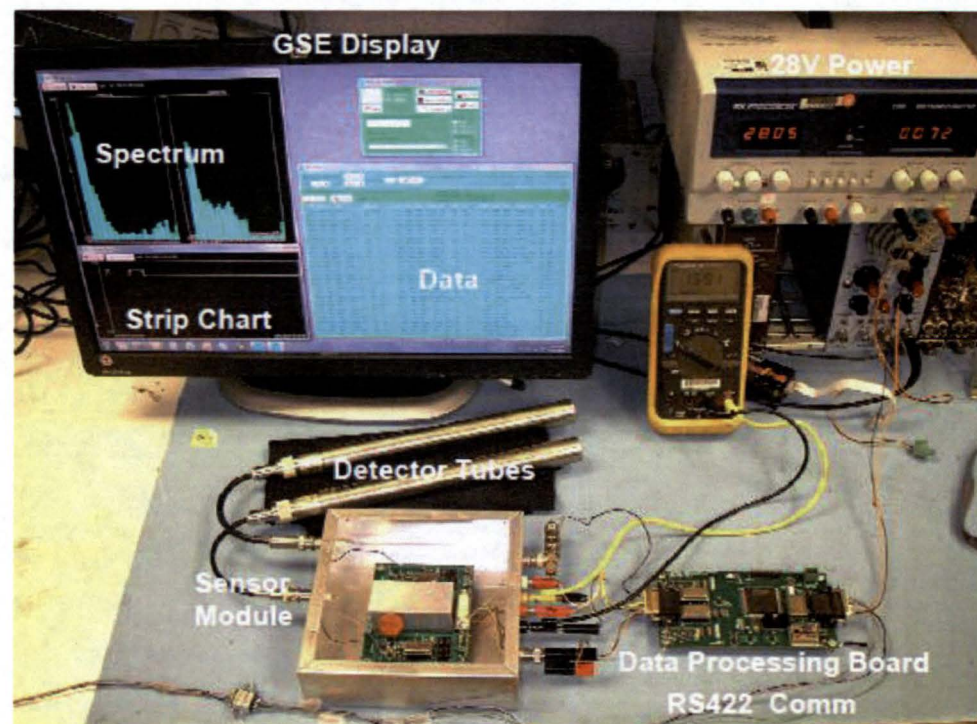



Near Infrared Volatile Spectrometer Subsystem (NIRVSS) Functions & Design Constraints






NIRVSS Functions

- Quantify amount of water in lunar regolith at a minimum abundance of 0.5% by mass
- Identify surface bound H₂O/OH
 - Map at rover speed $\leq 10 \pm 1$ cm/s
- Bound understanding of mineral content in regolith
- Identify volatiles, including water content and form evolved during auger/drilling
- Bound volatile presence in top 20-30 cm of regolith during auger/drilling
- Enable observation under all lighting conditions
- Image drill area with sufficient Field Of View to observe 22 cm of tailings with resolution at $\sim 200 \mu\text{m}$ scale
- Operate ≥ 6 hrs in permanently shadowed area

NIRVSS Design Constraints

- Identifying volatile and mineralogical features in the near-infrared spectrum in the range of 1.8-3.2 microns with a spectral resolution of less than $0.05 \mu\text{m}$.
- Mount Near IR, Camera, and lamp to view auger/drill area
- Achieve SNR ≥ 100 at $3 \mu\text{m}$ while drilling
- Operate $+0$ to $+45^\circ\text{C}$
- Mass: 7.7 kg
- Power: 16.31 W ave (NIR, camera, & lamp)

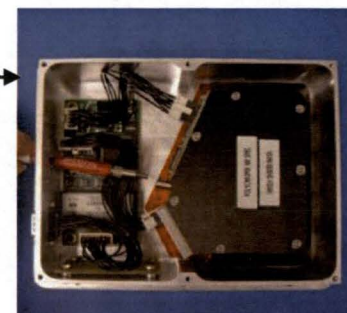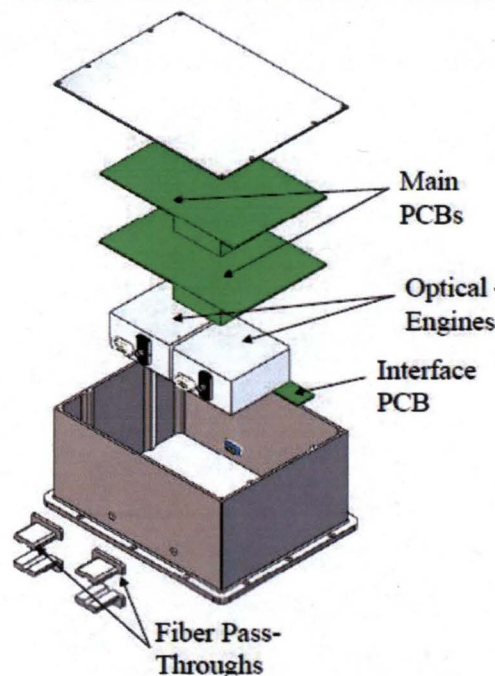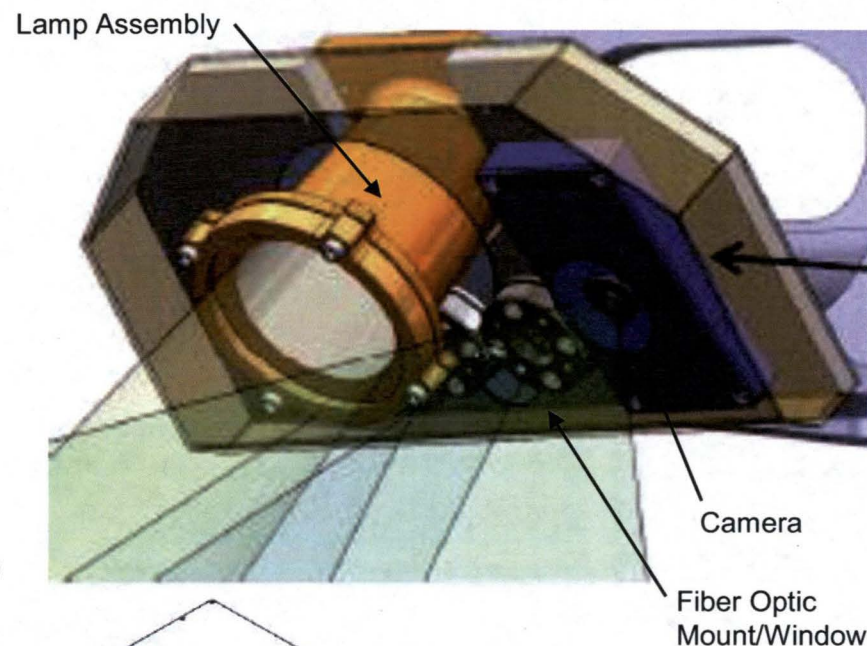

Spectrometer Assembly



Sample Acquisition Subsystem: DESTIN Functions & Design Constraints




DESTIN Functions

- Penetrate substrate 1m to collect, retrieve, transfer samples from >0.75 m
 - Unconsolidated regolith to consolidated
- Auger material to surface from depth up to 0.5 m
- Maintain sample phase, chemical state, and stratigraphy
- Sample Transfer Receptacle (STR)
 - Section samples to 12.5 cm length x 1.6 cm dia.
 - Deliver $\geq 90\%$ of sample to OVEN
 - Deliver sample to OVEN at <150 K
- Minimize sample cross-contamination
- Abandon drill rod if stuck and recover
- Autonomous operation in shadowed region for ≥ 6 hrs
- Measure temperature near core
- Measure: sample hardness, energy required for penetration, rate of cut, drill depth, instantaneous drilling power, weight on bit, torque, rpm

DESTIN Design Constraints

- Dimensions of all components $<1.35\text{m} \times 0.75\text{m} \times 1\text{m}$
- Safe position must not interfere with rover locomotion
- Operate when tilted up to 15 degrees
- No consumables
- 50 auger/drill operations
- Mass: <40 kg max.; 25 kg goal
- Power: <150 W ave.
- Static Force: 120 N max.

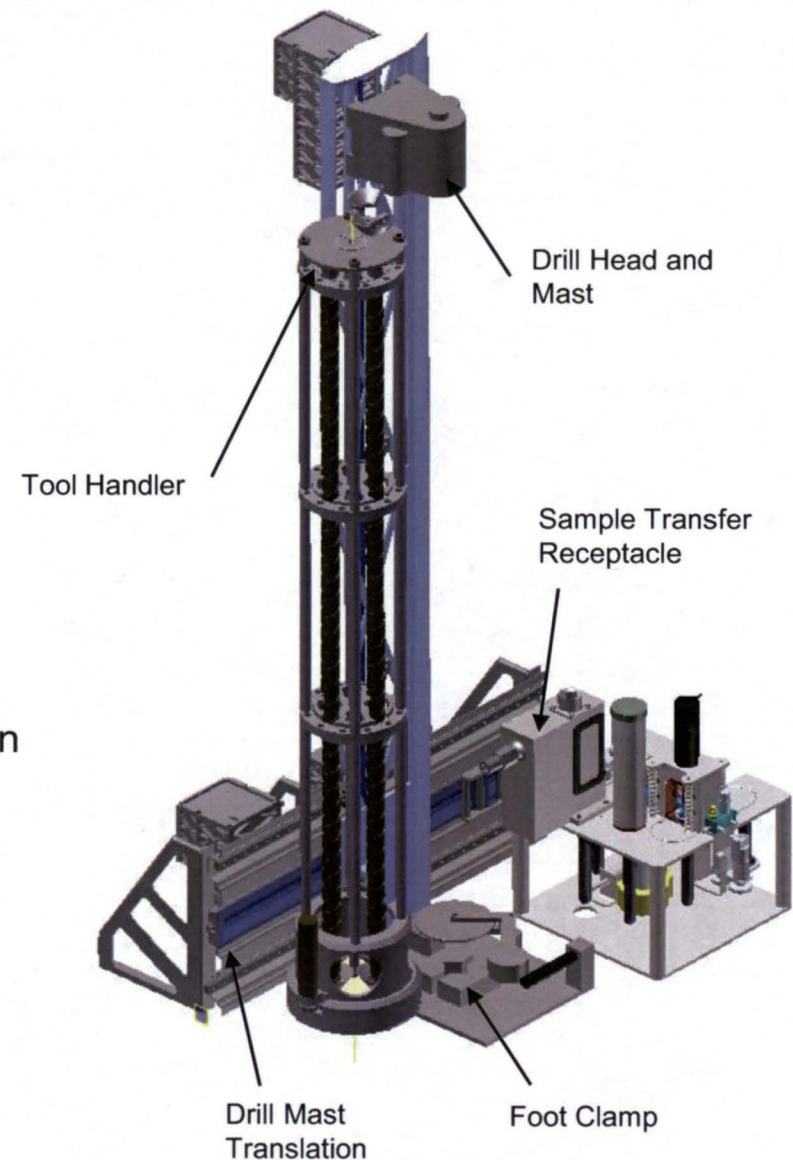



Oxygen and Volatile Extraction Node (OVEN) Functions & Design Constraints




OVEN Functions

- Heat multiple samples to 423 K (150 °C) for volatile extraction
- Heat a minimum of one sample to 1173 K (900 °C) in presence of hydrogen for oxygen extraction from regolith
- Accept sample from Sample Transfer Receptacle (STR)
 - May be solid core or granular material:
16 mm dia. by 125 mm ; 25 to 60 gms
- Accept samples at <150 K
- Maintain samples at <175 K prior to sealing.
- Measure sample mass before and after processing to +/- 0.1 gm accuracy
- Dump sample after processing to remove a minimum of 95% of sample material
- Selectively accept and reject (dump) sample without processing
- Transfer volatiles released (at 150 and 900 °C sample temperatures) to surge tank in analysis and water droplet demo.

OVEN Design Constraints

- Minimum leak rate during sample processing - psi over 4 hrs at 100 psi (at 150 and 900 °C)
- Heat 40 samples in 5-7 day mission
- Mass: 10 kg
- Power: <300 W ave.
- Minimize height of OVEN subsystem since STR drill tube must extend above OVEN

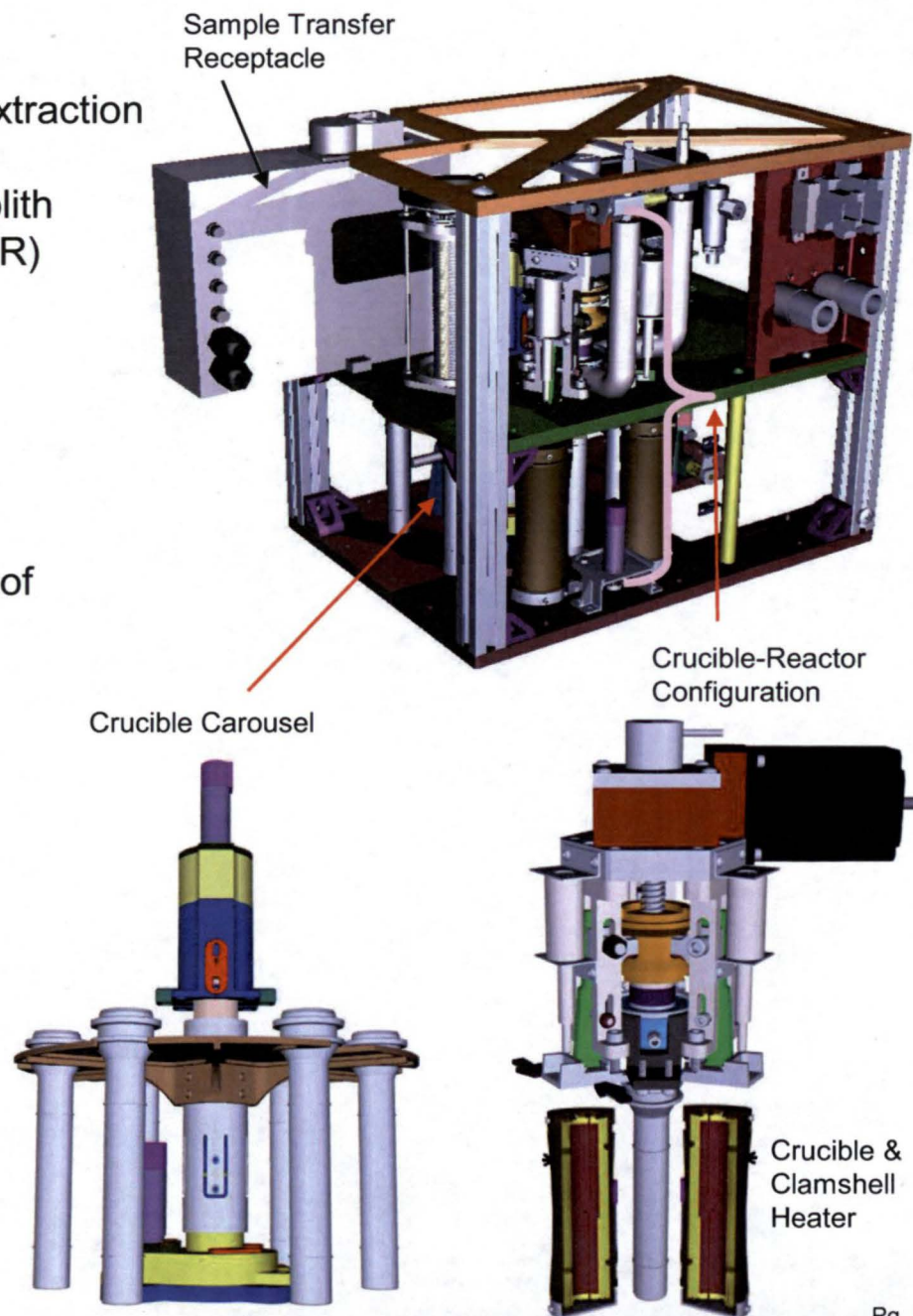





Lunar Advanced Volatile Analysis (LAVA) Functions & Design Constraints



LAVA Functions


- Accept gas samples from OVEN at up to 100 psi and 423 K (150 °C) to 1173 K (900 °C)
- Identifying and measuring the relative abundance of the volatile constituents of the lunar regolith.
 - Quantify amount of evolved water in from the lunar regolith at a vapor concentration of 0.1% to 95% with a standard deviation of 5% relative standard deviation or the absolute standard deviation of 0.2% water , whichever is greater
 - Measure volatile constituents of the lunar regolith to 70 amu including, CO, H₂O (g), H₂, [H₂S, NH₃, SO₂, C₂H₄. are Goals]
 - Measure the D/H and O^{16/18} isotope ratios (Goal)
- Collect and provide images of water collected through volatile extraction and hydrogen reduction of regolith
- Provide and regulate hydrogen gas to OVEN for Regolith Oxygen Extraction

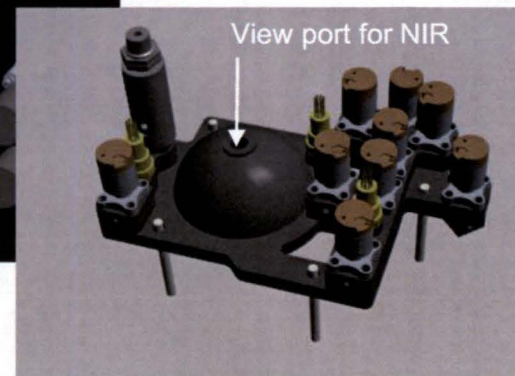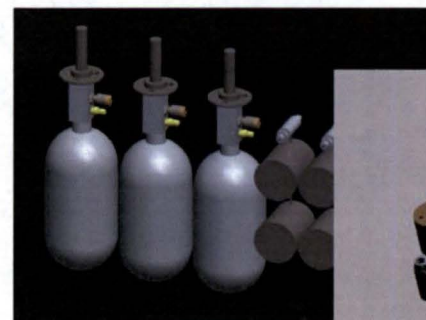

Surge Tank and Fluid Subsystem



Water Droplet Demo w/ Camera

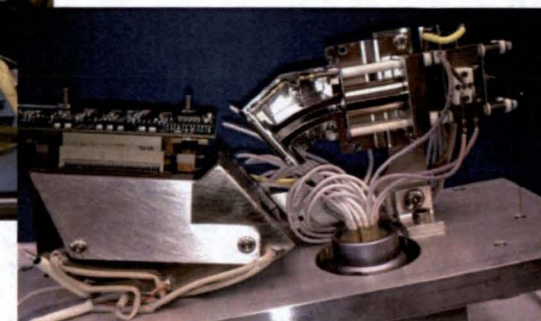

Gas Chromatograph/Mass Spectrometer

LAVA Design Constraints

- Complete GC-MS analysis in under 2 minutes
- Mass: 15 kg
- Power: <100 W ave.



RESOLVE Gen III



Purpose: Develop a flight-like unit that can fit on a rover and operate in the lunar environment

Sample Acquisition System

Auger/Core Drill Subsystem [CSA]

- Collect and transfer subsurface material down to 1 m below surface
- Maintain sample stratigraphy and volatiles (below 150 K)
- Meter samples for processing
- Auger material to surface for evaluation
- Measure geotechnical properties of regolith during drilling

Surface Mineral/Volatile Evaluation

Near Infrared Volatile Spectrometer Subsystem (NIRVSS) - ARC

- Measure surface bound OH/H₂O while traversing (at min. of 0.5% by mass)
- Detect form of water (ice/hydration) in auger tailings
- Detect water vapor in evolved gases
- Image surface and drill tailings

Resource Localization

Neutron Spectrometer Subsystem (NSS) -ARC

- Locate hydrogen and hydrogen bearing volatiles down to 1 meter below the surface while traversing (at min. of 0.5% by mass)

RESOLVE Mission Requirements

- Nom. Mission Life = 5+ Cores; 14 Days
- Mass = 170 kg rover/80 kg payload
- Ave. Power; 200-300 W

Volatile Content/Oxygen Extraction

Oxygen & Volatile Extraction Node (OVEN) - JSC

- Accept samples from Sample Acquisition System
- Heat samples from <150 K to 423K for volatile extraction
- Heat samples to 1173 K for oxygen extraction
- Transfer evolved gases to LAVA volatile analyzer

Volatile Content Evaluation

Lunar Advanced Volatile Analysis (LAVA) - KSC


- Accept evolved gas from OVEN; provide hydrogen for oxygen extraction
- Perform analysis in under 2 minutes
- Measure water content in evolved gas
- Characterize volatiles of interest (below 70 amu)
- Measure D/H and O^{16/18} isotopes
- Capture & image water evolved

Operation Control Flight Avionics - KSC

- Space-rated microprocessor
- Control subsystems and manage data

Surface Mobility [CSA]

- Traverse wide range of lunar surface/material conditions
- Tele-operation and autonomous traverse modes
- Carry RESOLVE payload; provide power, comm., and thermal management

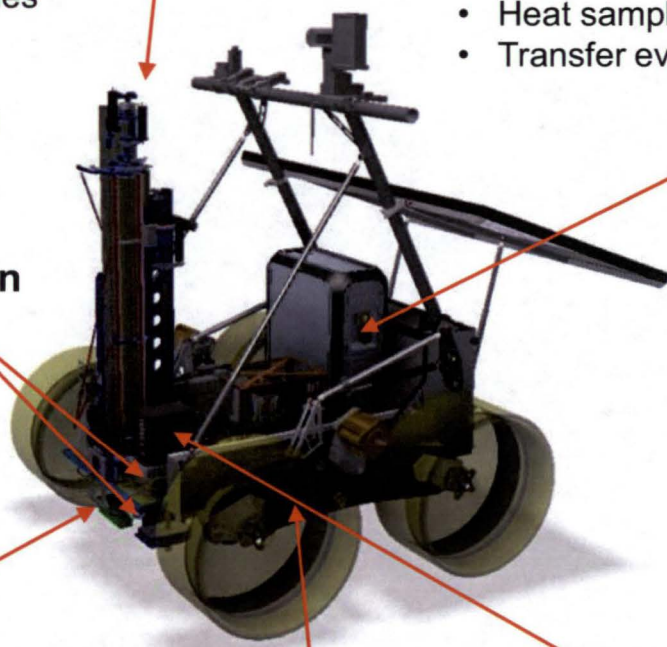



Vacuum Demonstration Unit GCMS Design: Integrated COTS Approach



Changes required for Vacuum Unit –Flight forward

-
- Maintain fast scan rate and analyte sensitivity
- Less mass
- Less power
- More rugged (Vibe tolerant)
- Dilution capability (prior unit saturated at about 5% water)
- Software control from Xiphos- Payload control
- Integrated Avionics controlled by payload
- Thermal Vacuum (materials and electronics)



Vacuum Chamber – Bell Jar



Inner Diameter - 13"


Total Height - 27"

Varian Turbo-V80 Turbo Vacuum Pump

Current configuration:



7 Total Ports – 5 Available for feed-throughs


- Pressure Transducer & Ion Gauge
- 20-pin Electrical Port
- Double 1/4" Swagelok Port
- Single 1/8" Swagelok Port

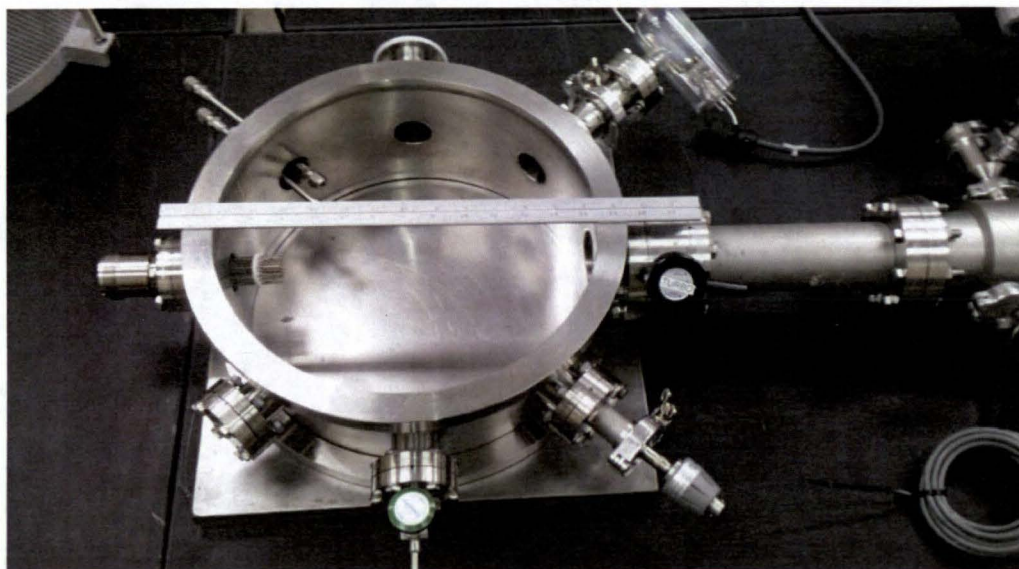

RESOLVE Sample Delivery System (SDS)


- Temperature and pressure-controlled sample chamber that provides a consistent sample environment for the GC or MS
- Accurately manages temperature (0-150°C) and sample pressure (0-35 psia), with the goal of reducing sample carry over, particularly water
- Can control sample gas and dilution gas delivery to target desired concentration ranges

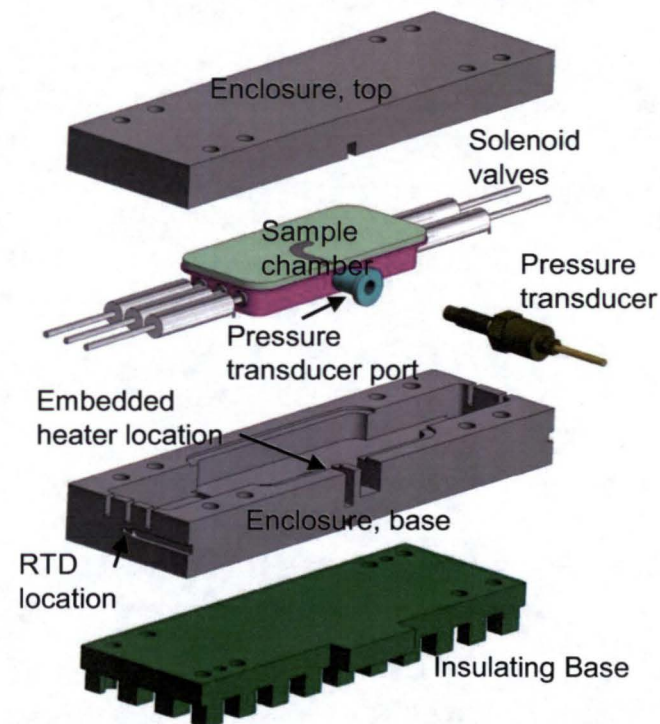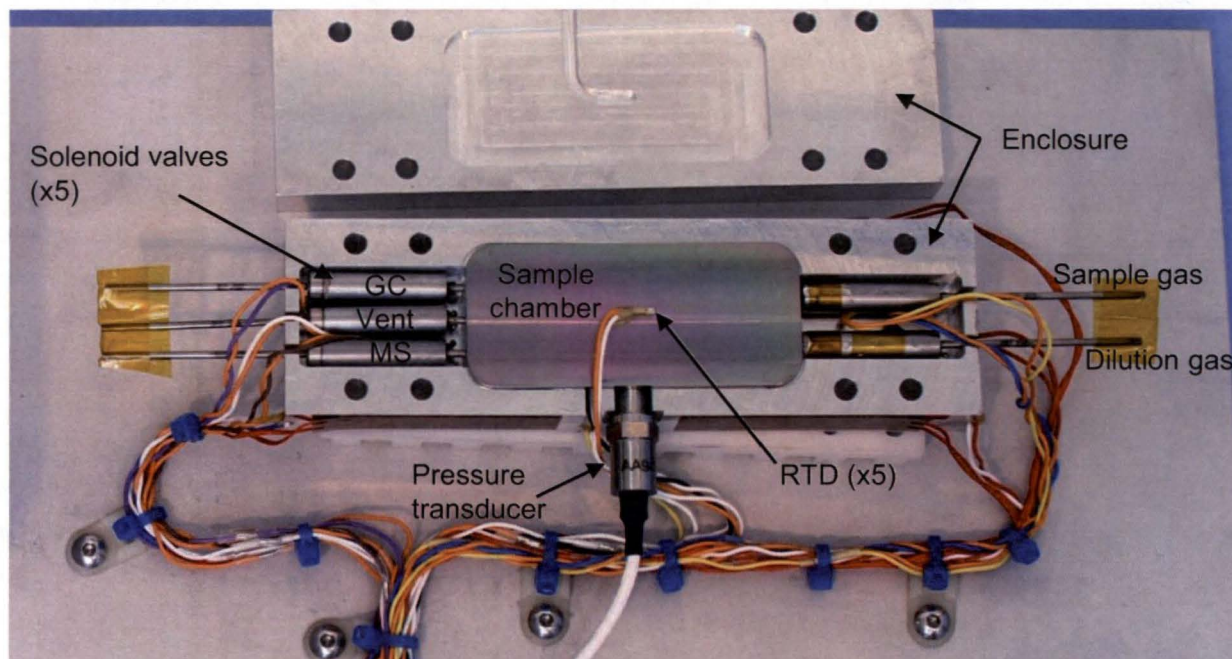

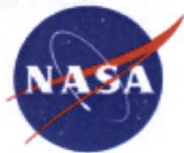

SDS Performance

10psig Sample, 480 ppm Water in Argon

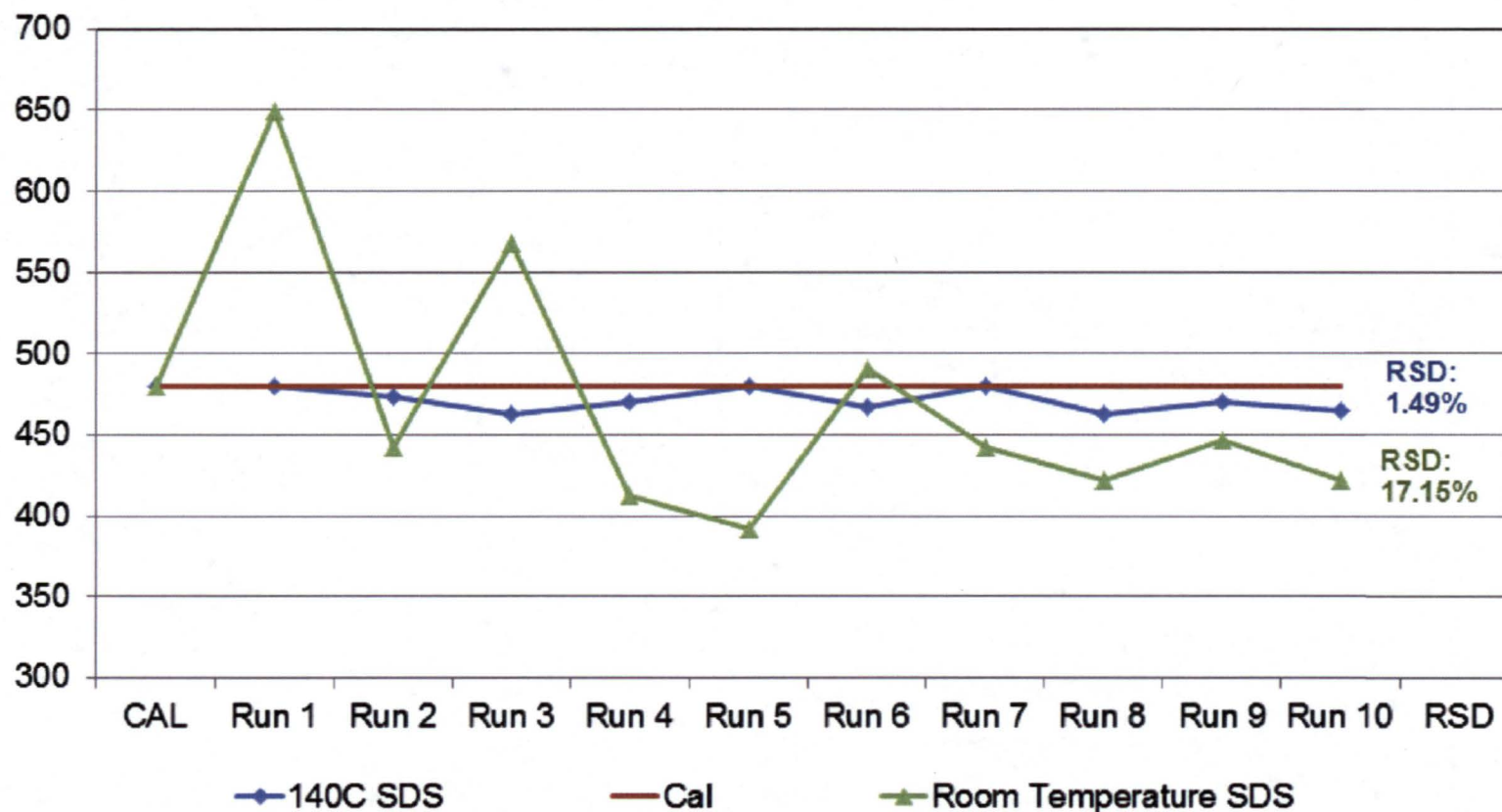






- SDS conditioning proved to be key for stable water quantification using GC
- High temperature conditioning sample with SDS (e.g. 140°C) provided much better water signal stability and repeatability



Key Performance Parameters (KPP)



Performance Parameter	State of the Art	Threshold Values	R&TD Goals	Actual Values / Current Best Estimate
System Mass	15 kg	15 kg 	11 kg	14.94 kg ETU (14.69 kg for flight) (MEL input, 5/1/13)
Average Power ¹	80 W	100 W	80 W 	75 W checkout (PEL, 2/28/13)
Peak Power	163 W	200 W 	160 W	200 W (PEL, 2/28/13)
Water Vapor Concentration ²	N/A ³	0.5-95%	0.1-99% 	0.1-99% (Test Data 6/3/2013)
Mass range (MS systems)	Ion trap 12-150amu Mag Sector 2-130amu Quad 1-60amu (1.8sec/mass scan) All Scanning	 Demonstrated data collection of a full mass spectrum at a sample rate of ≥ 6 Hz for 1-65 AMU	Demonstrated data collection of a full mass spectrum at a sample rate of ≥ 6 Hz for 1-80 AMU	6.7 Hz, 1-70 AMU (Test data, 5/3/13)


¹ State of the art values are based on laboratory instruments, not flight instruments

² Gas phase concentration, does not refer to percent mass in soil

³ Not stated for existing flight instruments, depends on background levels and carry over

GC Evolution (Inficon)

COTS (1 st Gen)	2nd Iteration
Cold spots (TCD and other Connectors)	Thermal control of heated zones
Operated by COTS software	Operated by Q6
Column oven inefficient and heavy	Decreased mass and power usage
Portable design	More rugged design

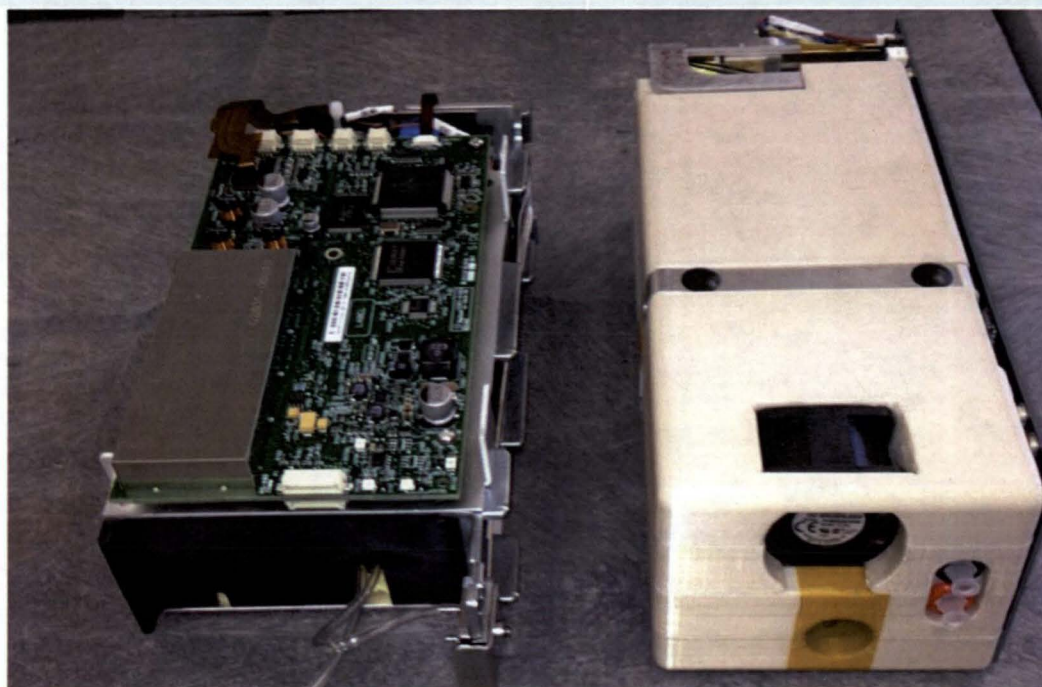

First Generation GC

Second Generation GC




MS Evolution (OIA)



COTS	1 st Iteration	2 nd Iteration
Integrated vacuum chamber	No vacuum chamber	No vacuum chamber
1.1 Kg magnet	1.1 Kg magnet	0.6 Kg magnet
Two detector boards	Two detector boards	Merged detector boards
Mounted to wall of vacuum chamber (3.9 Kg)	Mounted to lighter material base (2.0 Kg)	Mass further reduced (1.5 Kg)
No dust shield	Full/partial dust shield	Full dust shield
No internal thermal isolation	Internal thermal isolation	Internal thermal isolation
Standard electronics	Standard electronics	Some modified electronics
Electronics in ambient air	Electronics in ambient air	Electronics in vacuum

MS Design Iteration 1

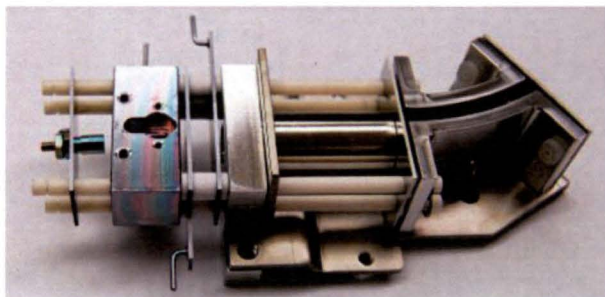


- COTS Instrument
 - HV tester (hardware)
 - IonCam w/ Silconert

- Fabrication
 - Bell Jar MS


- Software
 - Three ICDs

- Design Trades (to be detailed later)



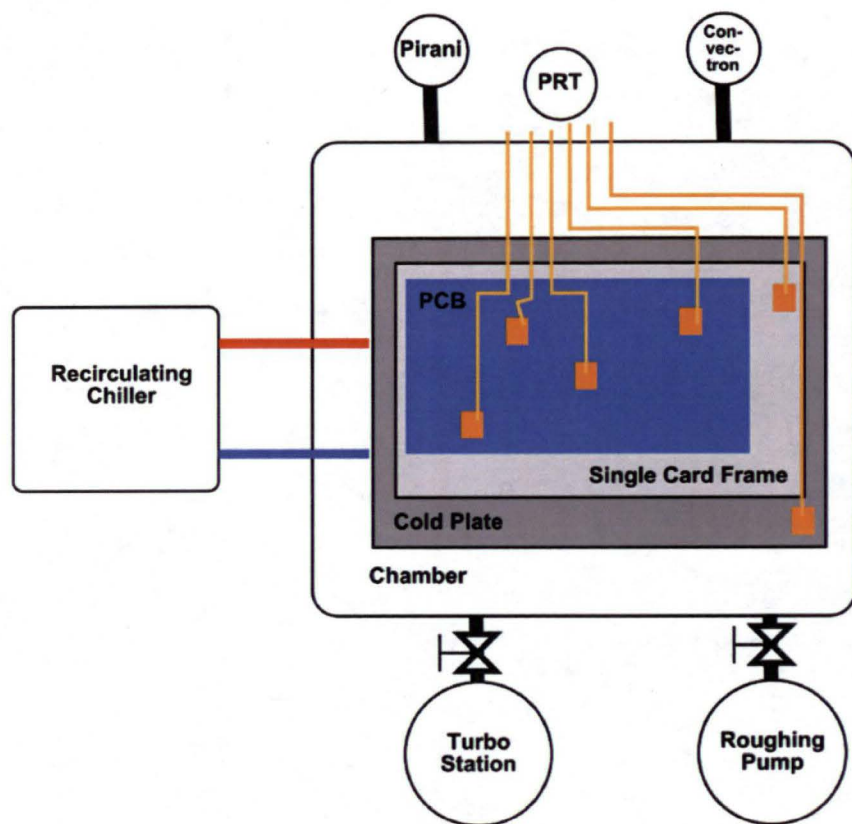

Future testing for other COTS boards will utilize PRTs in this manner.

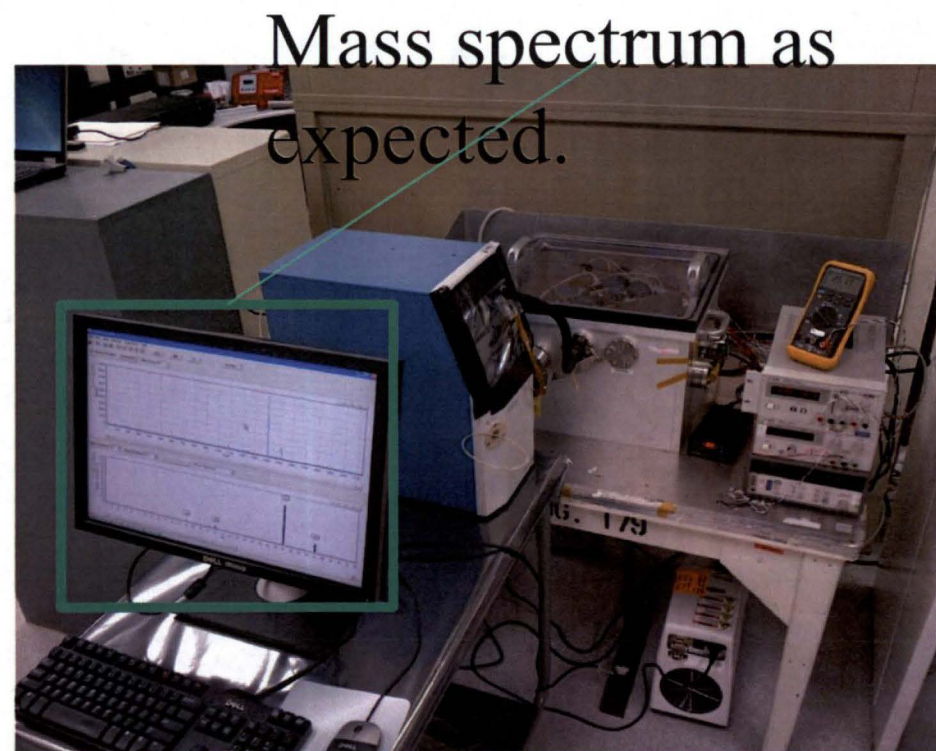

IonCam operating with the IHV board in thermal vacuum.



Thermal Imaging



- Thermal imaging was used on the board operating in ambient air to determine location of heat sources.

Majority of
heat from
Blackfin
processor



IHV Board in Mounting

DIG Board

Cover was removed and temperatures were at steady-state: thermocouples were not attached. Camera setting was auto range. Vacuum-side was on.




MS Electrical: Thermal Management



Each electronics card is mounted to a chassis

Current status: Card chassis are in fabrication. Final design of mounting frame in progress (may change from what is shown). IR Camera work in progress.

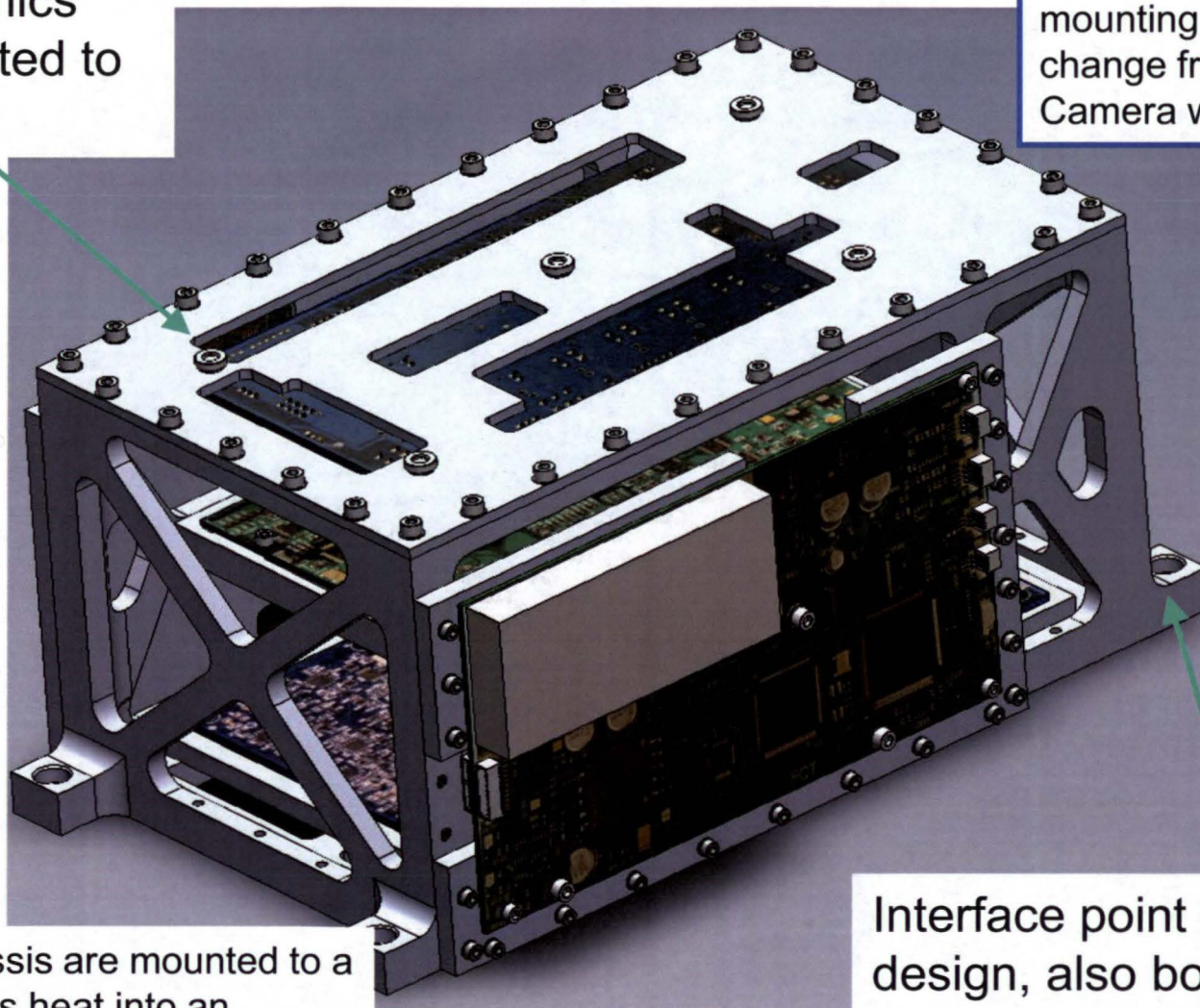


The various chassis are mounted to a frame which sinks heat into an interface with RESOLVE.

Interface point (4 in current design, also bottom surface of frame).




SDS/GC/MS Block Diagram

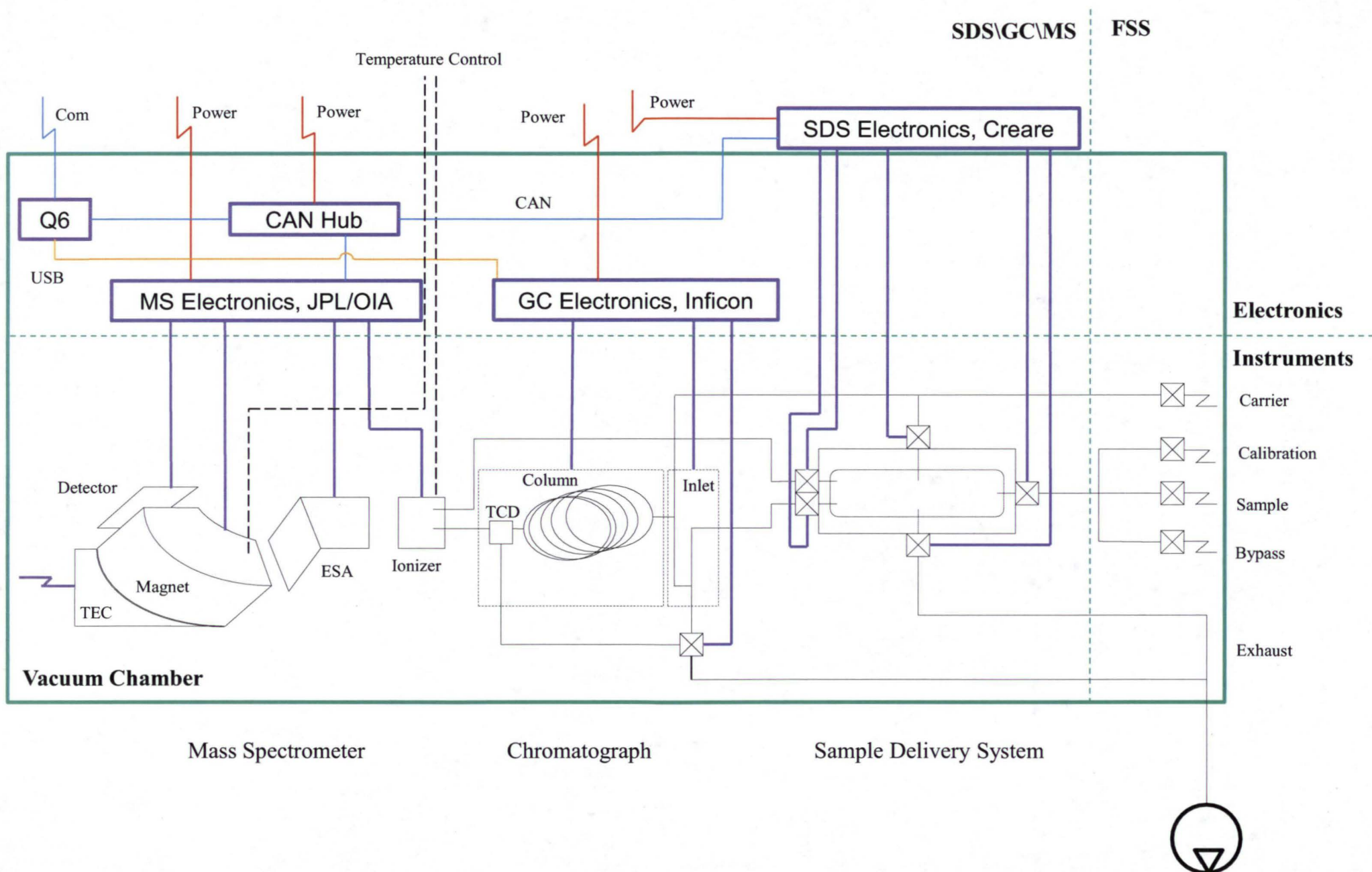



LAVA Software Overview – ETU Architecture

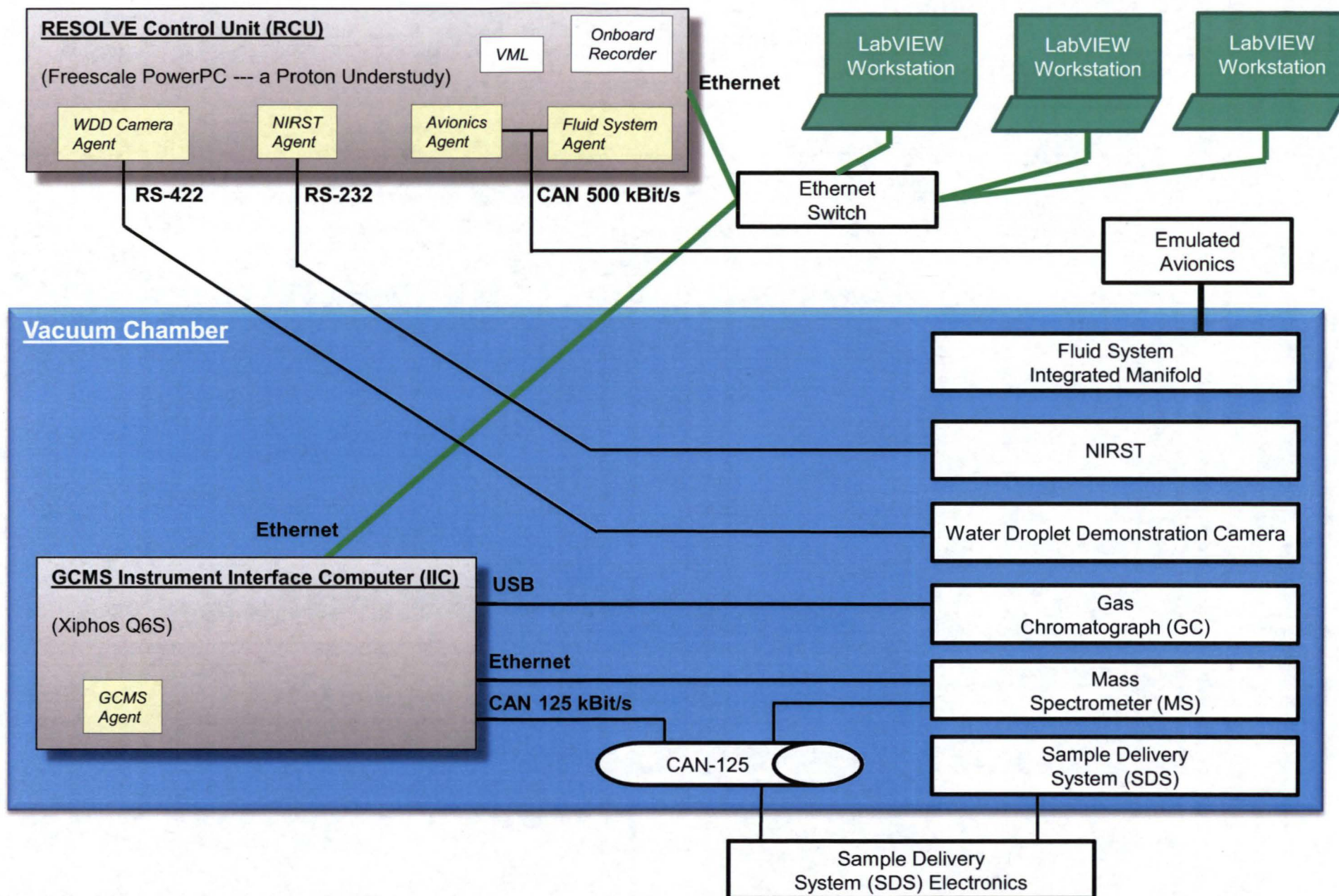



Key Performance Parameters (KPP)



Performance Parameter	State of the Art	Threshold Values	R&TD Goals	Actual Values / Current Best Estimate
System Mass	15 kg	15 kg	11 kg	14.94 kg ETU (14.69 kg for flight) (MEL input, 5/1/13)
Average Power ¹	80 W	100 W	80 W	75 W checkout (PEL, 2/28/13)
Peak Power	163 W	200 W	160 W	200 W (PEL, 2/28/13)
Water Vapor Concentration ²	N/A ³	0.5-95%	0.1-99%	0.1-99% (Test Data 6/3/2013)
Mass range (MS systems)	Ion trap 12-150amu Mag Sector 2-130amu Quad 1-60amu (1.8sec/mass scan) All Scanning	 Demonstrated data collection of a full mass spectrum at a sample rate of ≥ 6 Hz for 1-65 AMU	Demonstrated data collection of a full mass spectrum at a sample rate of ≥ 6 Hz for 1-80 AMU	6.7 Hz, 1-70 AMU (Test data, 5/3/13)


¹ State of the art values are based on laboratory instruments, not flight instruments

² Gas phase concentration, does not refer to percent mass in soil

³ Not stated for existing flight instruments, depends on background levels and carry over



Vacuum Demonstration Unit GCMS Design: Integrated COTS Approach

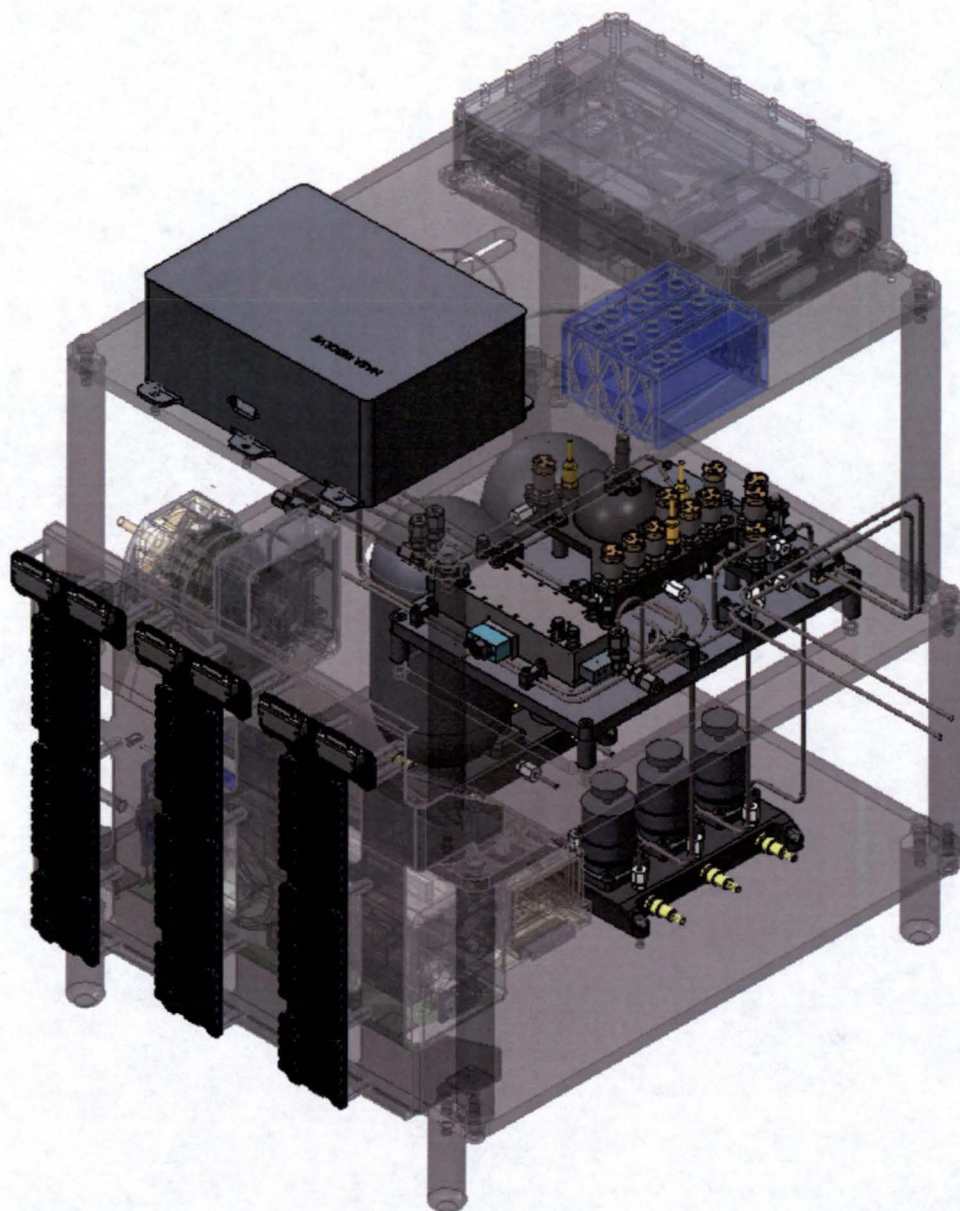



RESOLVE AES/OCT Project Status



Gen IIIA: Field Development Unit (FDU)


▪ FDU System Requirements Review	03/03/11	Completed
▪ FDU 30% Design Review	05/25/11	Completed
▪ FDU 90% Design Review	08/26/11	Completed
▪ FDU 90% Delta Design Review	10/28/11	Completed
▪ Field Demo Subsystem HW Initial Delivery to KSC	02/27/12	Completed
▪ Field Demo HW Integration onto Rover Complete	06/29/12	Completed
▪ Field Demo HW Delivered to Field Test Location	07/09/12	Completed
▪ Demonstrate Integrated RESOLVE ops on Rover in Field Test	07/27/12	Completed
▪ <i>AES Project Continuation Review</i>	09/18/12	Completed

Gen IIIB: Engineering Test Unit (ETU)

▪ ETU SRD Initial Delivery	12/16/11	Completed
▪ Complete ETU System Requirements Review	08/29/12	Completed
▪ ETU 30% Design Review	12/14/12	Completed
▪ ETU 90% Design Review	07/26/13	Completed
▪ <i>AES Project Continuation Review</i>	09/13	
▪ <i>OCT Project Evaluation/Continuation Review</i>	09/13	
▪ ETU Subsystem Environment Testing Complete	05/12/14	




LAVA Software Overview – Flight Architecture

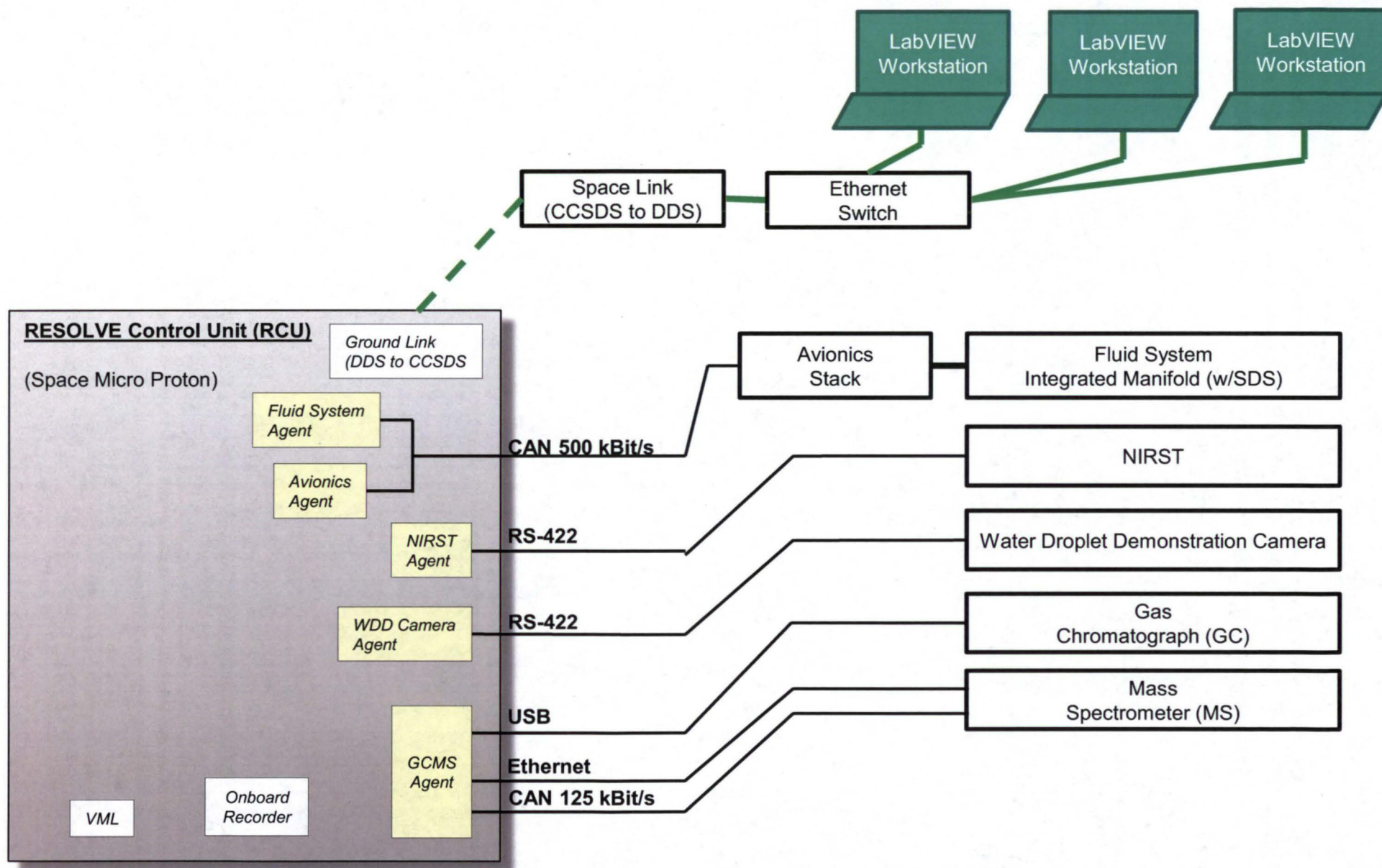



Brassboard GC/MS for RESOLVE

Ken Wright (Inficon); Andres Diaz (NTCR); Bob Kline-Schoder, Paul Sorensen, Brandon Smith (Creare)

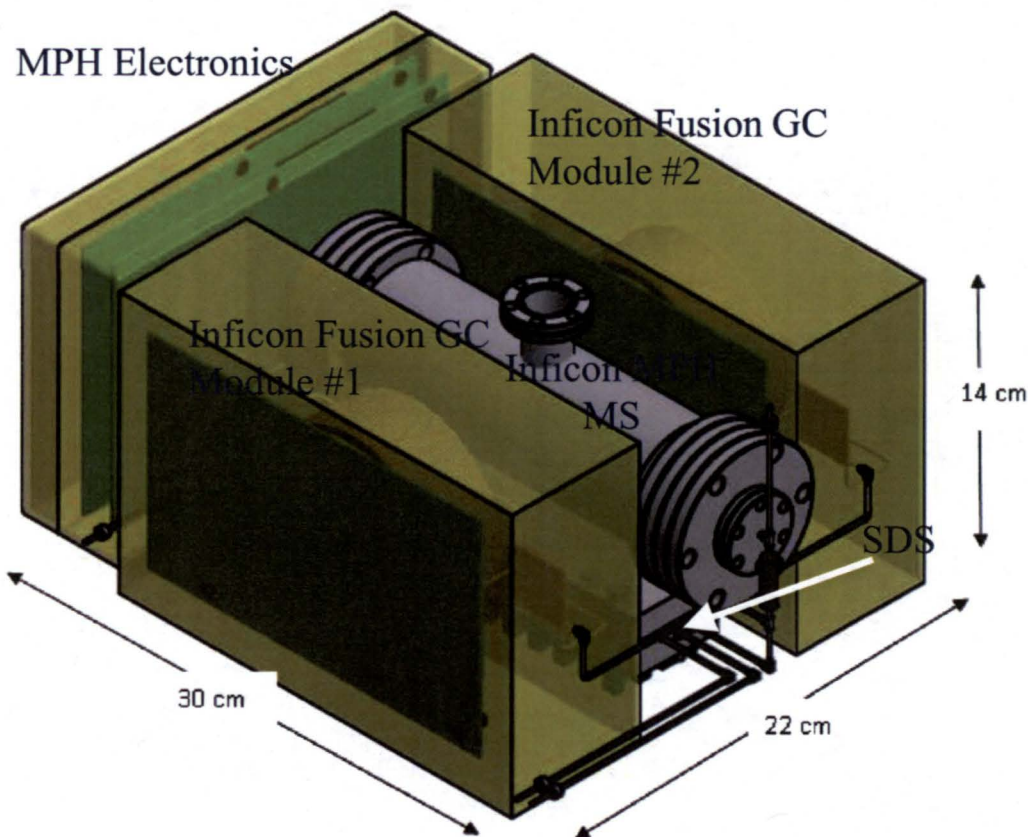

INFICON's Transceptor MPH:

- Ultra Fast Measurements (0-65 m/z @ 8.5Hz).
- High Sensitivity and Low Noise = Lower Minimum Detectable Partial Pressure ($<1\text{E-}15$ torr).
- Smaller and Lighter ($>30\%$ weight reduction)
- Nine Decades of Dynamic Range.
- New High Performance and Field Replaceable EM/FC Detector (Developed with Detector Technology Inc.).
- Field Replaceable Dual Filament Assemblies (W or Y2O3/Ir).
- Easy to Use Programming Interface - JSON over HTTP or LabVIEW.



Continuing Challenges

Budgetary Cuts
Requirements evolving
Commercial Partners Deliveries/IP
issues



Flight Forward Plans




- Stand down on purchases until payload matures
- Ground Calibration
- Failure Mode Testing
- Material off-gassing testing
- Vibe tests
- Verification and Validation
- Compare results from SBIR Phase II



BACKUP

Complete RESOLVE Mission Traverse

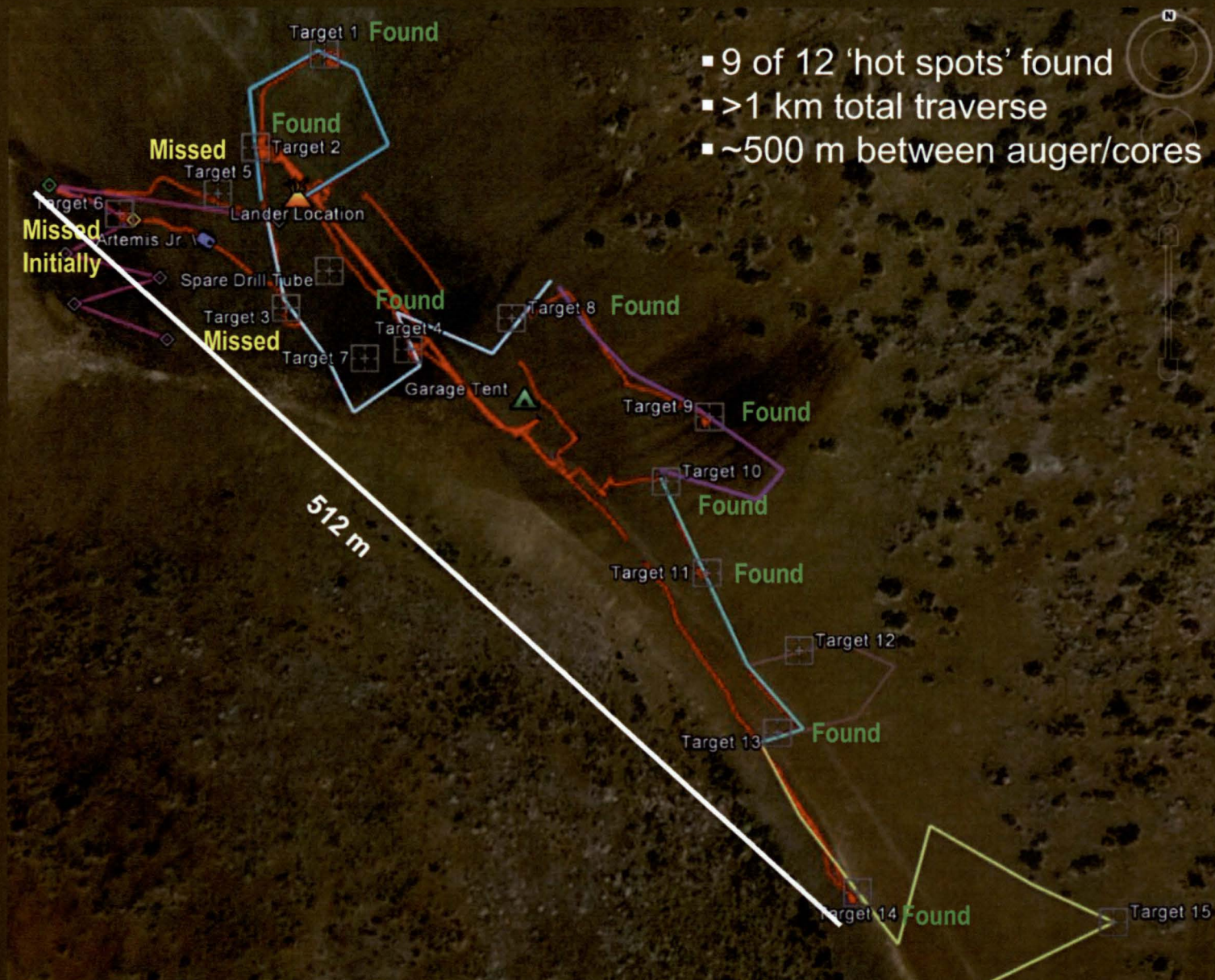



Key RESOLVE Mission Design Trades



Mission Attributes	Base	Mid	Full
Location	Long duration sunlight	Min. Sun/Shadowed	Permanent Shadow
Sample Site Selection	Surface features/minerals	Neutron Spec on Rover	Neutron Spec with GPR
Subsurface Sample Acquisition	Arm/scoop	<i>Auger w sample transfer</i>	Core Drill/Push Tube w sample transfer
Sample of Interest	Rock/regolith	Ice	Polar volatiles
Sample Depth	<i><0.75 m</i>	1.0 m	2.0 m
Sample Measurement	Downhole Optical for ice	Oven w Tunable Diode Lasers	Oven with GC/MS and Near IR
Sample Preparation	None	Crushing	Thin Section
Mineral Characterization	None	Single instrument - Near IR	Multiple Instruments
Regolith/Dust Physical Characterization	None	Camera & Drill Response	Microscope & Geotechnical Instruments
Volatile/Product Collection	None	Water	Water and gas volatiles
Oxygen Extraction from Regolith	None	H₂ Reduction w Same Oven	Separate demo
Temperature/Radiative Environment Characterization	None	External temp sensor	Instrumented Radiator
Mobility	None - Lander	Hopper	Rover
Power	<i>Non-recharge battery</i>	Battery/Solar Array	Nuclear
Communications	Direct to Earth-rover	<i>Direct to Earth-lander; rover relay</i>	Comm Relay Satellite

Blue Bold = Baseline

Red Italics = Backup



Modified High Voltage Board




- General requirements
 - Retain current board dimensions and part location of larger parts
 - Boards should be manufactured under J-STD-001 E
 - Follow NASA recommendations for vacuum parts as much as practical
 - Do not install test points
- Modifications
 - Attempt to protect low voltage DAC (U33)
 - The DAC is usually destroyed when HV arcing occurs
 - Improve emission current measurement and control
 - Designing a new power control
 - Design of a new emission current control
 - New firmware to support changes
 - See “docushare” for design documentation (BOM, schematics) and control loop simulation



IHV Board - Filament Control Function



For tuning and start up the electron emission is controlled by a power set-point. After the Electron Impact ionizer is tuned and stable, the control can be switched to emission control.

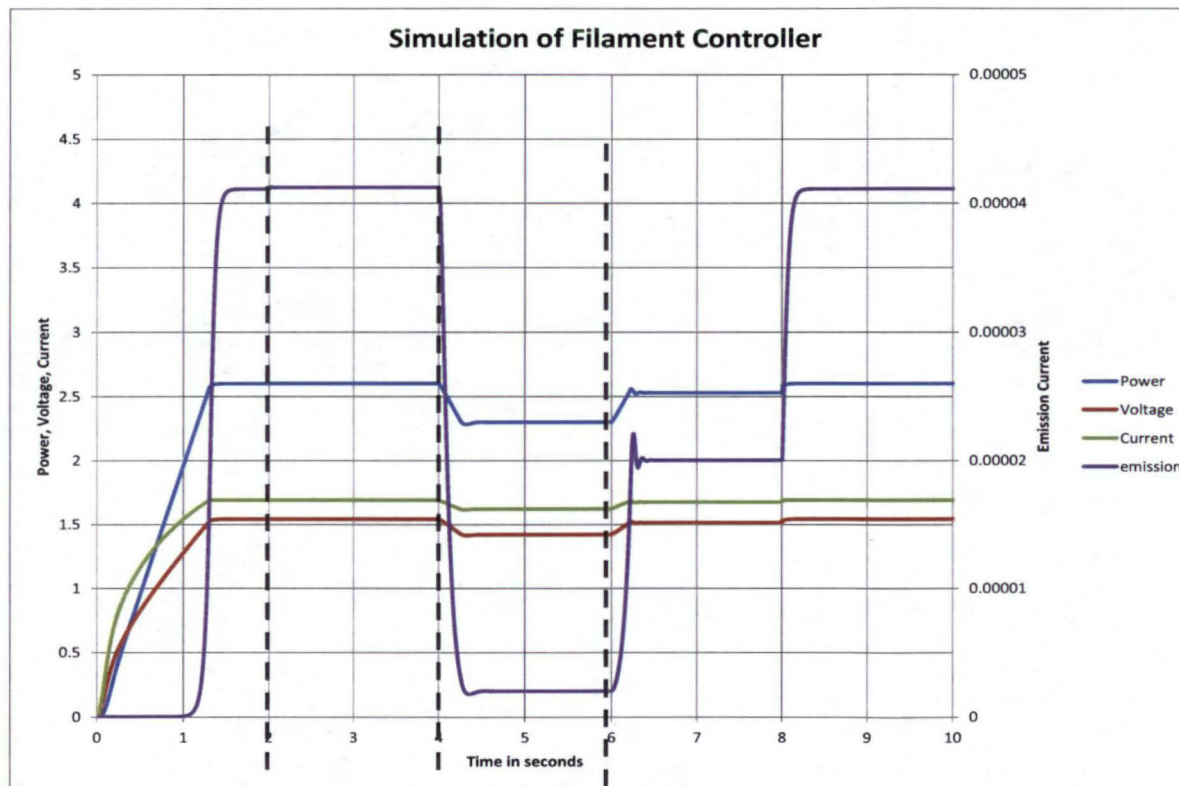

0 to 2 seconds:

Power mode set to 2.6 W

Software ramp 2 W/ s

2 to 4 seconds:

Switch to emission control (at 42 uA). No instabilities, small randomness bump added. **4 to 4 to 6**

seconds:

In emission control set to 2 uA.

6 to 10 seconds:

In emission control set to 20 uA

8 to 10 seconds:

Power mode set to 2.6 W



IHV Board – Timeline & Risks



August 5 - Deliver 1 modified IHV to JPL for vacuum testing

September 2 – Deliver 1 modified IHV to KSC as part of Phase II MS

- July 24 receive assembled board at OI
 - EMCO DC/DC has long lead time
- General risks
 - Manual routing of the IHV board invites errors
 - Scale factors and component values were estimated and might need adjustment
 - J-STD-001 E will not allow patches but will require spinning the board
- Off-set for emission measurement is not controlled




Merged Detector Board



- Redesign IonCCD detector
 - Merge Detector Board and Personality Board
 - Reduce surface area, size, and weight
 - Eliminate board to board connector, improve routing
 - Incorporate DCRs
 - Improve mounting and ruggedization of detector
 - See “docushare” for design documentation (BOM, schematics)
- Ceramic PCB substrate
 - Better temperature equalization across the IonCCD die
 - Reduces VOC outgassing (low coefficient of thermal expansion)
- Standards to be used
 - AMS2404 for Nickel plating
 - AMS2422 and MIL-G-45204 for Gold plating
 - J-STD-001, rev. E and IPC-A-610 for assembly



Merged Detector Board - Dimensions

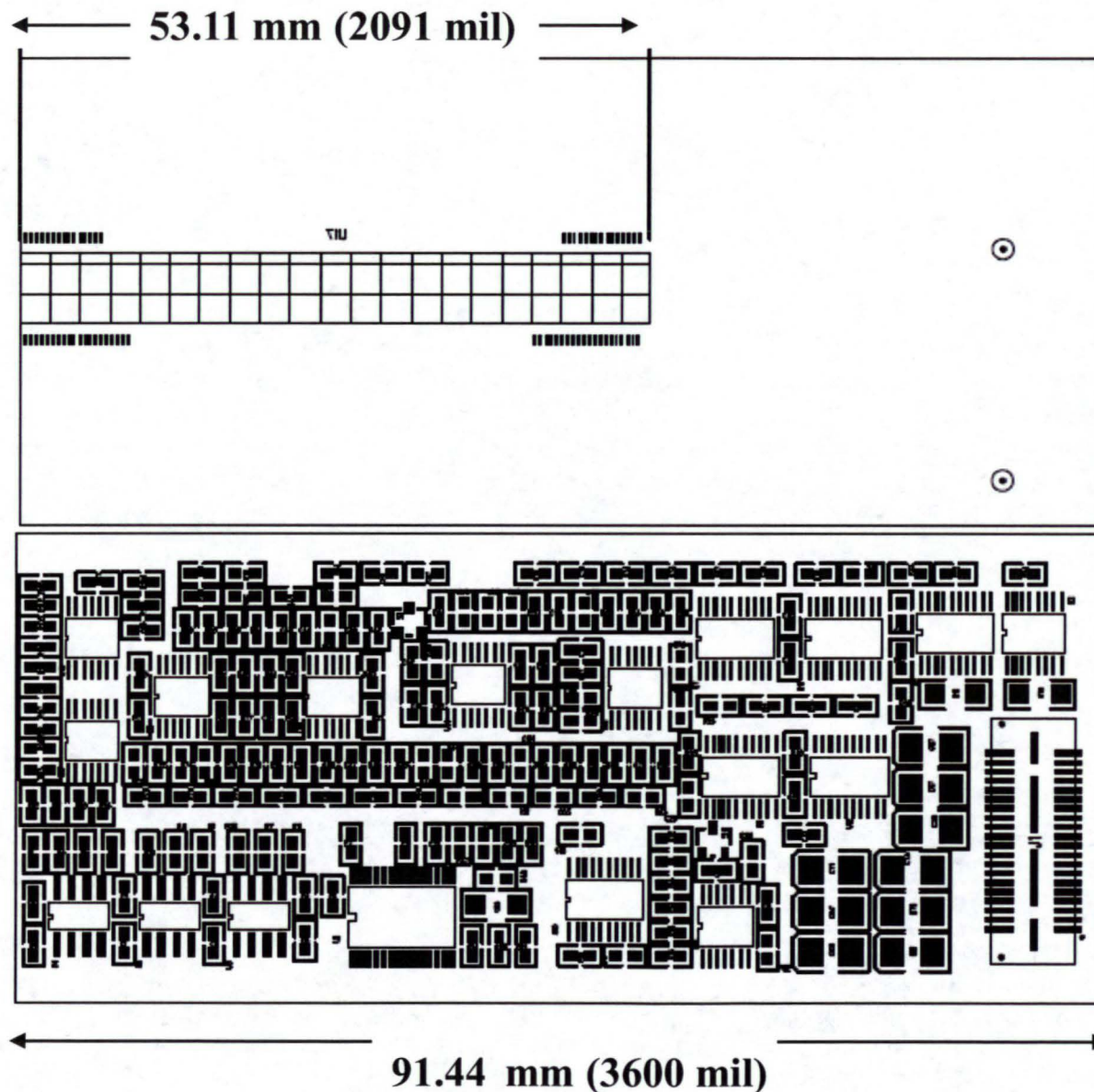

L: 91.44 mm

W: 40.64 mm

H: 2.50 mm

Weight: 50 g max

Front with IonCCD




Back with components



Mass Spectrometer




- Mass = 1.5 kg
- Overall Dimensions = 61 mm x 100 mm x 256 mm

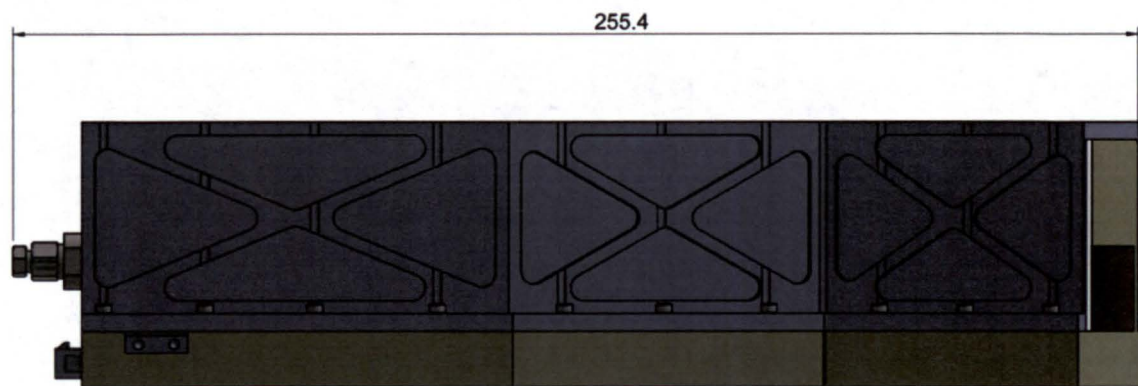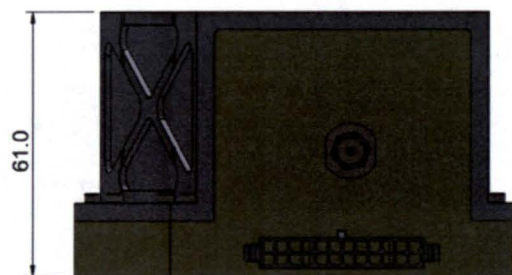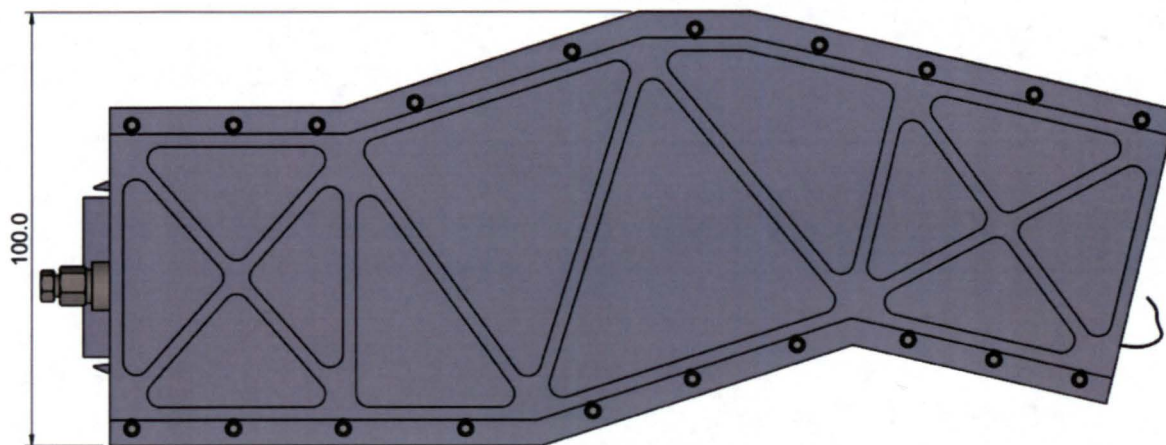




Mass Spectrometer

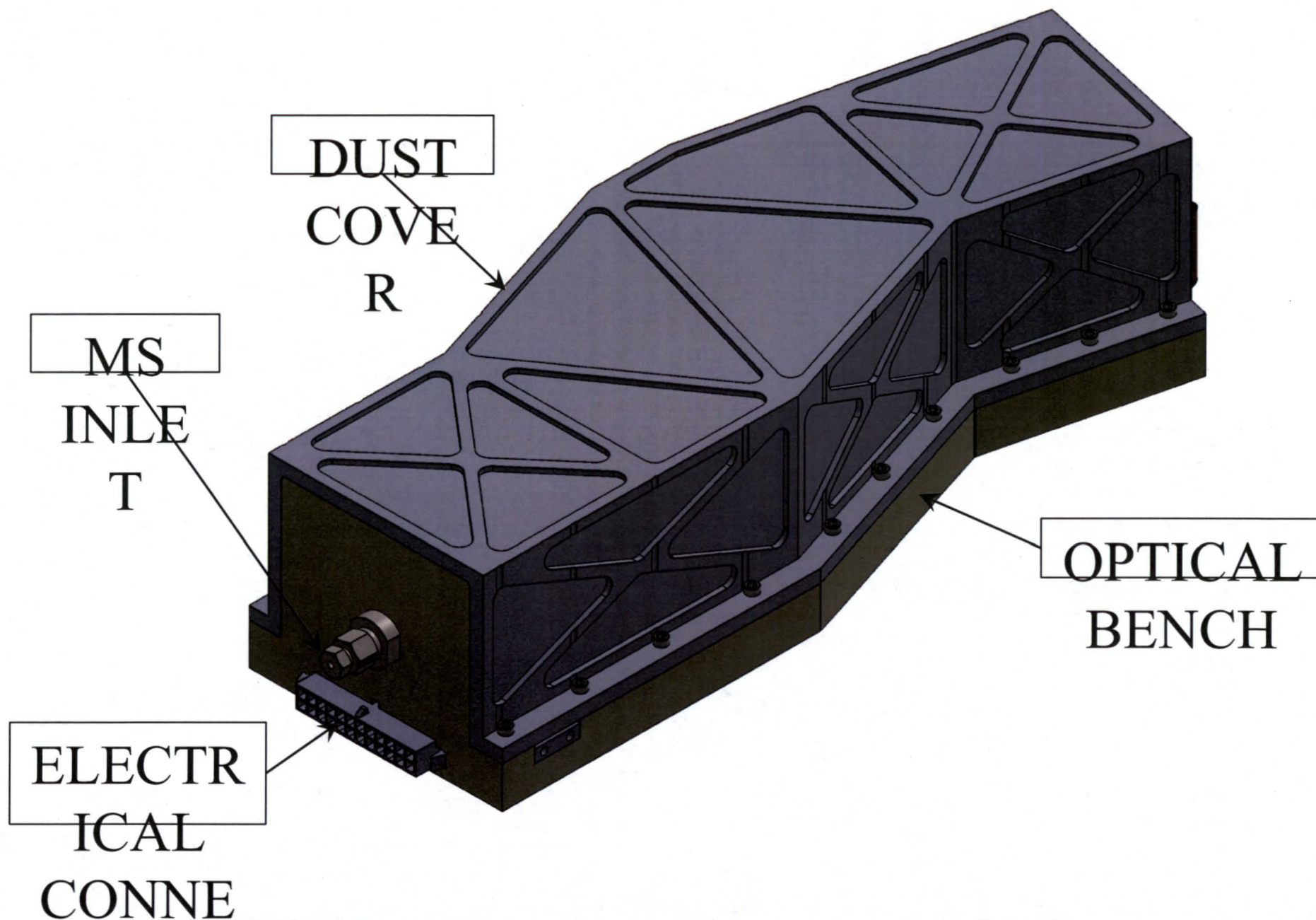




Mass Spectrometer

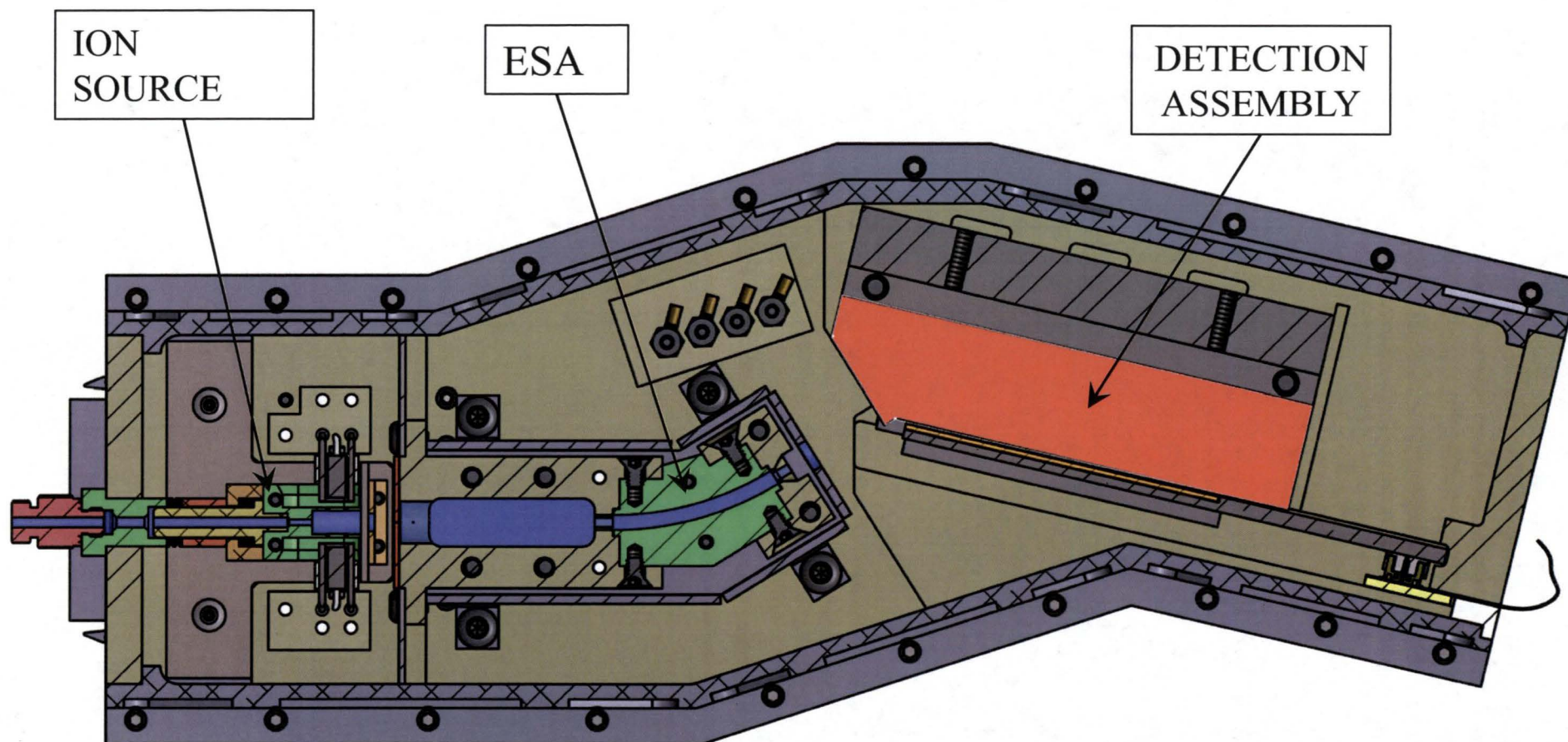



Mass Spectrometer

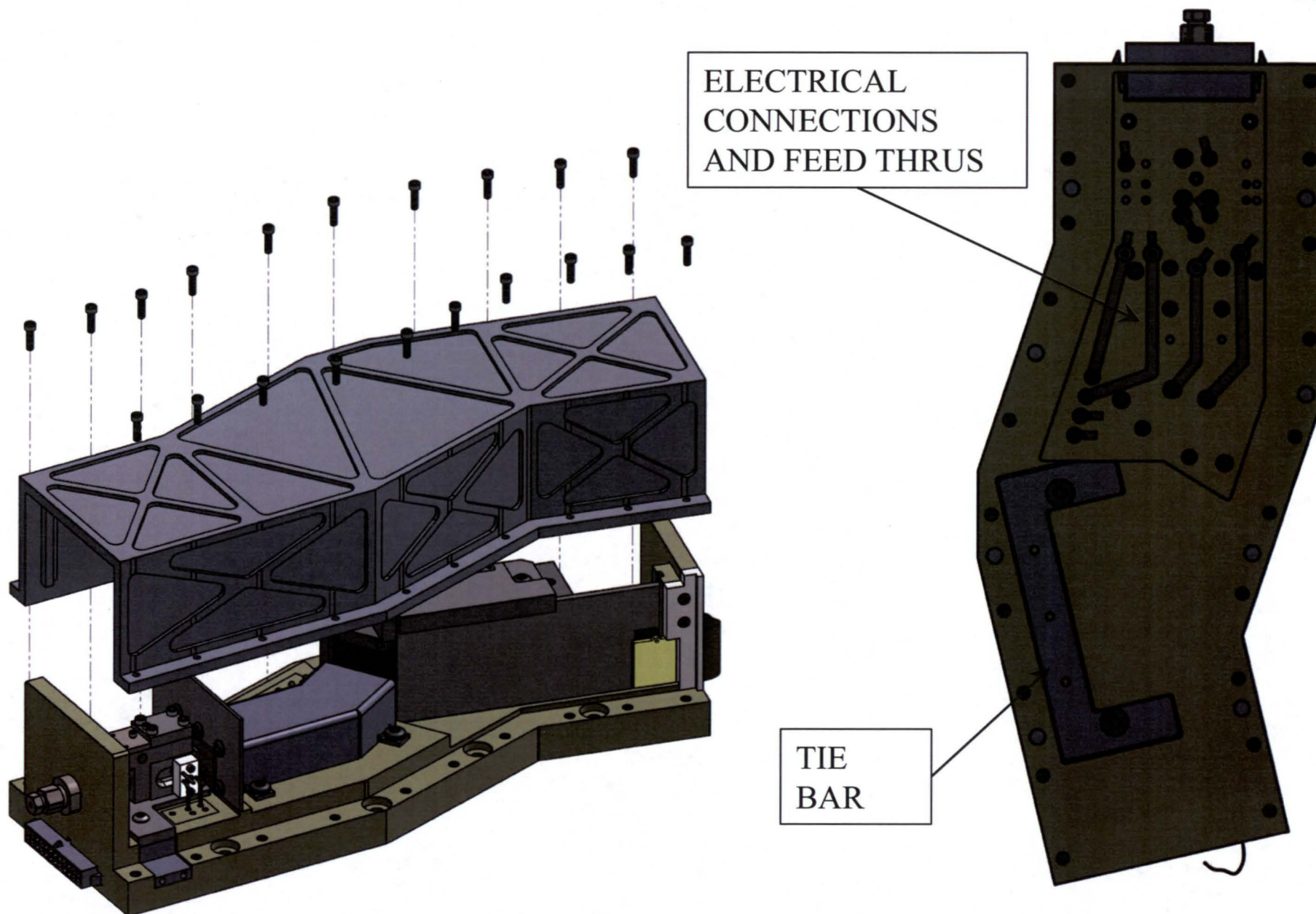




Ion Source With Column Transfer



- Mass = 0.067 kg
- Thermal
 - Heat Source: 15 W cartridge heater
 - Max Temperature: 200 C
- Column Transfer
 - Beswick compression fitting: 1/16" compression with 0.067" dia bore
 - Ferrule: Restek 1/16" dual column, Vespel/Graphite
 - Column transfer fitting heated indirectly from ion source heater
 - Alumina 92% electrical insulator used to isolate fitting from 1000 VDC on ionizer block
- Ionizer Block
 - Ionization volume unchanged from existing IONCAM systems




Ion Source With Column Transfer

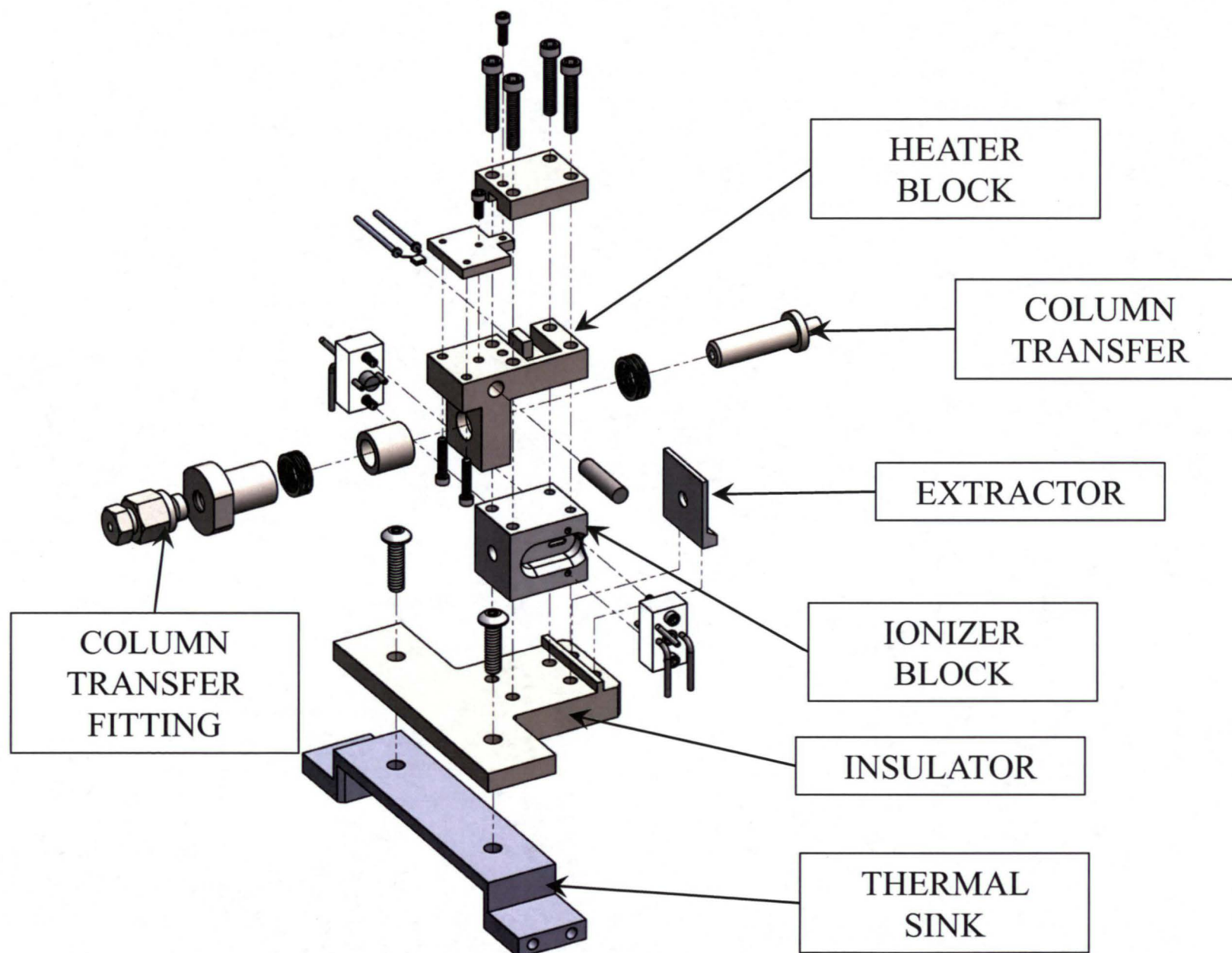



Ion Source With Column Transfer

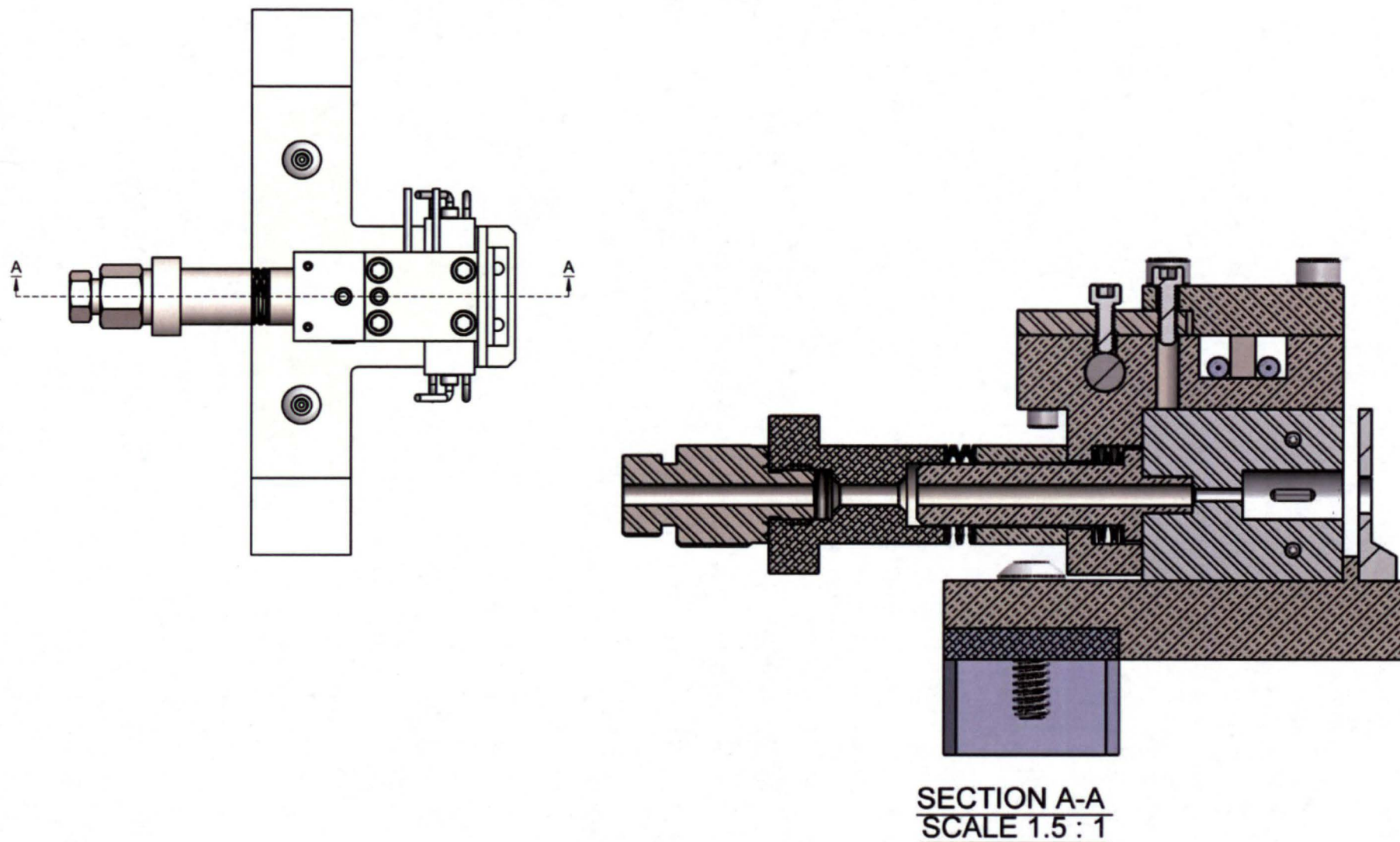




ESA

CHAMB
ER
DIVIDER


ESA
INSERTS

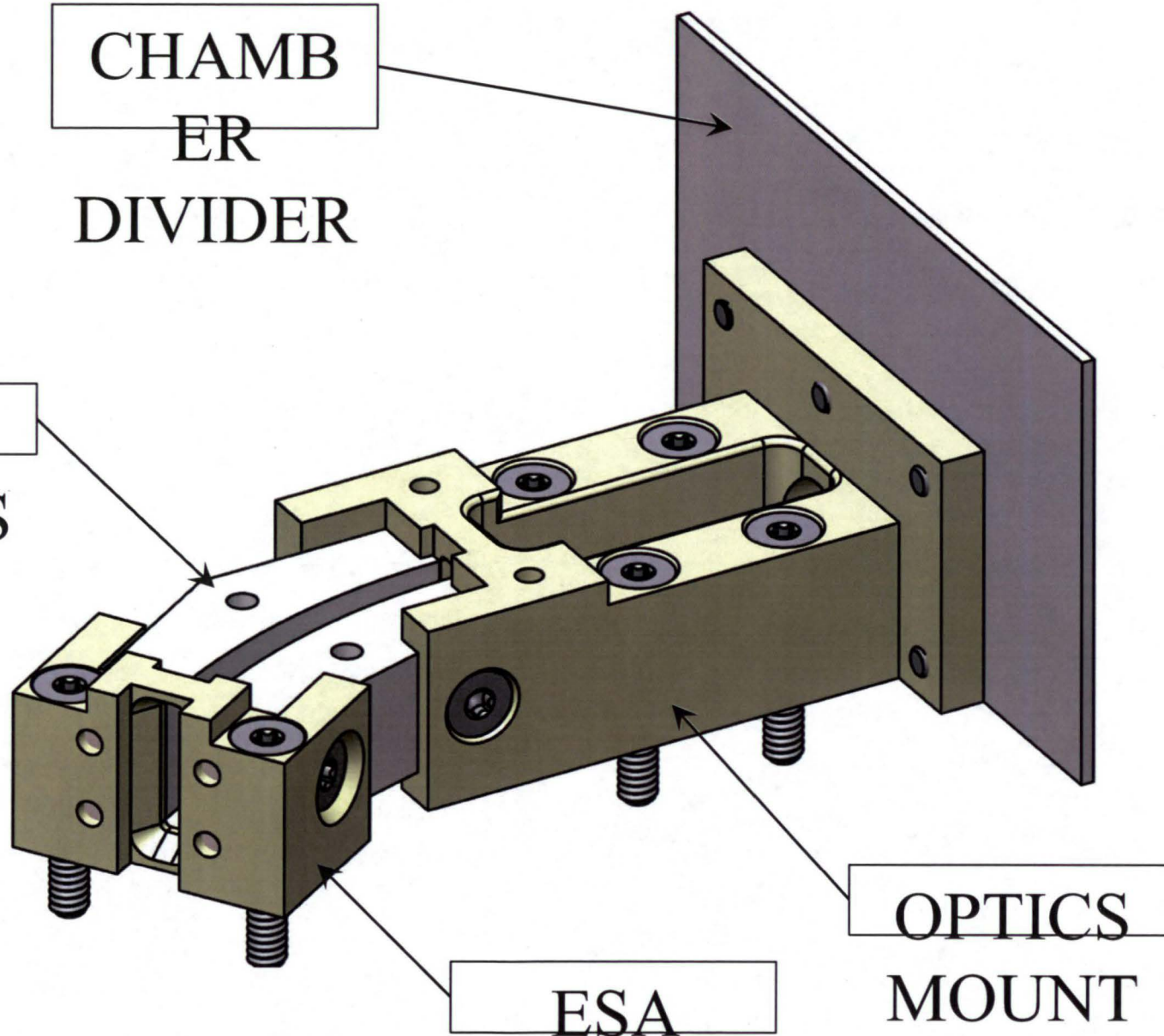

ESA
BRACKET



Detector Assembly

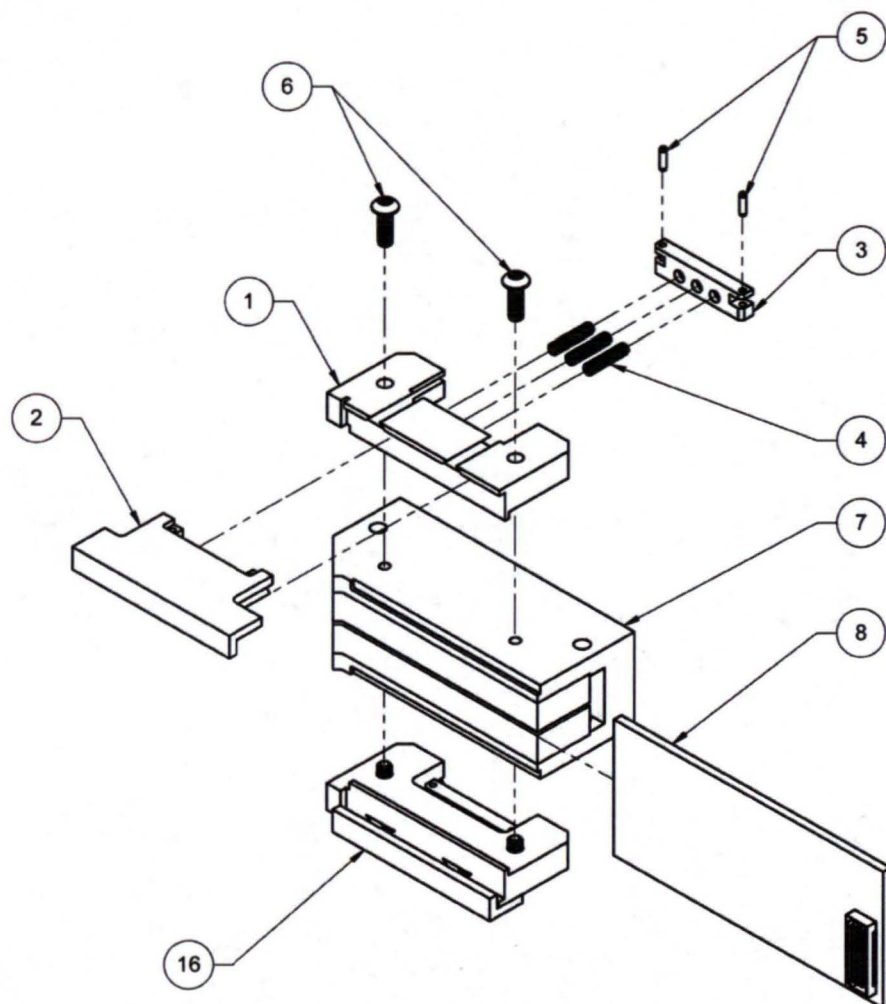

ITEM NO.	QTY.	DESCRIPTION	Material
1	1	TOP CLAMP SLIDE BASE	6061-T6
2	1	TOP CLAMP SLIDE	6061-T6
3	2	CLAMP SPRING RETAINER	6061-T6
4	6	SPRING, COMPRESSION	MUSIC WIRE OR 302SS
5	4	ROLL PIN	4xx SS
6	4	SCREW, #6-32	18-8 SS
7	1	YOKE	ASTM A36 Steel
8	1	CCD BOARD ASSY	ALUMINA
16	1	BOTTOM CLAMP SLIDE	6061-T6
17	1		





- IonCCD, Magnet, Upper and Lower Clamp
- Mass: 0.611 kg
- Field: 7400G +/-20G



Compact HDMR Magnet






- Mass = 0.503 kg (Δ 0.222 kg)
- Field = 7400G \pm 20G
- Magnetic Gap = 1 mm

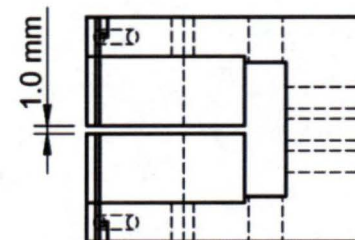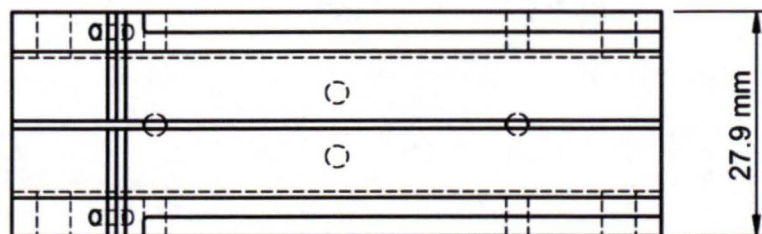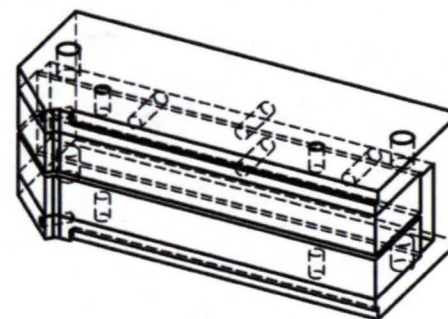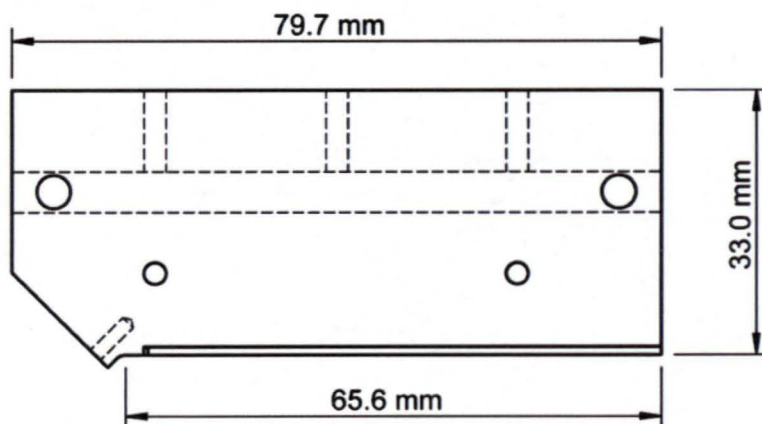

Optical Bench

- Mass = 0.350 kg
- Material: PEEK 30% GF
- Plating
 - Base layer: Copper
 - Top Layer: Nickel
- Monolithic Construction

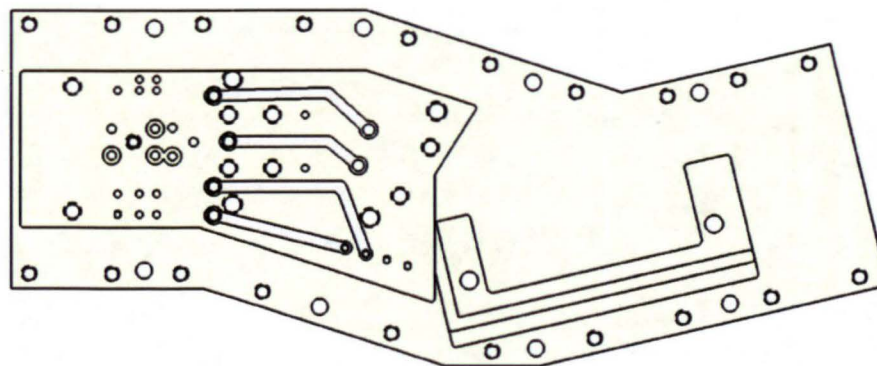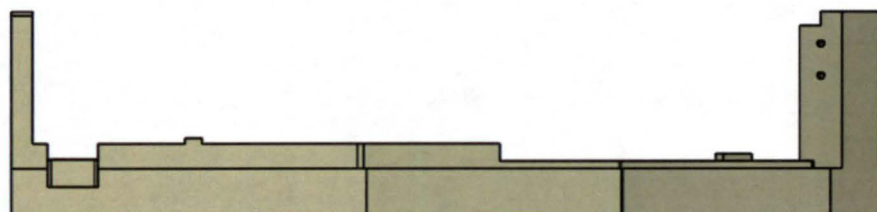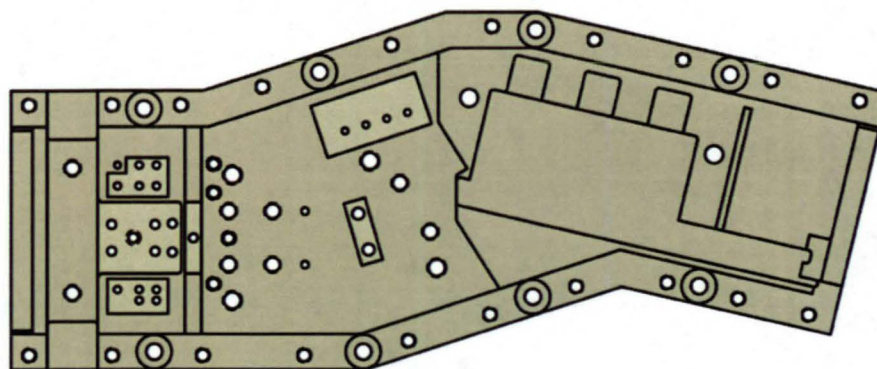




Optical Bench



- Traces for some electrical connections will be machined into plating
- Plating will provide continuous ground path similar to the all metal vacuum chamber of the IONCAM



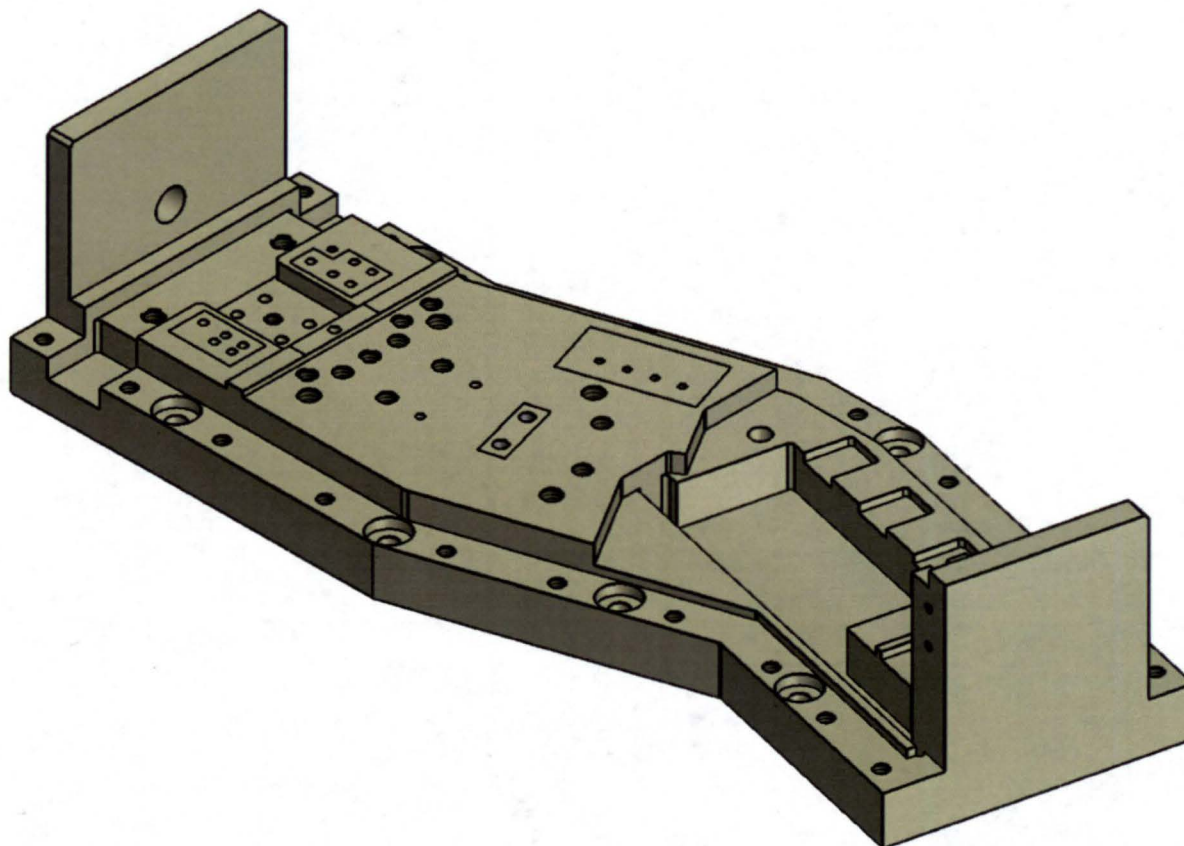




Dust Cover




- Mass = 0.337 kg
- Material: 6061 Al
- Unvented
- Monolithic construction

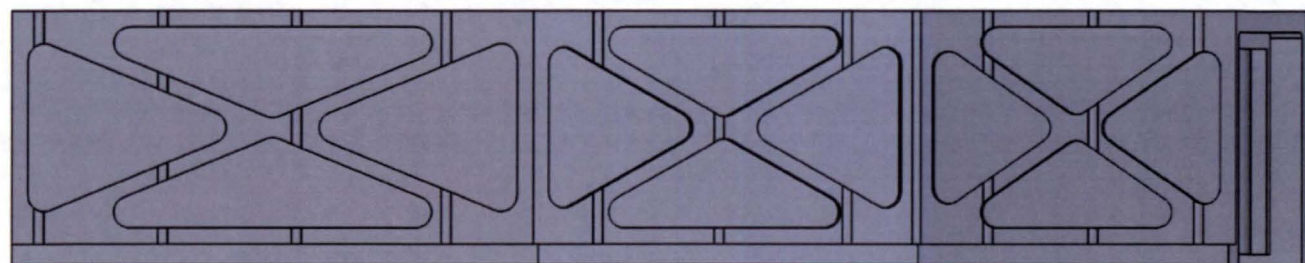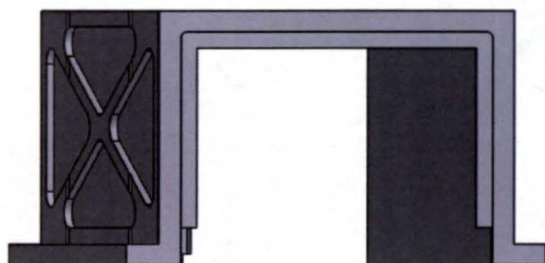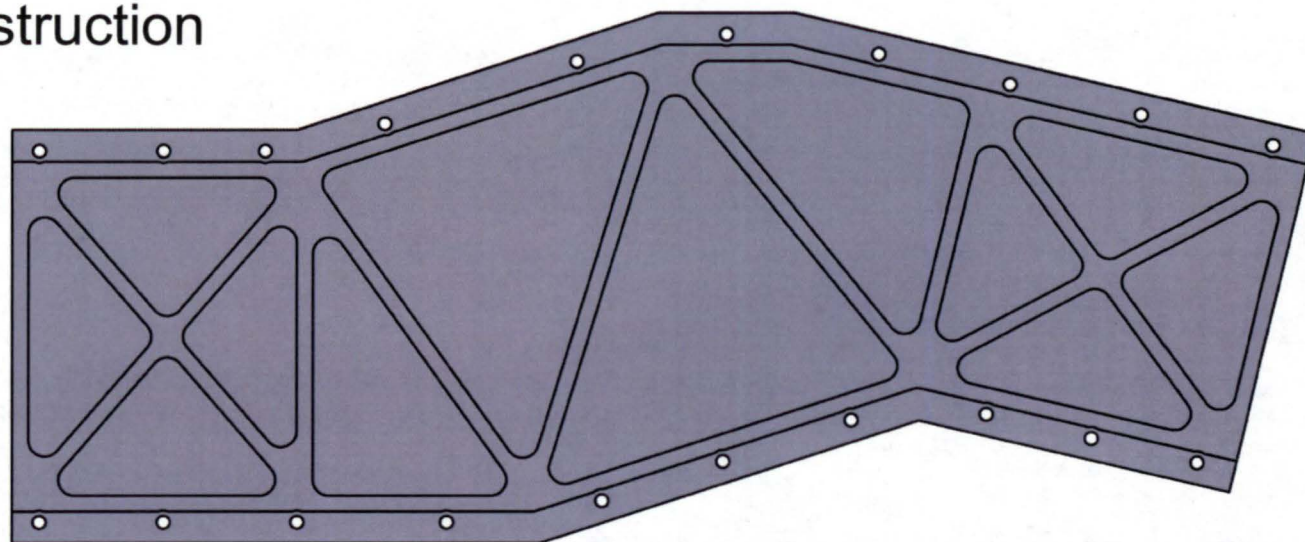




Mounting

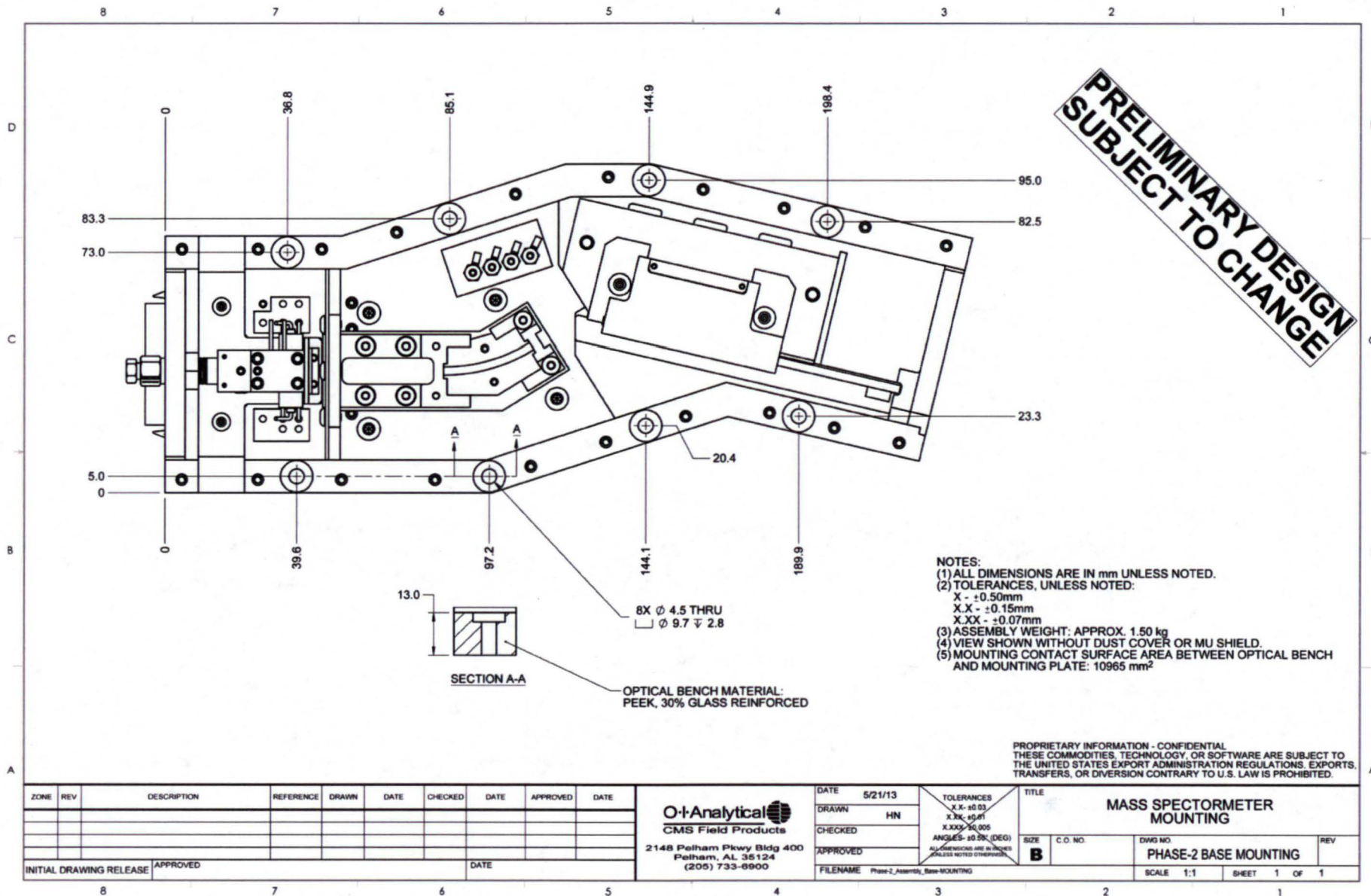



Electrical Connection

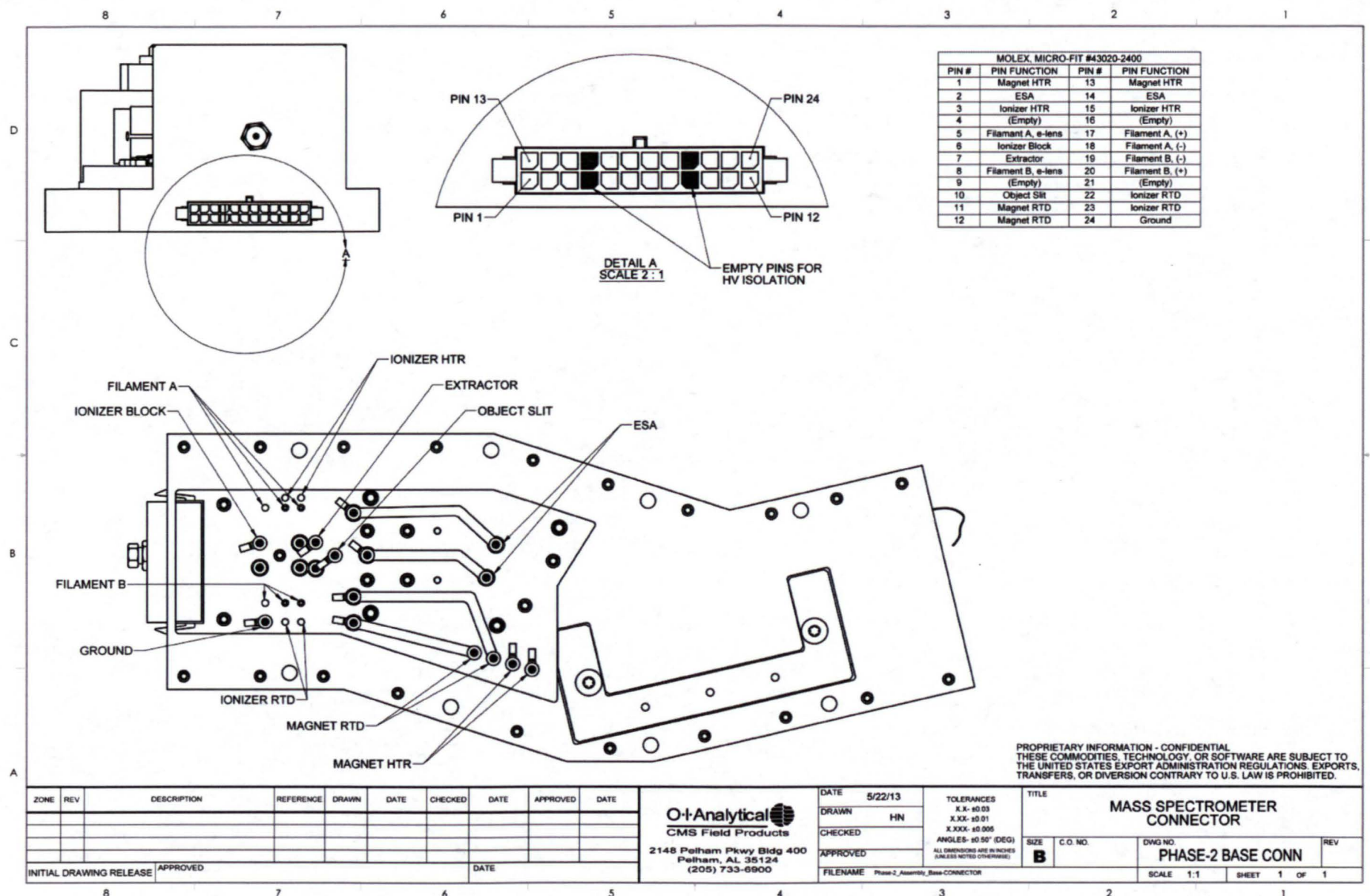



Thermal Study +40°C Ambient


- 40°C Loads
 - Static
 - Heat Loads
 - 3W filament wire
 - 0.2W at CCD face (0.1W for CCD and 0.1W for board components)
 - 13.75W for cartridge heater
 - Radiation
 - Large Peek faces; Inside of the dust cover; Ion source mounting block and heater body
 - Temperature
 - 40°C on the bottom face of the bench
 - -10°C on the bottom face of the magnet tie bar (for certain simulations)
 - Transient
 - Heat Loads
 - 3W filament wire
 - 0.2W at CCD face (0.1W for CCD and 0.1W for board components)
 - 15W for cartridge heater, which is controlled by a thermostat set at the RTD location to keep the temperature between 205 and 210
 - Radiation
 - Same as above
 - Temperature
 - 40°C on the bottom face of the bench
 - 40°C initial temperature



Thermal Study +40°C Ambient



- With tie-bar at -10°C

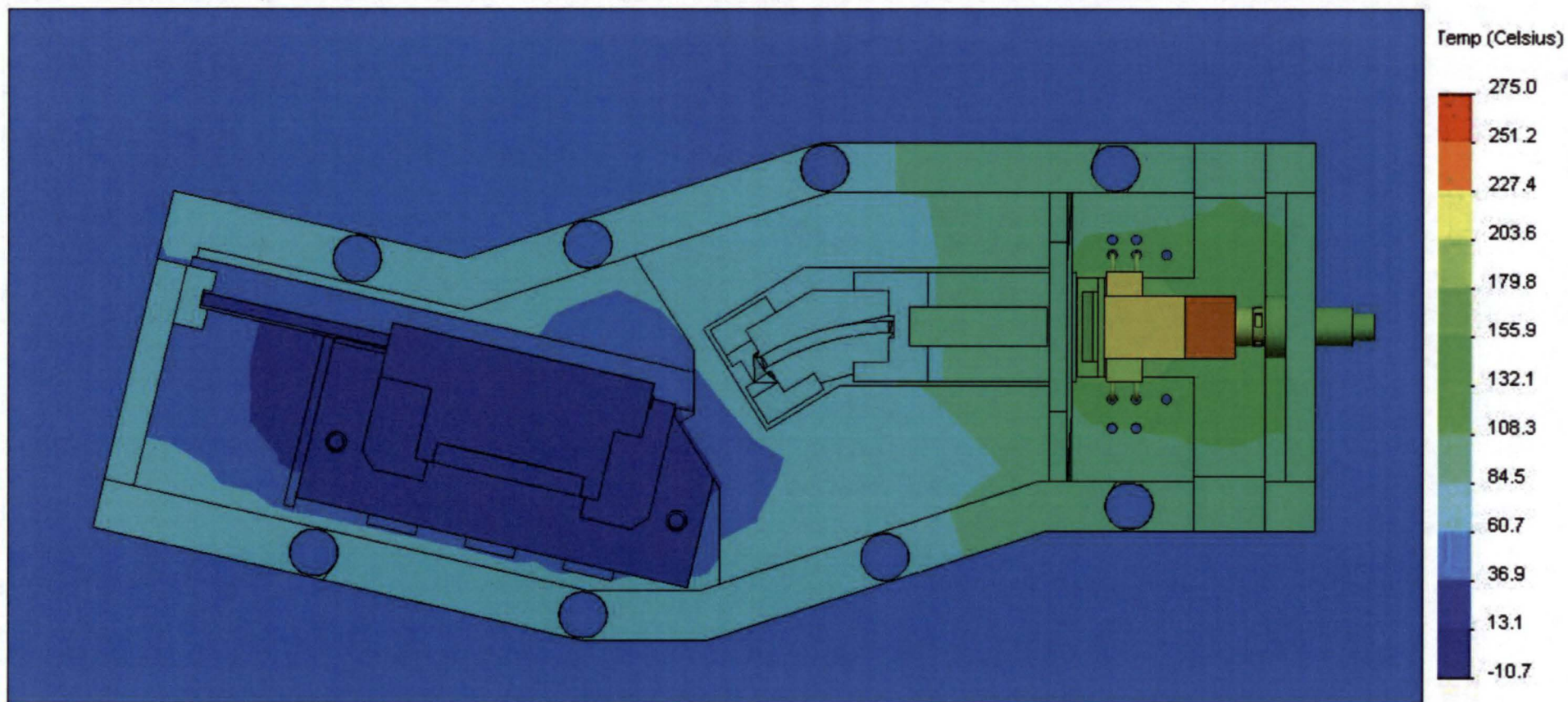




Thermal Study +40°C Ambient



- With tie-bar at -10°C

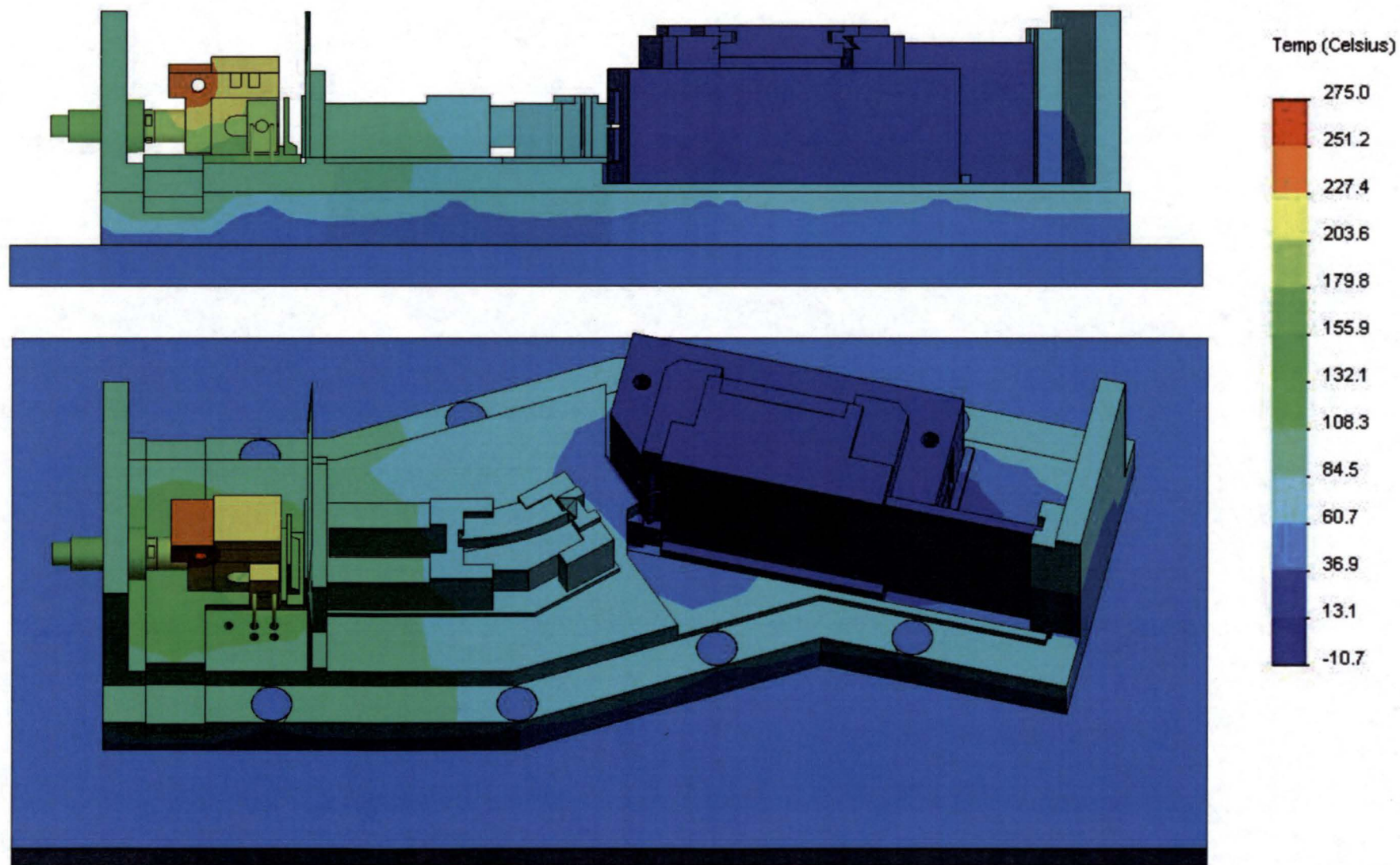




Thermal Study +40°C Ambient



- With tie-bar at -10°C

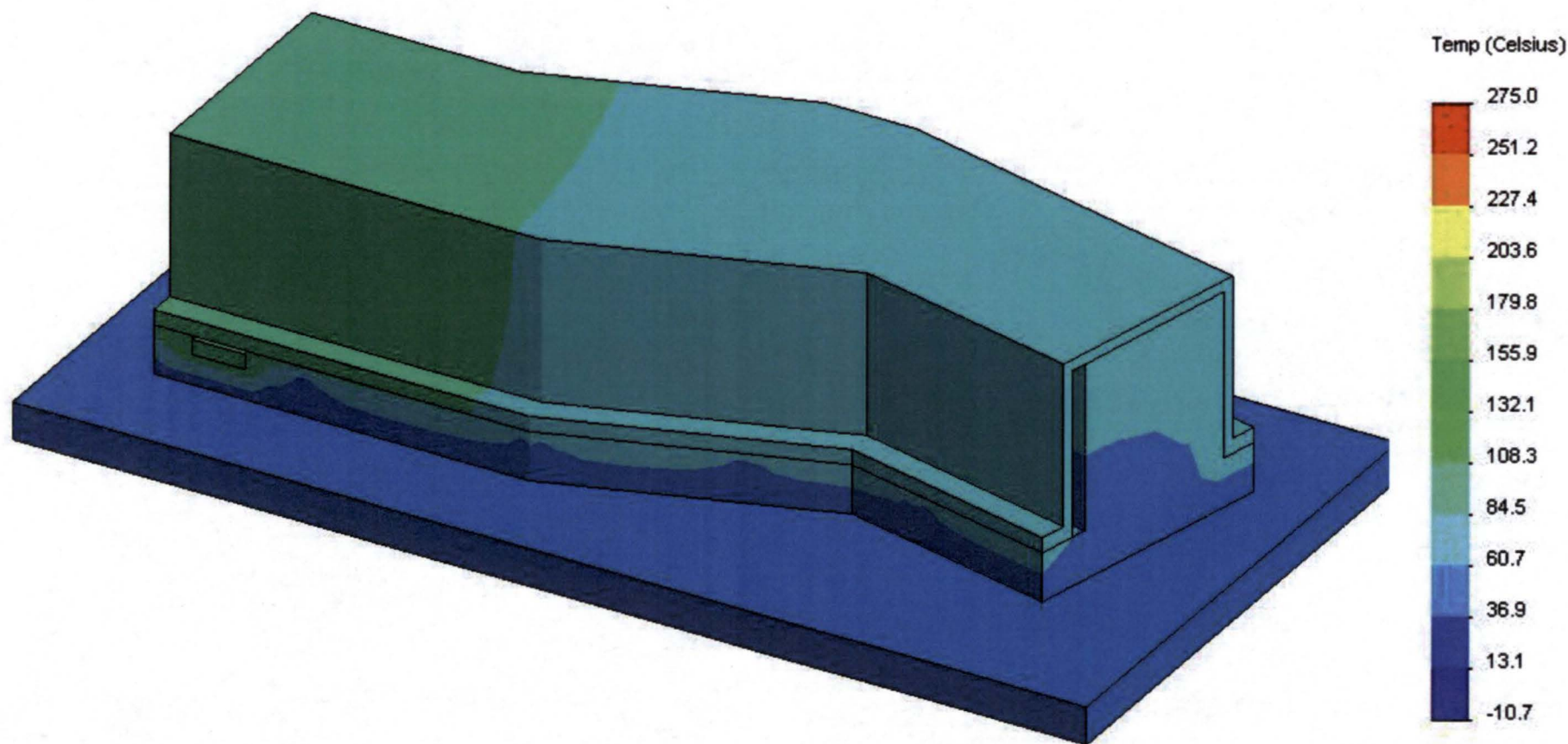




Thermal Study +40°C Ambient



- With NO tie-bar heat sink

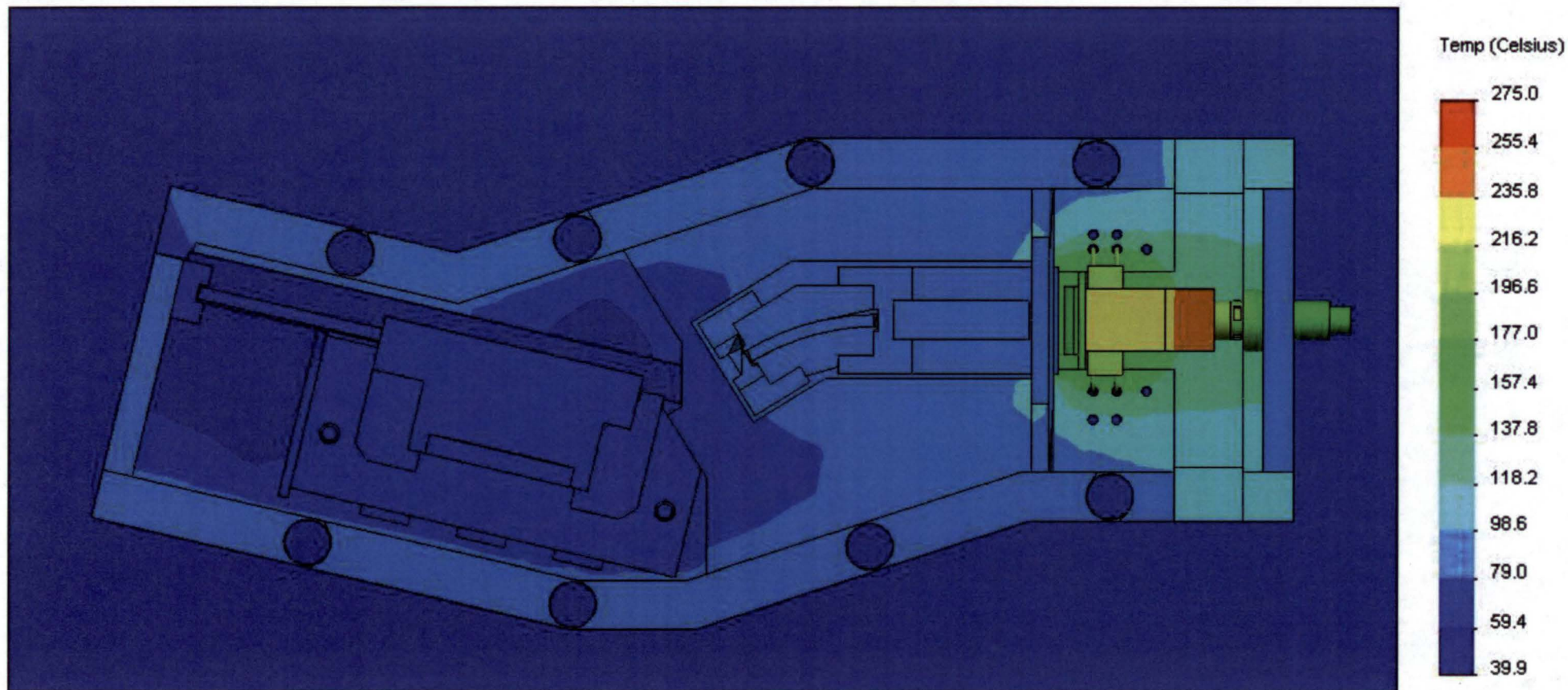




Thermal Study +40°C Ambient



- With NO tie-bar heat sink

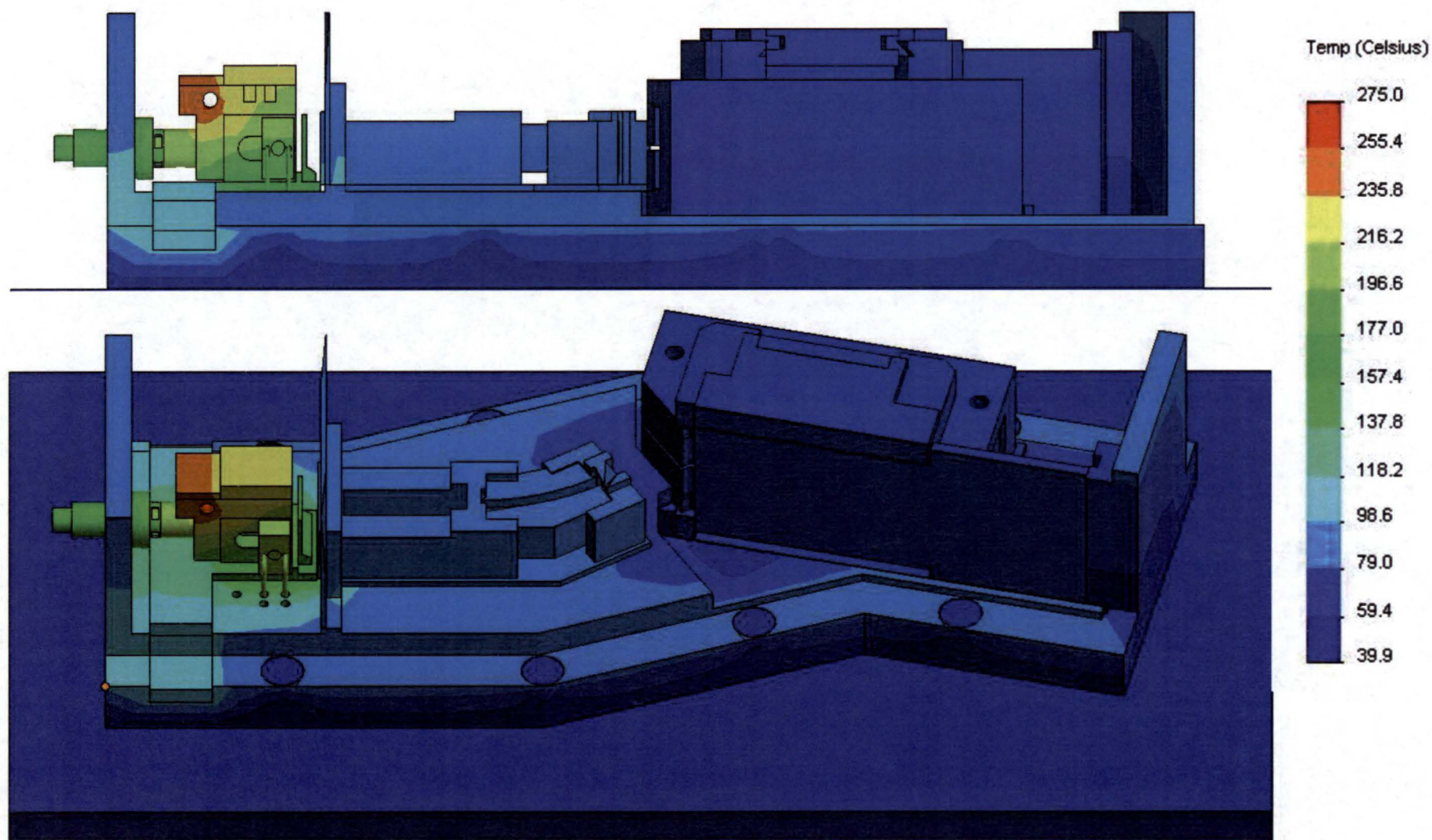




Thermal Study +40°C Ambient



- With NO tie-bar heat sink
- Ionization volume temperature virtually unaffected by presence of tie bar under the Detector Assembly

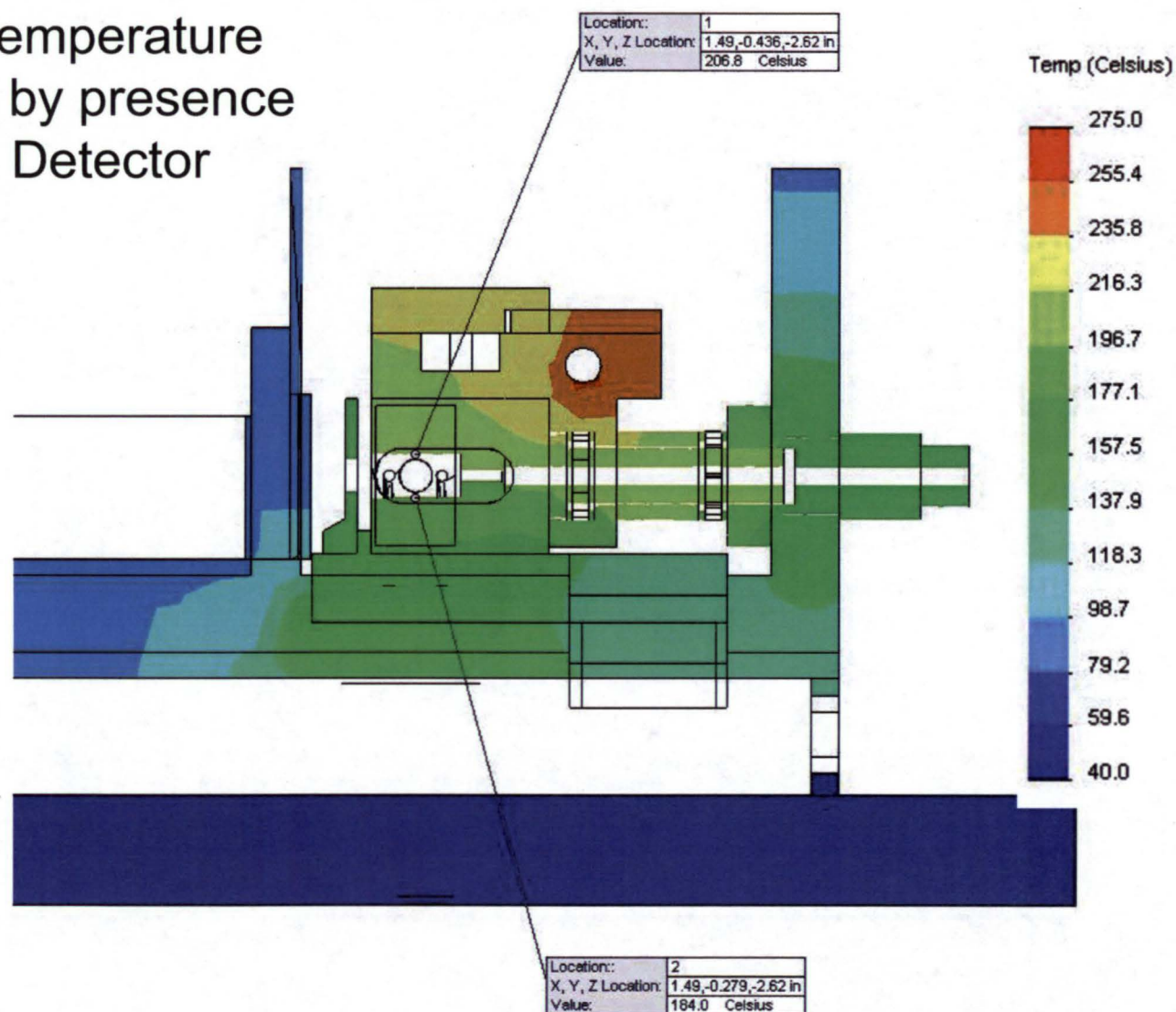



Thermal Study +40°C Ambient




Study name: Transient03
Plot type: Thermal Thermal2


17:10 Monday June 03 2013

X Title: Time (sec)
Y Title: Temp (Celsius)

Point	X	min	Y1 (Location 1)	
1	0	0	40	4000
2	240	4	99.037	9903.7
3	480	8	131.02	13102
4	720	12	149.57	14957
5	960	16	161.59	16159
6	1200	20	170.08	17008
7	1440	24	176.49	17649
8	1680	28	181.56	18156
9	1920	32	185.71	18571
10	2160	36	189.2	18920
11	2400	40	192.17	19217
12	2640	44	194.75	19475
13	2880	48	196.99	19699
14	3120	52	198.96	19896

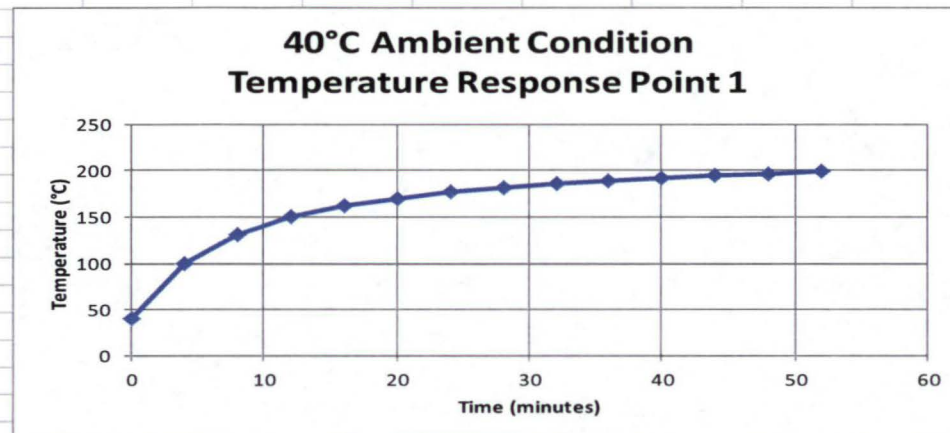

	X	min	Y1 (Location 1)	
1	0	0	40	4000
2	240	4	83.756	8375.6
3	480	8	111.09	11109
4	720	12	127.94	12794
5	960	16	139.19	13919
6	1200	20	147.27	14727
7	1440	24	153.44	15344
8	1680	28	158.35	15835
9	1920	32	162.41	16241
10	2160	36	165.83	16583
11	2400	40	168.75	16875
12	2640	44	171.29	17129
13	2880	48	173.51	17351
14	3120	52	175.46	17546

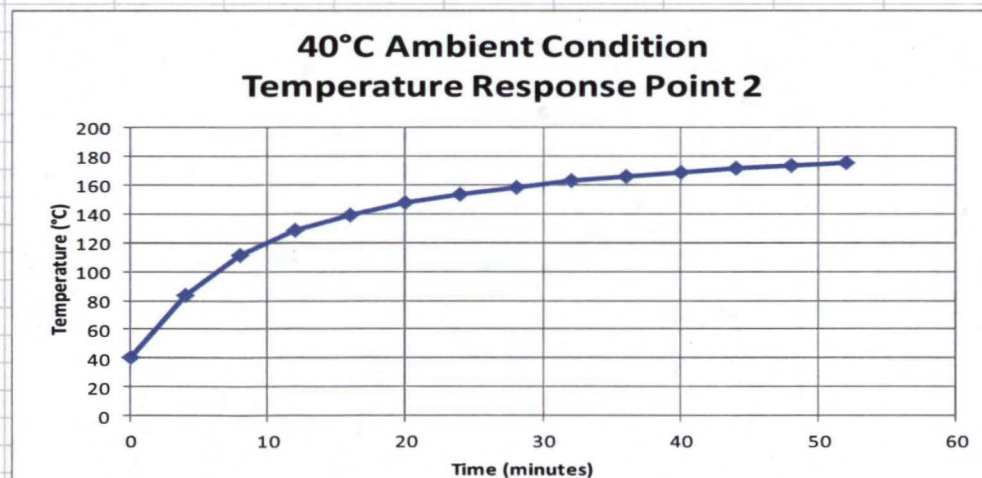



Thermal Study +10°C Ambient




- 10°C Study
 - Static
 - Heat Loads
 - 3W filament wire
 - 0.2W at CCD face (0.1W for CCD and 0.1W for board components)
 - 15W for cartridge heater
 - Radiation
 - Large Peek faces; Inside of dust cover; Ion source mounting block and heater body
 - Temperature
 - 10°C on the bottom face of the bench
 - -10°C on the bottom face of the magnet tie bar (for certain simulations)
 - Transient
 - Heat Loads
 - 3W filament wire
 - 0.2W at CCD face (0.1W for CCD and 0.1W for board components)
 - 15W for cartridge heater, which is controlled by a thermostat set at the RTD location to keep the temperature between 205 and 210
 - Radiation
 - Same as above
 - Temperature
 - 10°C on the bottom face of the bench
 - 10°C initial temperature



Thermal Study +10°C Ambient



- With tie-bar at -10°C

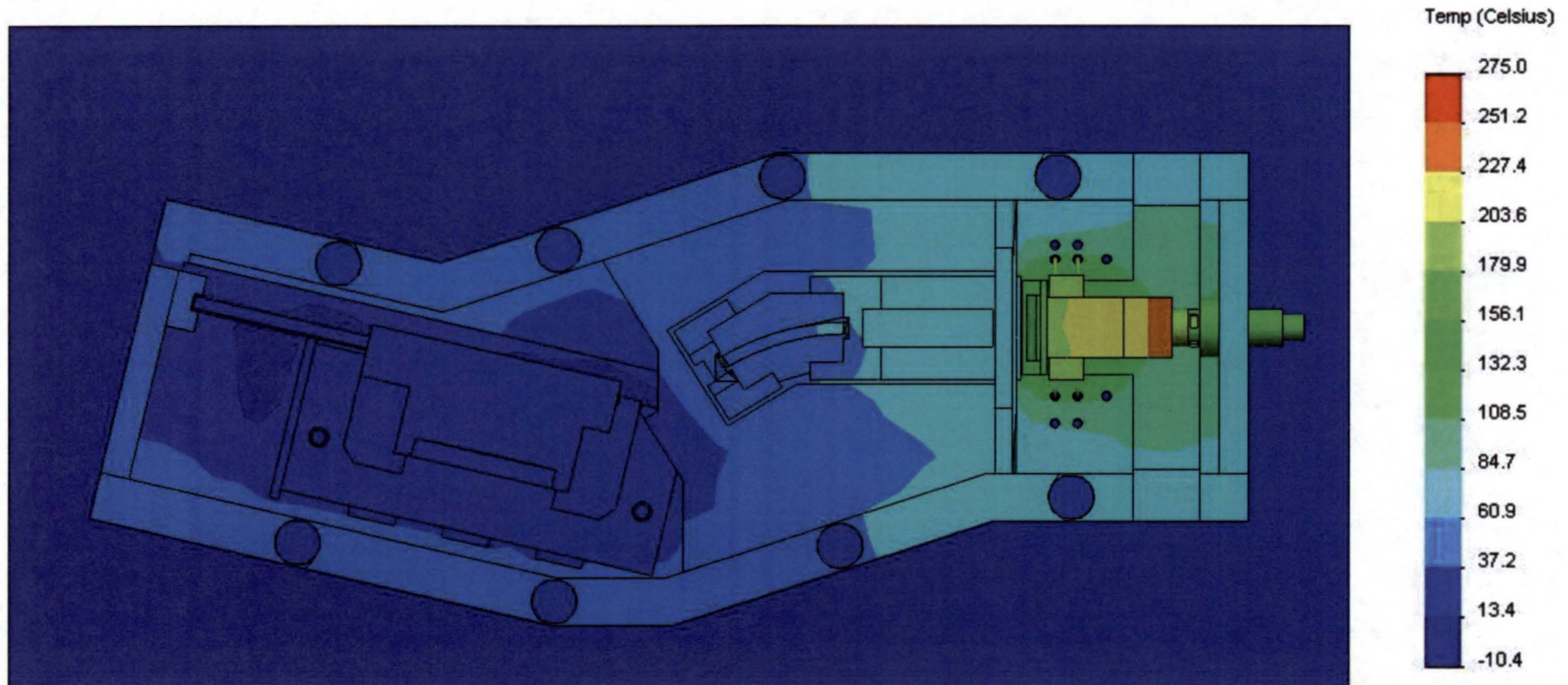




Thermal Study +10°C Ambient



- With tie-bar at -10°C

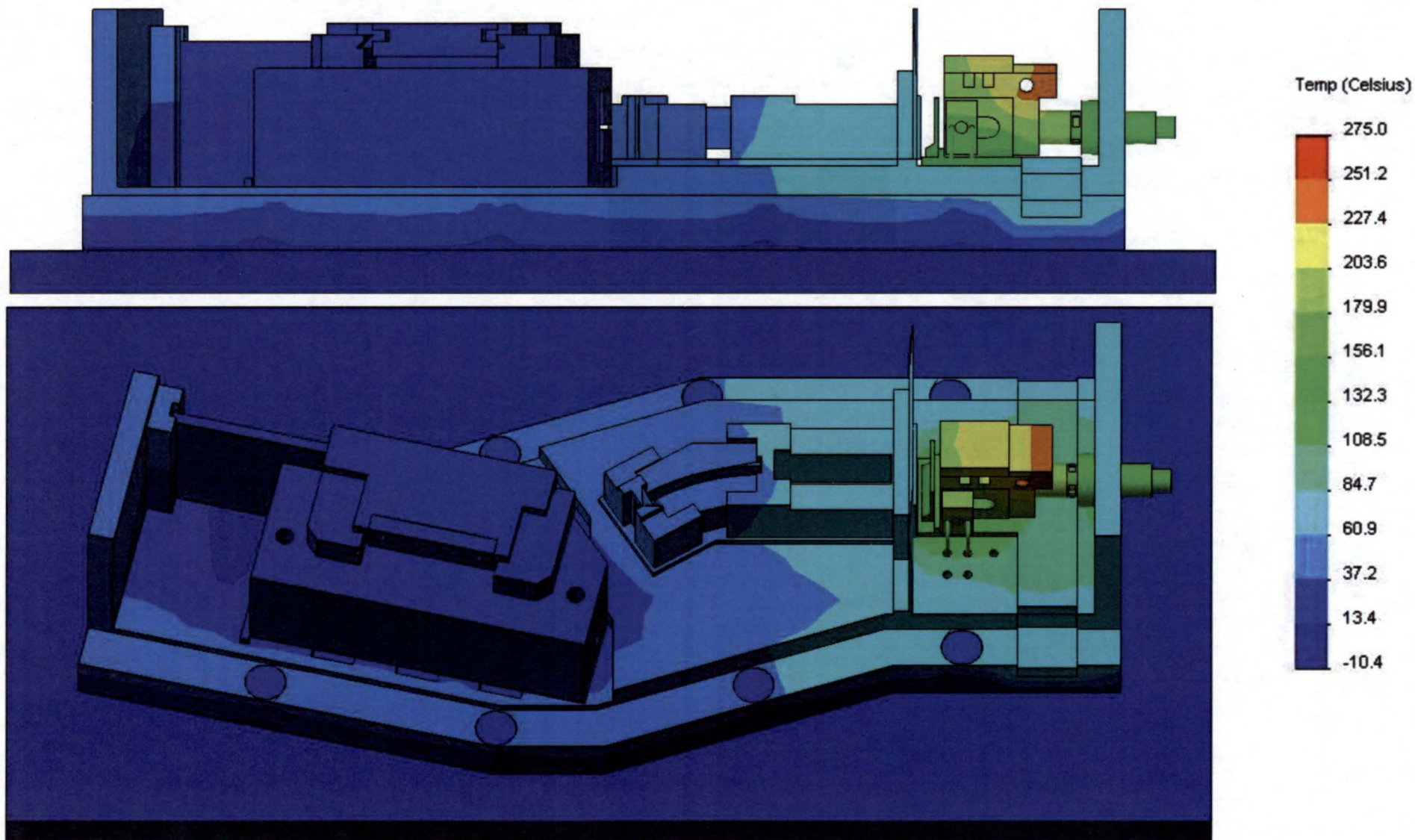




Thermal Study +10°C Ambient



- With tie-bar at -10°C

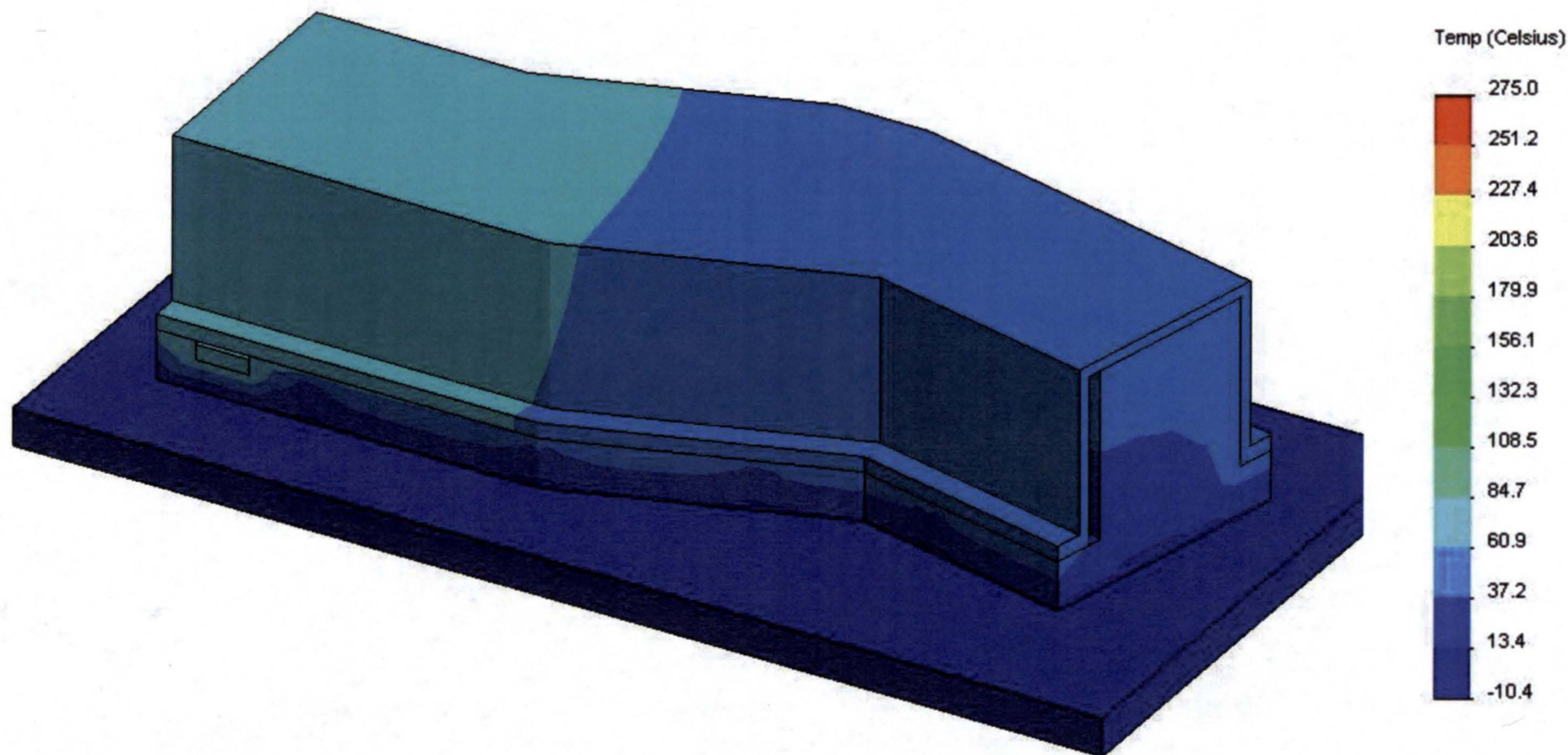




Thermal Study +10°C Ambient



- With tie-bar heat sink
- Ionization volume temperature virtually unaffected by presence of tie bar under the Detector Assembly

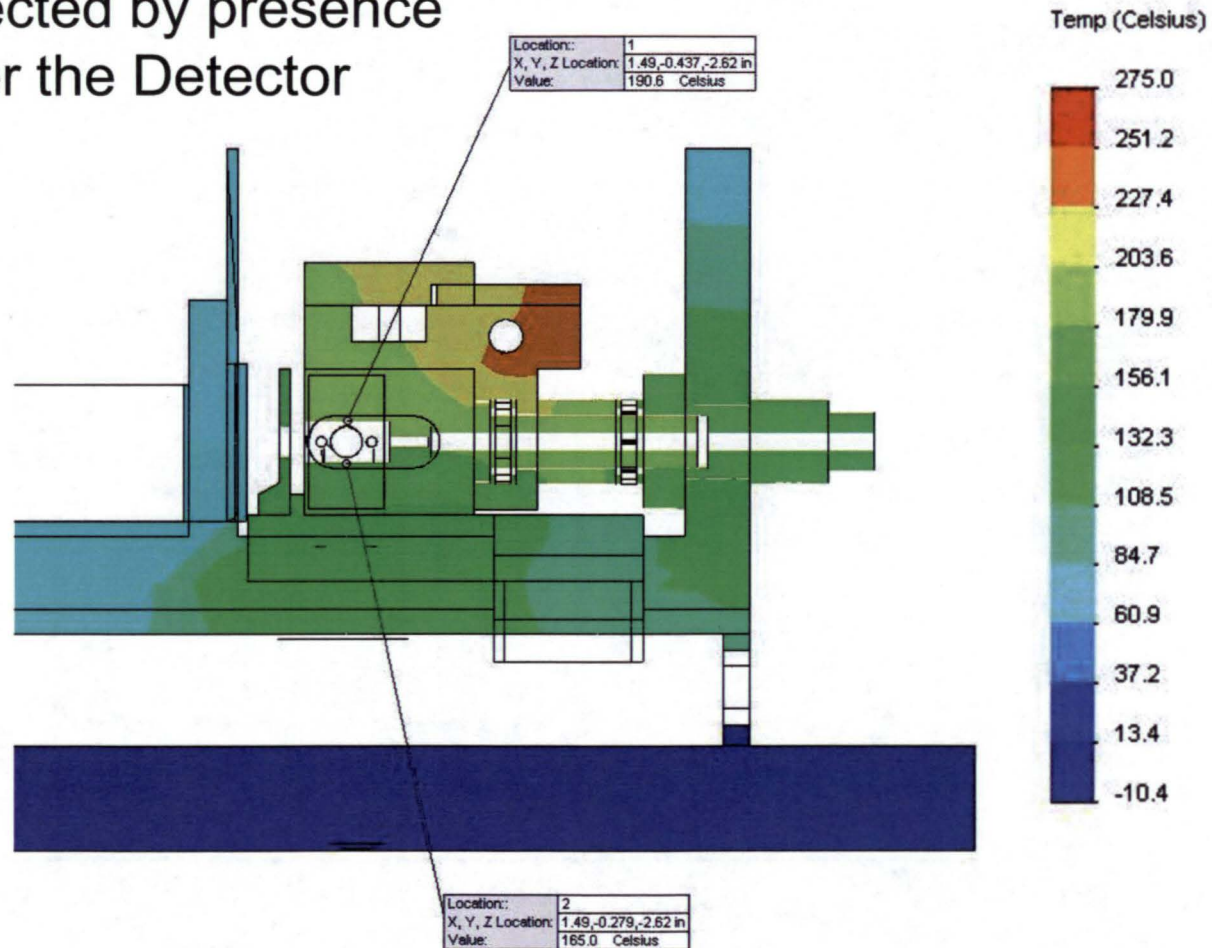




Thermal Study +10°C Ambient



- With NO tie-bar heat sink

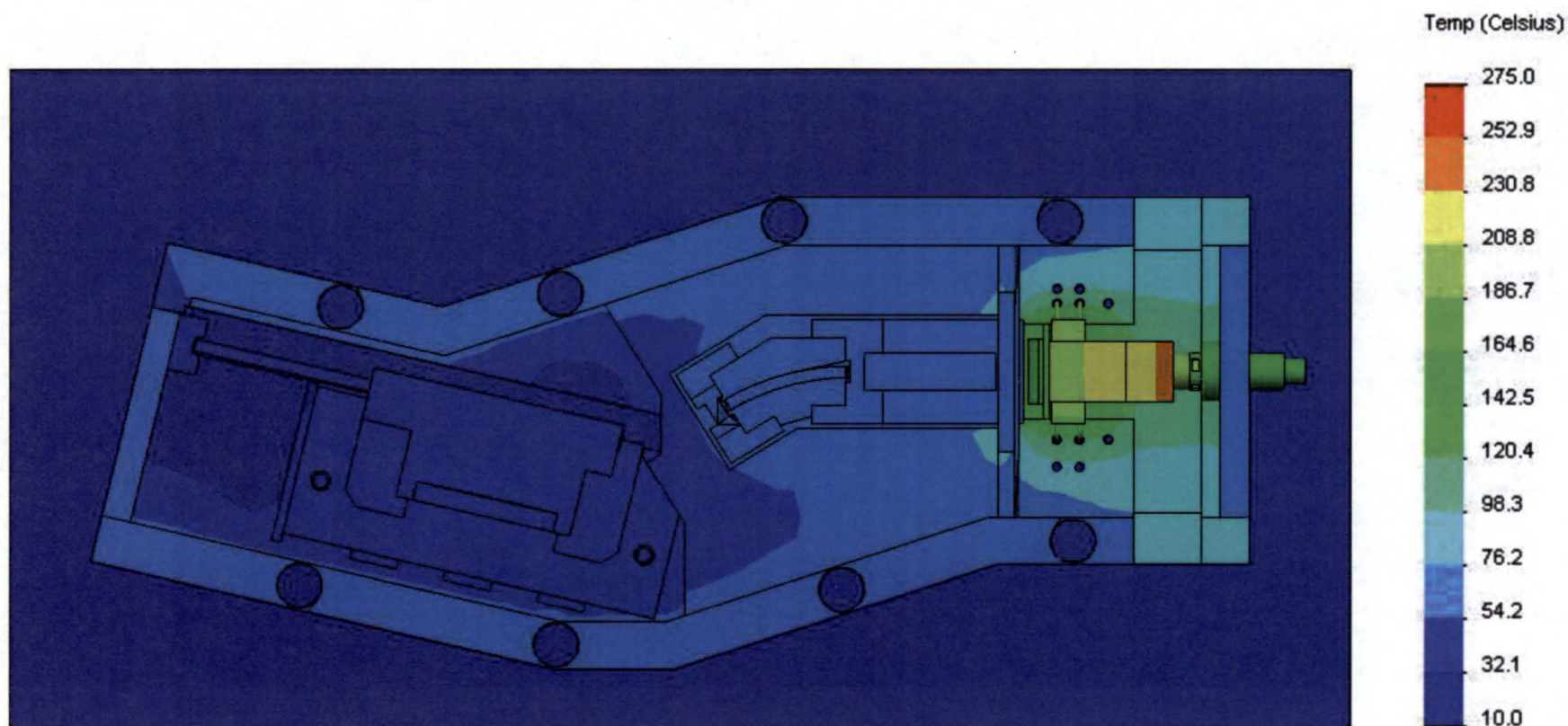




Thermal Study +10°C Ambient



- With NO tie-bar heat sink

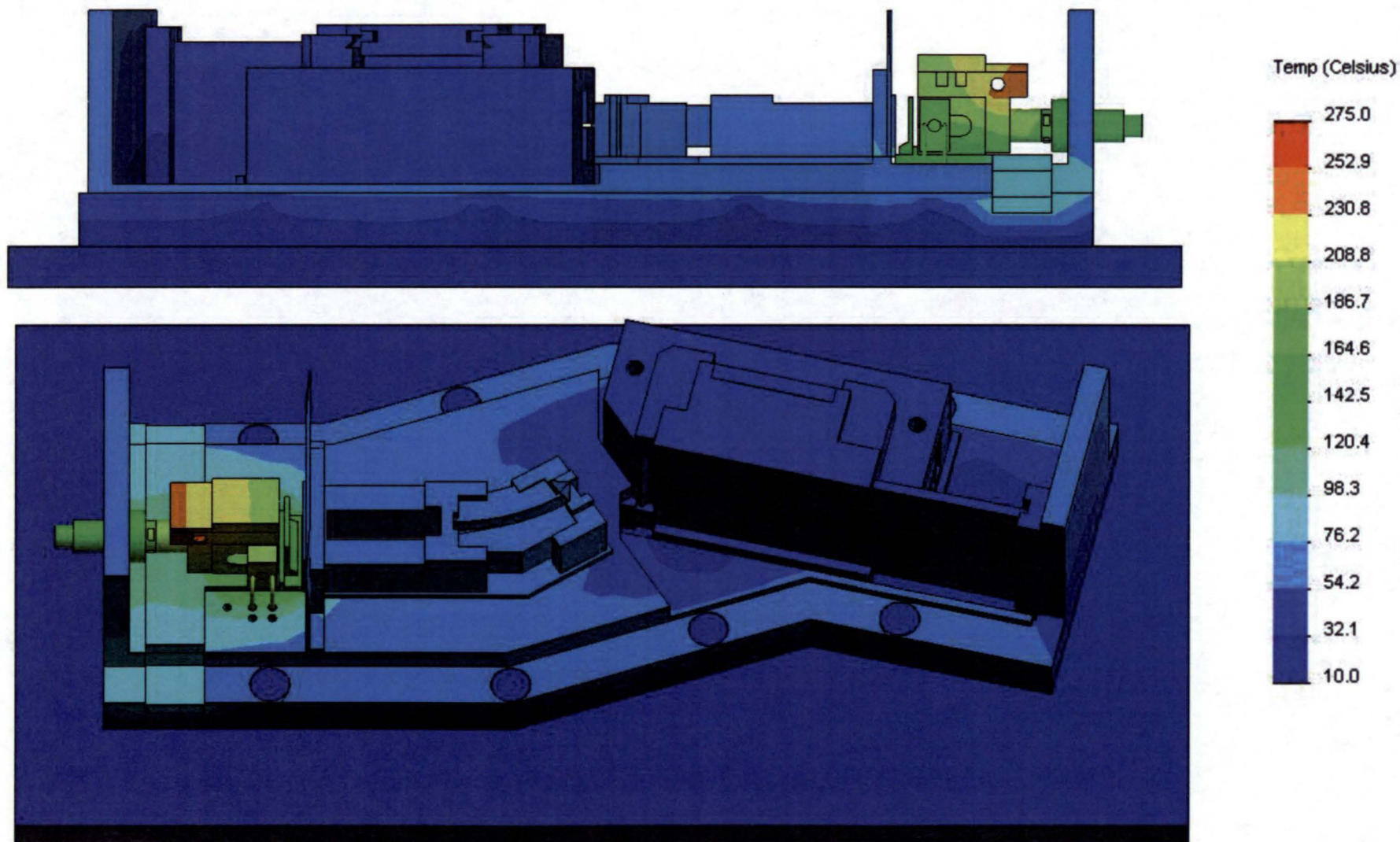




Thermal Study +10°C Ambient



- With NO tie-bar heat sink

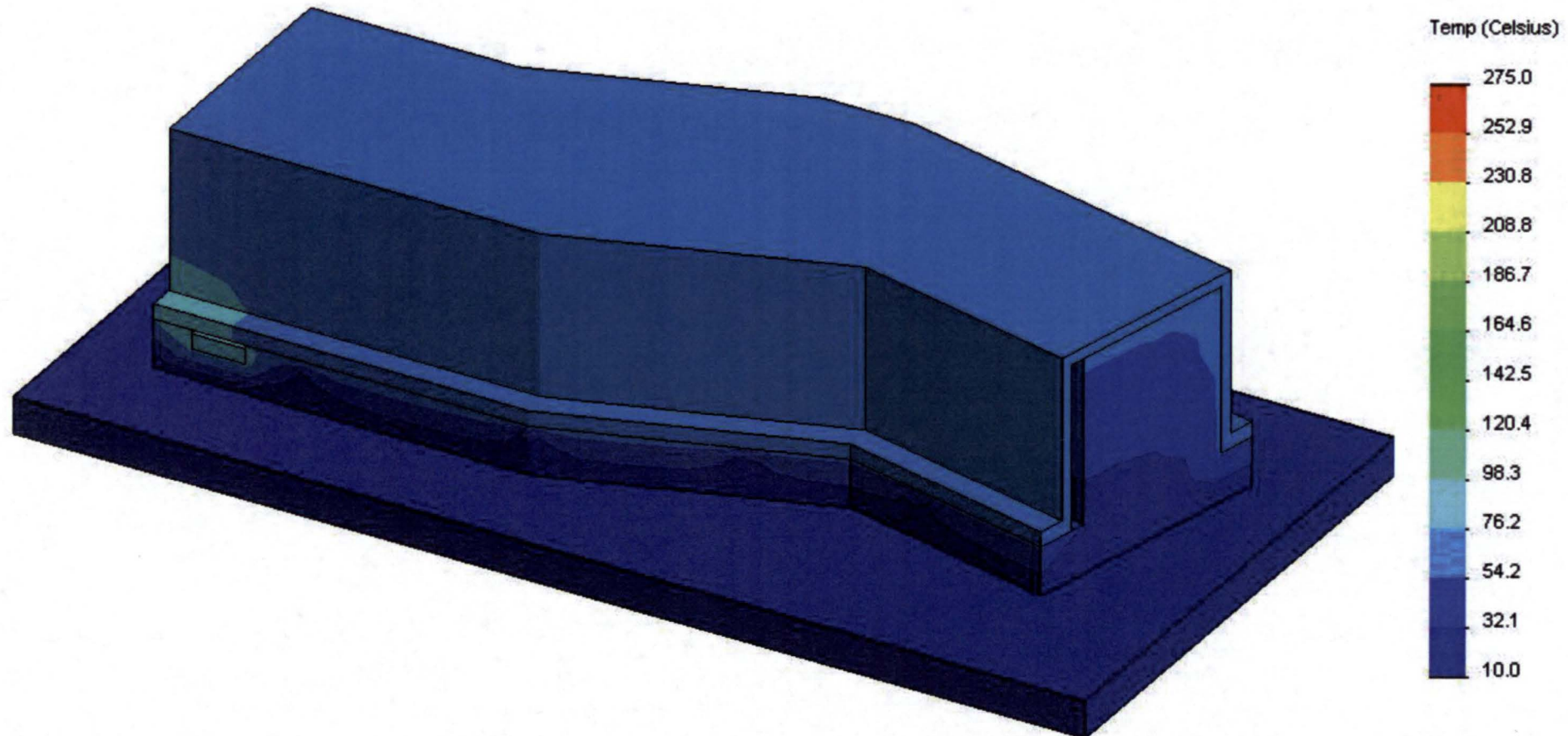




Thermal Study +10°C Ambient




Study name: Transient1-10
Plot type: Thermal Thermal2
19:35 Monday June 03 2013

X Title: Time (sec)
Y Title: Temp (Celsius)

Point	X	Y1 (Location 1)		
		0	10	1000
1	240	69.829	6982.9	
2	480	103.05	10305	
3	720	122.39	12239	
4	960	134.91	13491	
5	1200	143.75	14375	
6	1440	150.41	15041	
7	1680	155.69	15569	
8	1920	160.02	16002	
9	2160	163.66	16366	
10	2400	166.78	16678	
11	2640	169.48	16948	
12	2880	171.85	17185	
13	3120	173.92	17392	
14	3360	175.76	17576	
15	3600	177.4	17740	
16	3840	178.85	17885	
17	4080	180.15	18015	
18	4320	181.32	18132	
19	4560	182.37	18237	
20	4800	183.31	18331	

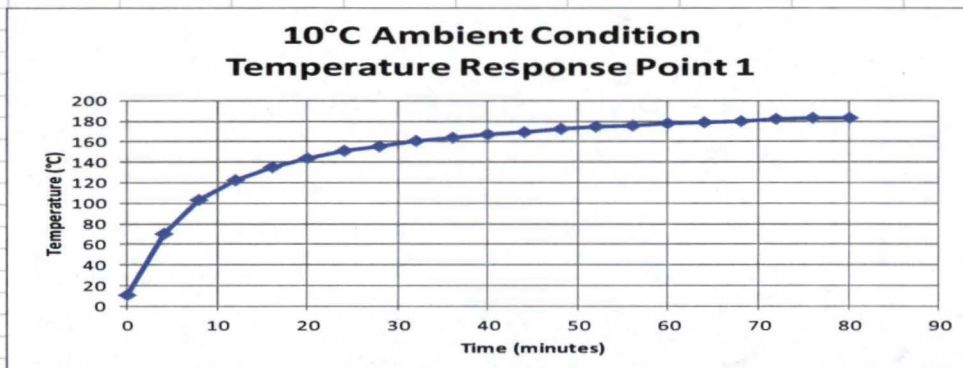

0	0	0	10	1000
1	240	53.787	5378.7	
2	480	81.695	8169.5	
3	720	99.036	9903.6	
4	960	110.63	11063	
5	1200	118.98	11898	
6	1440	125.35	12535	
7	1680	130.45	13045	
8	1920	134.66	13466	
9	2160	138.21	13821	
10	2400	141.27	14127	
11	2640	143.93	14393	
12	2880	146.26	14626	
13	3120	148.31	14831	
14	3360	150.13	15013	
15	3600	151.74	15174	
16	3840	153.19	15319	
17	4080	154.48	15448	
18	4320	155.63	15563	
19	4560	156.67	15667	
20	4800	157.61	15761	

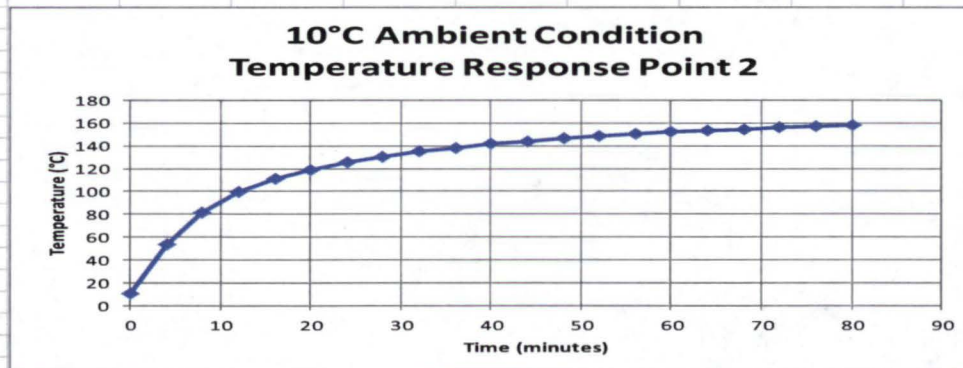



Thermal Study -20°C Ambient




- -20°C Study
 - Static
 - Heat Loads
 - 3W filament wire
 - 0.2W at CCD face (0.1W for CCD and 0.1W for board components)
 - 19W for cartridge heater
 - Radiation
 - Large Peek faces; Inside of dust cover; Ion source mounting block and heater body
 - Temperature
 - -20°C on the bottom face of the bench
 - -10°C on the bottom face of the magnet tie bar (for certain simulations)
 - Transient
 - Heat Loads
 - 3W filament wire
 - 0.2W at CCD face (0.1W for CCD and 0.1W for board components)
 - 20W for cartridge heater, which is controlled by a thermostat set at the RTD location to keep the temperature between 205 and 210
 - Radiation
 - Same as above
 - Temperature
 - -20°C on the bottom face of the bench
 - -20°C initial temperature



Thermal Study -20°C Ambient



- With tie-bar at -10°C

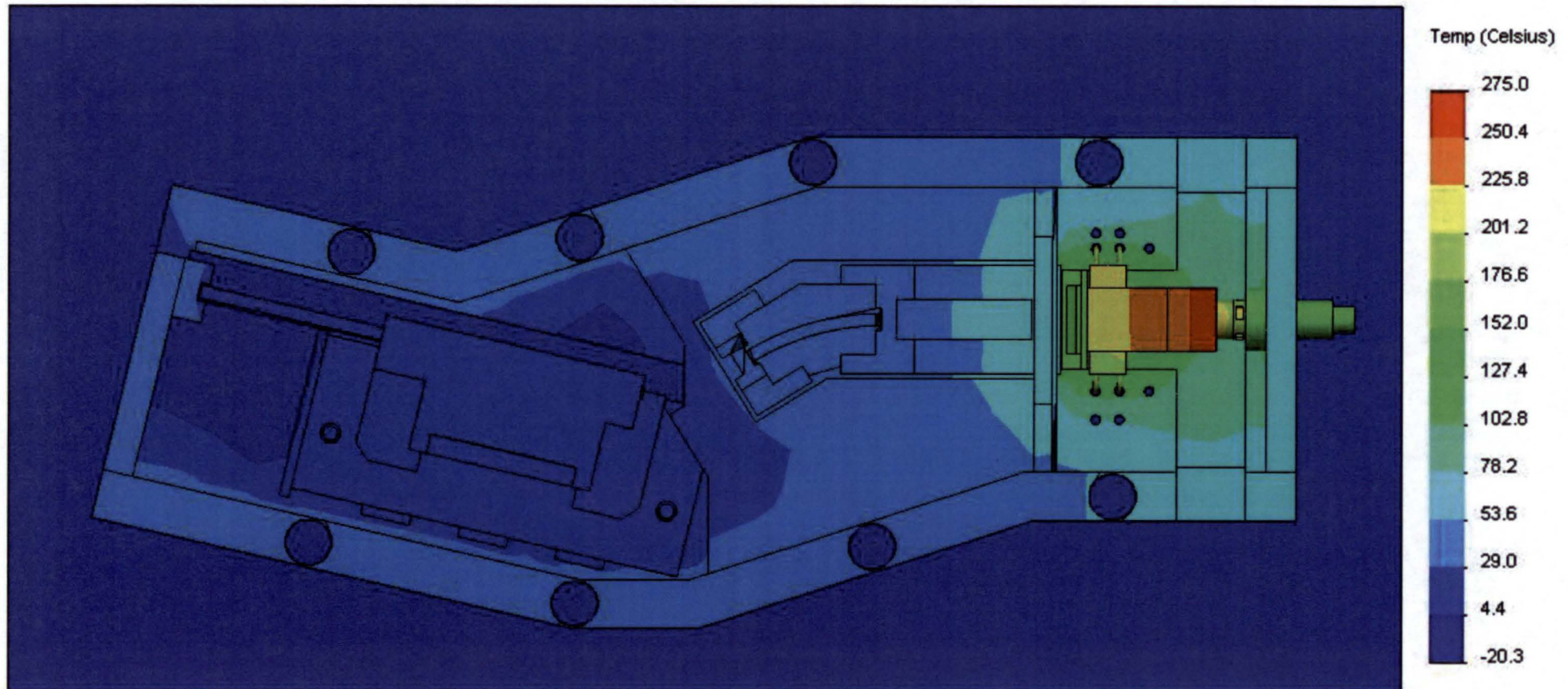




Thermal Study -20°C Ambient



- With tie-bar at -10°C

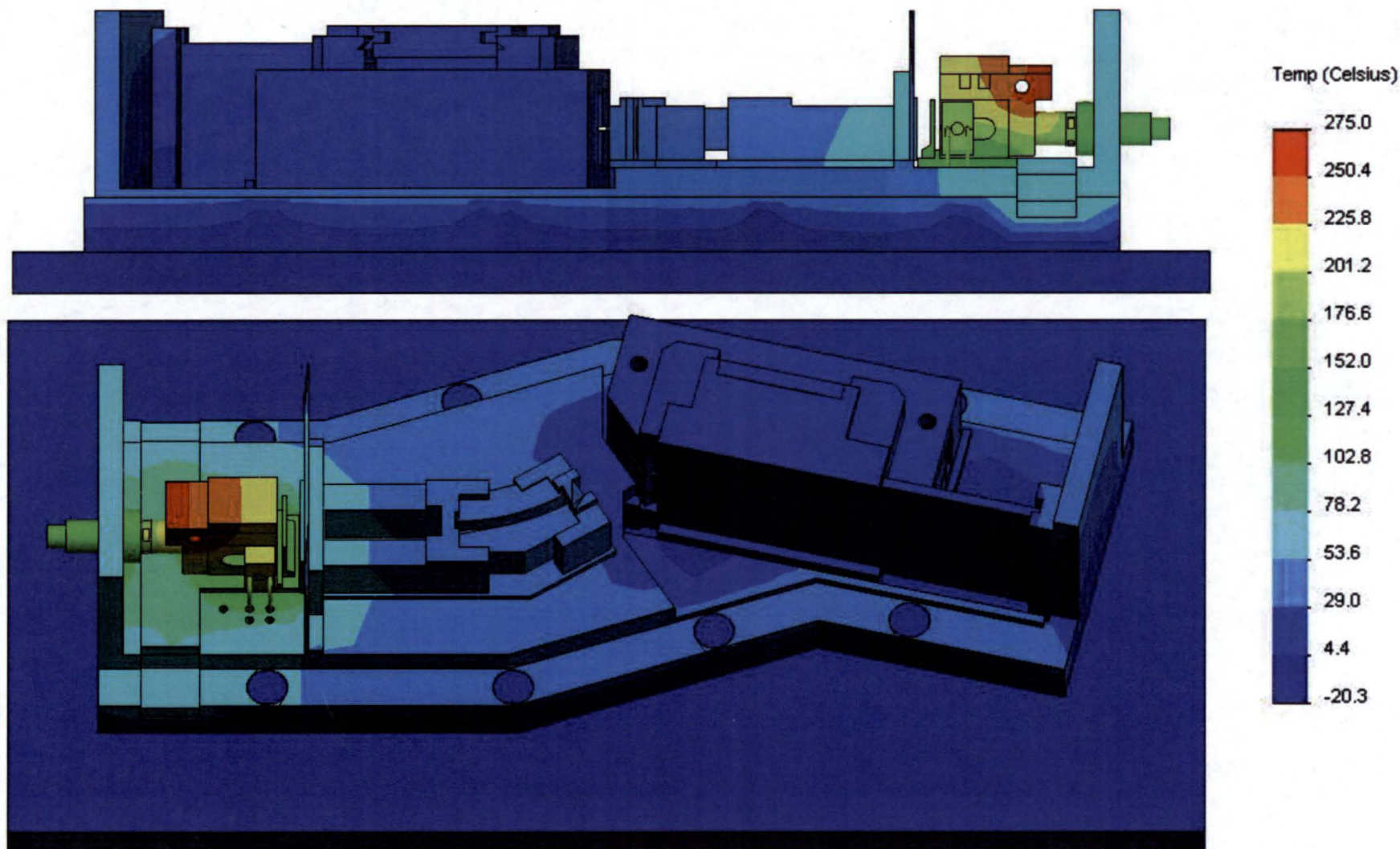




Thermal Study -20°C Ambient



- With tie-bar at -10°C
- Ionization volume temperature virtually unaffected by presence of tie bar under the Detector Assembly

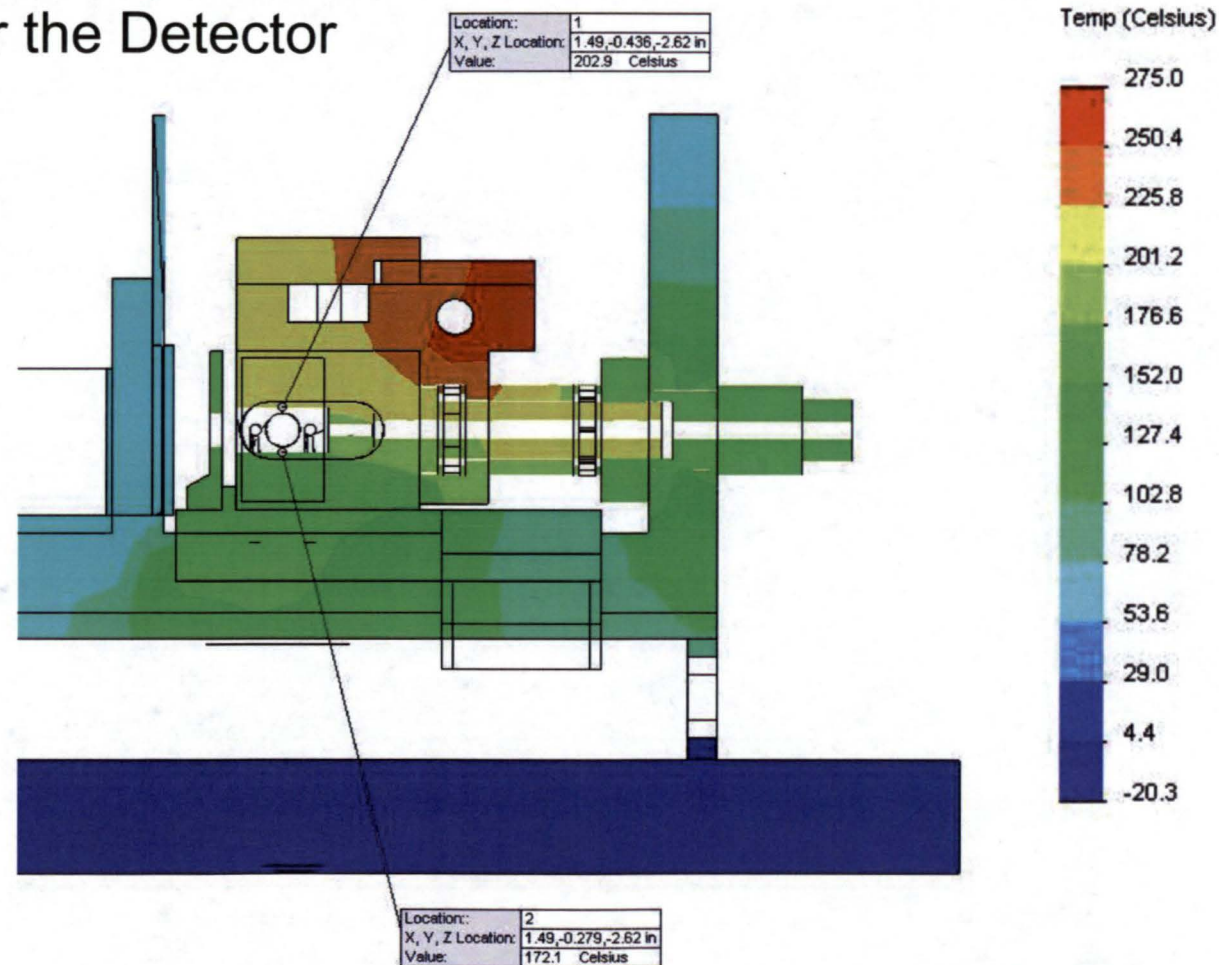




Thermal Study -20°C Ambient



- With NO tie-bar heat sink

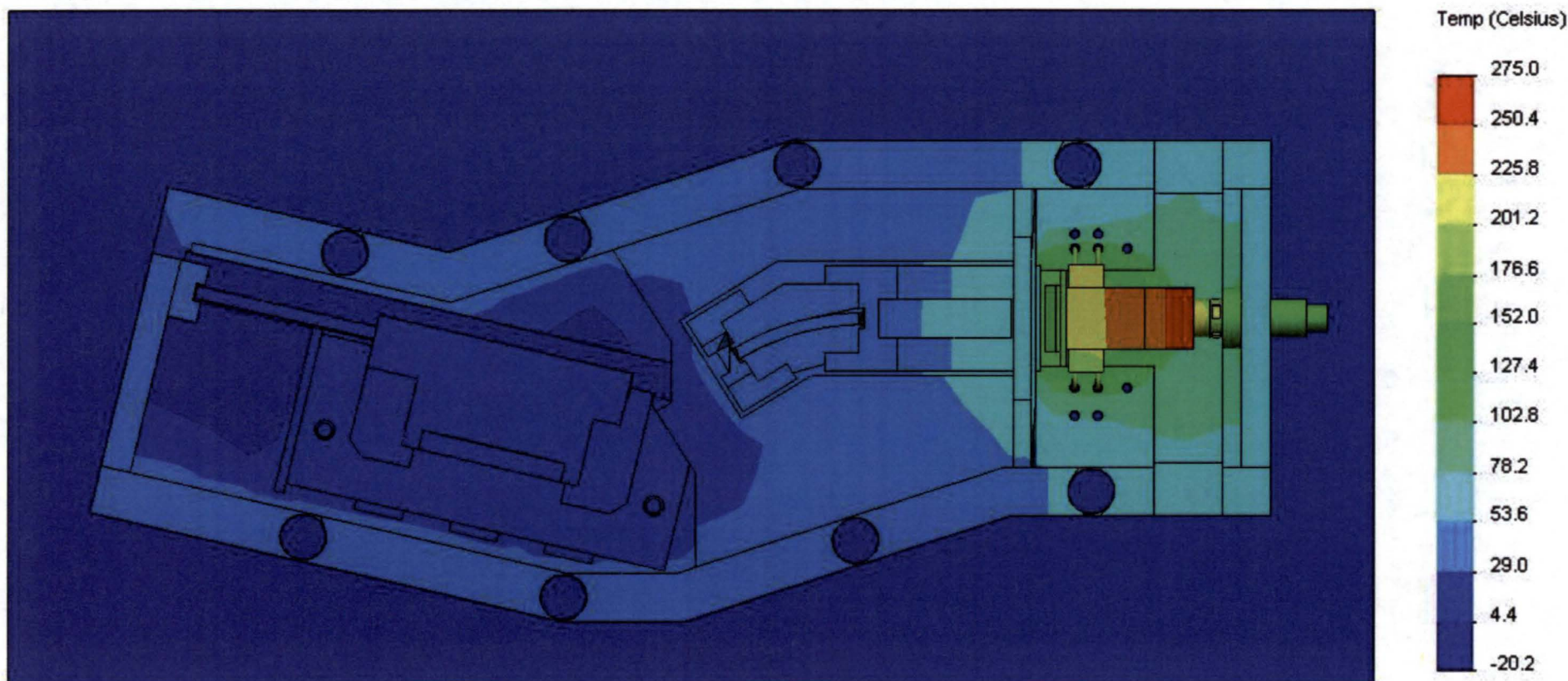




Thermal Study -20°C Ambient



- With NO tie-bar heat sink

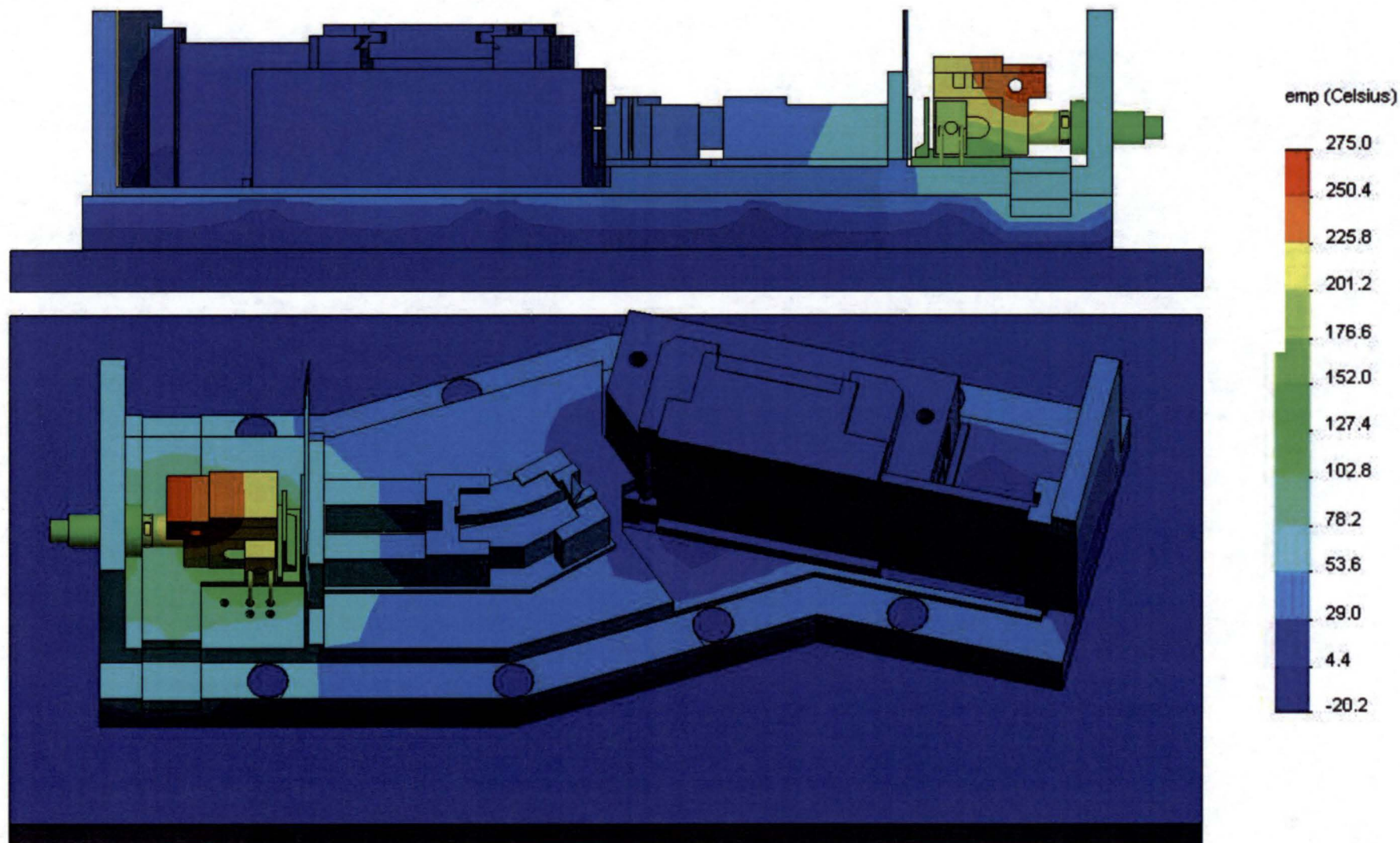



Thermal Study Heat Dissipation



+40 deg C

No Tie Bar		
-Mounting screws: 6.9 W		
	Average:	7069.7
W/m^2		
	Max:	26292
W/m^2		
	Min:	386.58
W/m^2		
-Magnet Tie Bar: 2.1 W		
	Average:	1465.5
W/m^2		
	Max:	10280
W/m^2		
	Min:	0.38284

W/m^2		
-Aluminum Heat Sink Bar: 3.0 W		
	Average:	30829
W/m^2		
	Max:	1.19e+05
W/m^2		
	Min:	4384.6
W/m^2		
Tie Bar		
-Mounting Screw: 5.3 W		
	Average:	5425.1
W/m^2		
	Max:	24545
W/m^2		
	Min:	7563.9

-20 deg C

No Tie Bar		
-Mounting screws: 9.5 W		
	Average:	9745.3 W/m^2
	Max:	36035 W/m^2
	Min:	533.63 W/m^2
-Magnet Tie Bar: 2.8 W		
	Average:	1980.4 W/m^2
	Max:	11563 W/m^2
	Min:	.93043 W/m^2
-Aluminum Heat Sink Bar: 3.9 W		
	Average:	40863 W/m^2
	Max:	1.575e+05
W/m^2		
	Min:	5848W/m^2

Tie Bar		
-Mounting Screw: 9.3 W		
	Average:	9469.7 W/m^2
	Max:	35640 W/m^2
	Min:	522.4 W/m^2
-Magnet Tie Bar: 2.8 W		
	Average:	525.93 W/m^2
	Max:	1989.7 W/m^2
	Min:	1.0277 W/m^2
-Aluminum Heat Sink Bar: 3.9		
	Average:	41018 W/m^2
	Max:	1.5811e+05
W/m^2		
	Min:	5848.9 W/m^2

W/m^2